



Determination of Patulin in Processed Foods Using QuEChERS Extraction and UHPLC-MS/MS Analysis

UCT Part Numbers

ECMSSC-MP

Mylar pouch containing 4g MgSO₄
and 1g NaCl

CUMPSC1875CB2CT

2mL dSPE tube with 150mg
MgSO₄, 50mg PSA, 50mg C18 and
7.5mg GCB

or

ECQUEU615CT

15mL dSPE tube with 900mg
MgSO₄, 150mg PSA and 45mg GCB

SLDA50ID21-18UM

Selectra® DA UHPLC column
(50 × 2.1 mm, 1.8 μm)

SLDAGDC20-18UM

Selectra® DA guard cartridge
(10 × 2.1 mm, 1.8 μm)

SLGRDHLDR

Guard cartridge holder



Summary:

Patulin (Figure 1) is a naturally occurring mycotoxin that is produced by several species of fungi, such as *Aspergillus*, *Penicillium* and *Byssochlamys*. It typically grows on fruit, including apples, pears, peaches and grapes, but has also been reported in vegetables and cereal grains. Patulin has been implicated as a possible carcinogen and teratogen, although an official designation has not yet been made. The main risk arises when unsound fruit is used for the production of juices and other processed food products.

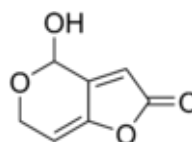


Figure 1. Structure of patulin.

The World Health Organization, U.S. Food and Drug Administration (FDA) and European Union (EU) have suggested a maximum limit of patulin in apple juice and apple juice ingredients at 50 μg/kg. Furthermore, the EU has set a limit of 25 μg/kg in solid apple products and 10 μg/kg in baby food (EC 1881/2006).

This application note outlines a QuEChERS procedure for the detection of patulin in processed food products. An apple-based baby food product was used as a representative sample matrix. The inclusion of graphitized carbon black (GCB) in the dispersive-SPE (dSPE) cleanup step produces a clean sample extract. Analysis is carried out by UHPLC-MS/MS using a Selectra® DA column. The unique chemistry of the Selectra® DA column, which contains a polyaromatic stationary phase, provides a high degree of retention and selectivity for aromatic compounds and improved retention of polar compounds.



FOOD

QuEChERS Procedure:

Sample Extraction:

1. Weigh 10 g of sample into a 50 mL polypropylene centrifuge tube.
2. Add 10 mL of acetonitrile.
3. Add the contents of the **ECMSSC-MP** Mylar pouch and shake for a minimum of 1 minute (by hand or mechanically).
4. Centrifuge the samples at $\geq 3000 \times g$ for 5 minutes.

Sample Clean-up of 1 mL Extract:

1. Transfer 1 ml of supernatant to a 2 mL dSPE tube (**CUMPSC1875CB2CT**).
2. Vortex the sample for 30 seconds.
3. Centrifuge the samples at $\geq 3000 \times g$ for 5 minutes.
4. Transfer 500-600 μL of purified supernatant into an autosampler vial (filter if desired).

Sample Clean-up of 5 mL Extract (for increased sensitivity):

1. Transfer 6 ml of supernatant to a 15 mL dSPE tube (**ECQUEU615CT**).
2. Vortex the sample for 1 minute.
3. Centrifuge the samples at $\geq 3000 \times g$ for 5 minutes.
4. Transfer 5 mL of purified supernatant into a glass tube and evaporate the sample to ≈ 1 mL at 40°C under a gentle stream of nitrogen.
5. Transfer the sample into an autosampler vial (filter if desired).

LC-MS/MS Parameters:

Instrumentation	
HPLC system	Thermo Scientific™ Dionex™ Ultimate™ 3000 UHPLC
MS system	Thermo Scientific™ TSQ Vantage™ (MS/MS)
HPLC column	UCT Selectra® DA, 50 \times 2.1 mm, 1.8 μm (p/n: SLDA50ID21-18UM)
Guard column	UCT Selectra® DA, 10 \times 2.1 mm, 1.8 μm (p/n: SLDAGDC20-18UM)
Guard column holder	p/n: SLGRDHLDLDR
Column temperature	40°C
Flow rate	400 $\mu\text{L}/\text{min}$
Injection volume	5 μL

Mobile Phase Gradient		
Time (min)	% Mobile Phase A Water	% Mobile Phase B Methanol
0.0	95	5
2.0	5	95
3.5	5	95
3.6	95	5
6.0	95	5

Note: no mobile phase additive was used as it was found to give better signal intensity in ESI⁻ for patulin.



MRM transitions (ESI)				
Compound	t _R (min)	Precursor ion	Product ion 1	Product ion 2
Patulin	2.32	153.0	109.0	81.0

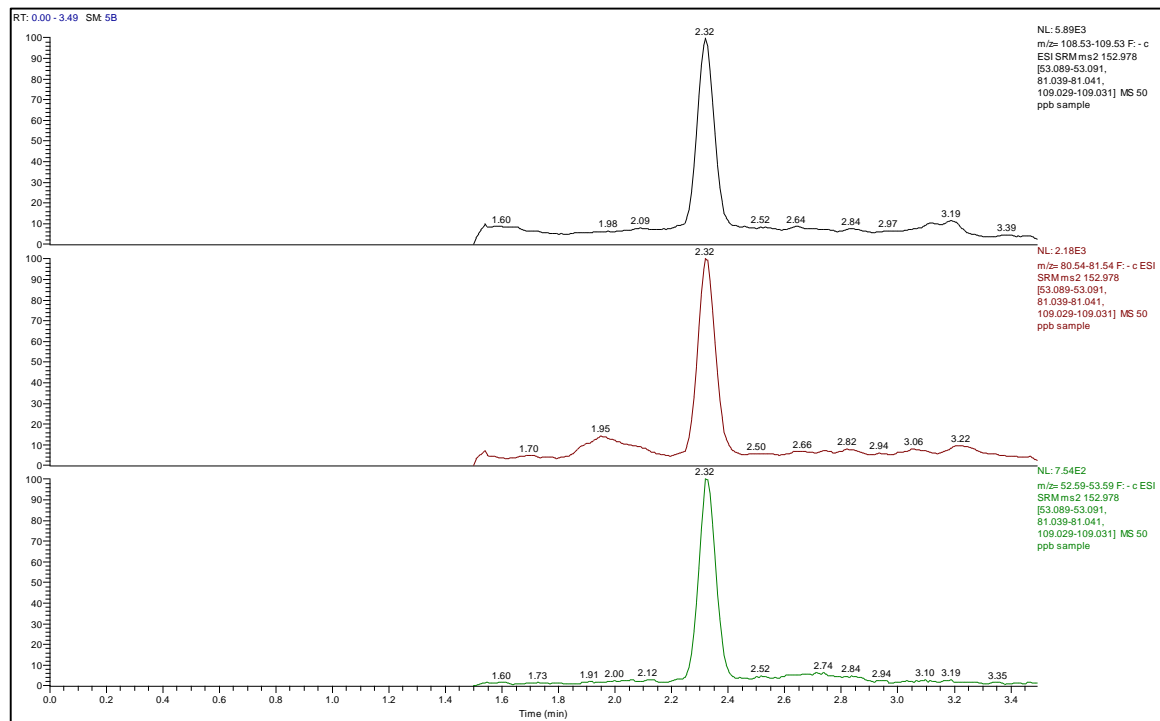



Figure 1. Chromatogram of an extracted baby food sample fortified with patulin at 50 µg/kg.

Results:

Recovery and Reproducibility		
	10 µg/kg ^a	50 µg/kg ^b
Sample 1	102.3	90.1
Sample 2	101.8	83.4
Sample 3	99.2	86.0
Sample 4	103.8	88.3
Sample 5	104.8	87.9
Sample 6	106.1	83.9
Mean Recovery (%)	103.0	86.6
RSD (%)	2.35	0.30

^a used 5 mL of sample extract and included a concentration step.

^b used 1 mL of sample extract and no concentration step.

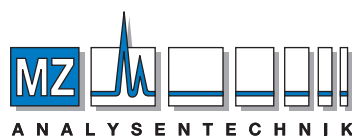

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de

6103-03-01

UCT, LLC • 2731 Bartram Road • Bristol, PA 19007 800.385.3153 • 215.781.9255

www.unitedchem.com Email: methods@unitedchem.com

©UCT, LLC 2014 • All rights reserved

