

ACQ-06 Column Switcher

—Provide Fast Purification

ACQ-06 column switcher allows chemists to switch automatically from one column to another, without evaporating solvents. It also offers multi-stationary phases for a complete chromatography solution.

- Time Saving
Start a new run by switching to a new column
- Efficient Purification
Further purification can be achieved through tandemed columns
- Good Compatibility
Compatible with the purification system from most of the manufacturers in the market




Structure



- 1 Indicator Lamps
Observe the status of the channels
- 2 Column Support
Suitable for columns of different specifications
Enable column tandem
- 3 Standard Adaptor
Compatible with flash columns from other manufacturers
- 4 Channel Switcher
Avoid cross contamination

Software

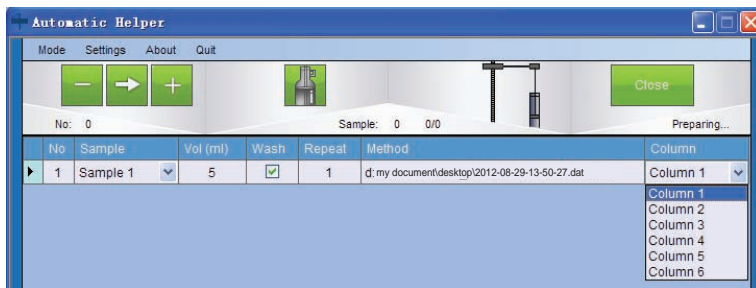

- Easy to Operation
Six channels enable the switching from column 1 to column 6
- Independent Software
Allow single use or connected with CHEETAH program

Still purify your sample by hand?

Free yourself by choosing a fully automatic purification kit !

Just build a connection between

- ATS Auto-Sampler
- ACQ-06 Column Switcher
- CHEETAH™ System

