

Conditions

Column: ACE Excel 1.7 C18
Dimensions: 50 x 3.0 mm
Part Number: EXL-171-0503U
Mobile Phase: 20 mM potassium phosphate
 pH 2.5 in MeCN/H₂O (30:70 v/v)
Flow Rate: 0.43 mL/min
Injection: 0.7 µL
Temperature: 20 °C
Detection: UV, 230 nm


Analytes

1. Phthalic acid
2. 4-Hydroxybenzoic acid
3. Sorbic acid
4. Methyl paraben
5. Ethyl paraben
6. Propyl paraben




AUTHORIZED DISTRIBUTOR


MZ-Analysentechnik GmbH
 Barcelona-Allee 17 • D-55129 Mainz
 Tel +49 6131 880 96-0
 Fax +49 6131 880 96-20
 e-mail: info@mz-at.de
 www.mz-at.de

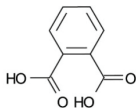



Phthalic acid

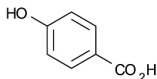


4-Hydroxybenzoic acid

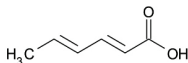


Sorbic acid

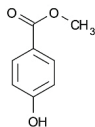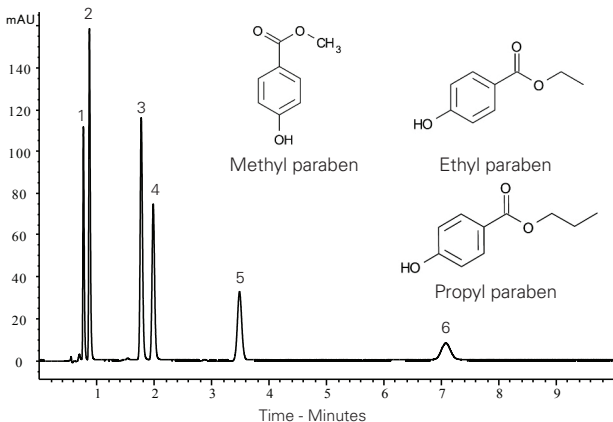

Methyl paraben

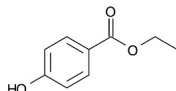

Ethyl paraben

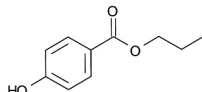

Propyl paraben