


Conditions

Column: ACE 5 C18
Dimensions: 250 x 4.6 mm
Part Number: ACE-121-2546
Mobile Phase: 50 mM KH_2PO_4 pH 5.7 in
 $\text{H}_2\text{O}/\text{MeOH}$ (70:30 v/v)
Flow Rate: 1 mL/min
Temperature: 22 °C
Detection: UV, 220 nm

Analytes


1. L-Ascorbic acid
2. Maleic acid
3. Acetylsalicylic acid
4. Benzoic acid
5. Salicylic acid

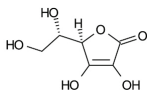

L-Ascorbic acid



Maleic acid

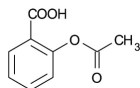

Acetylsalicylic acid



Benzoic acid



Salicylic acid

