



Methylmalonic Acid and Succinic Acid using HILIC Mode

ACE[®]
Ultra-inert
UHPLC & HPLC Columns


Application #AN4520

Conditions


Column: ACE 5 HILIC-B
Dimensions: 150 x 4.6 mm
Part Number: HILB-5-1546U
Mobile Phase: 10 mM ammonium formate pH 3.0 in MeCN/H₂O (90:10 v/v)
Flow Rate: 1.5 mL/min
Injection: 5 µL
Temperature: 25 °C
Detection: ELSD
Evaporator temperature: 30 °C
Nebuliser temperature: 30 °C
Gas speed: 1 SLM

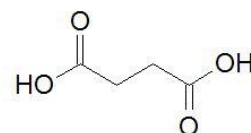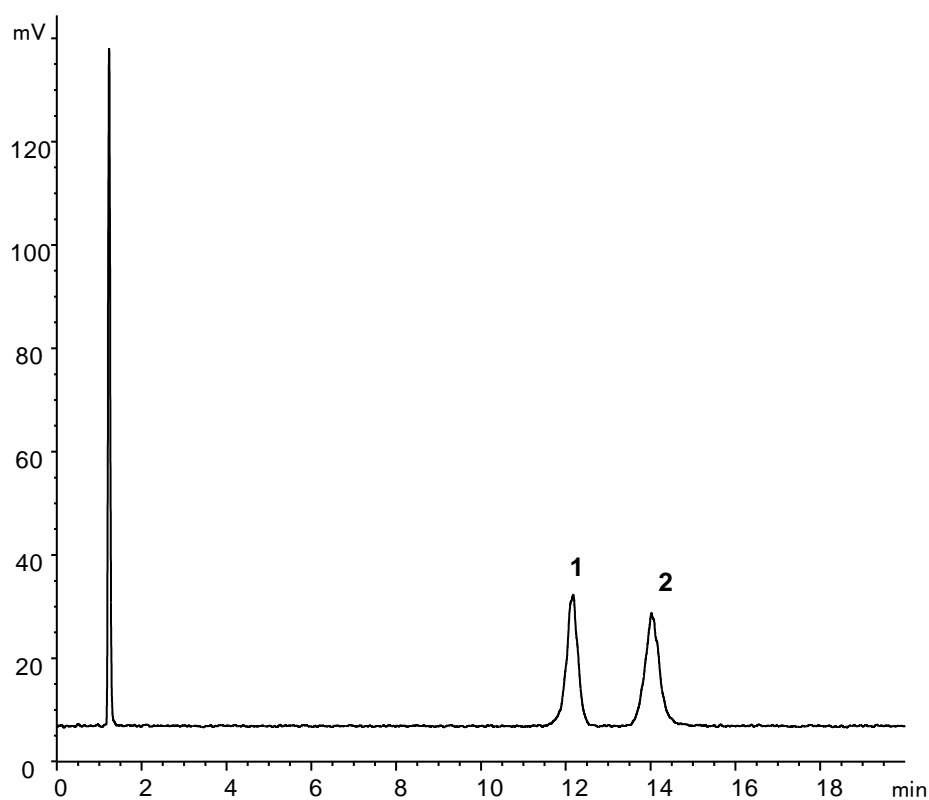

1. Succinic acid

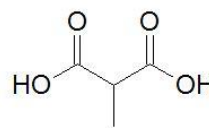

2. Methylmalonic acid

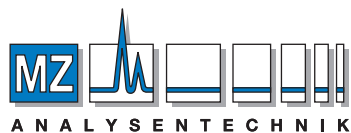

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Wöhlerstraße 2-6 • D-55120 Mainz
Tel +49 6131 68 66 19, Fax +49 6131 68 66 20
e-mail: info@mz-at.de, www.mz-at.de

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.

