



分岐鎖アミノ酸の分離

Branched-chain amino acids SunArmor NH2 5 µm, 250 x 4.6 mm i.d.


Temperature: 25 °C

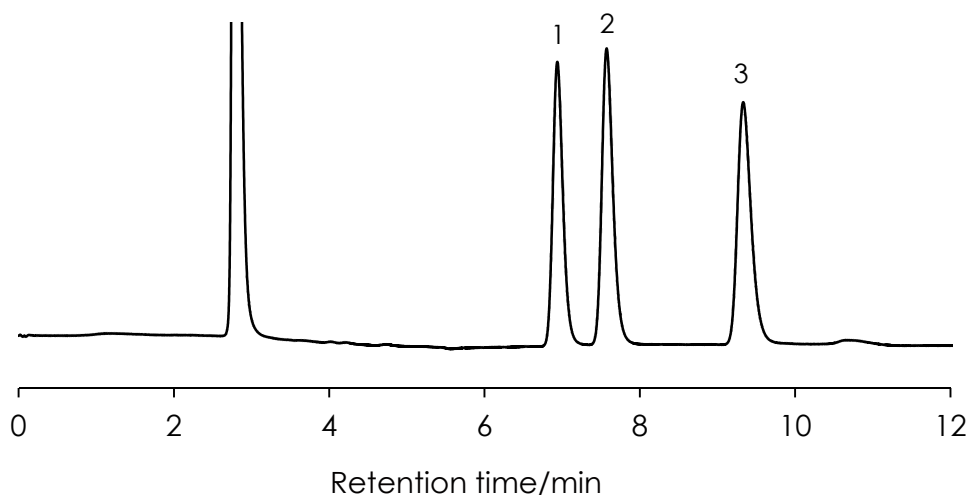


1. L-Leucine

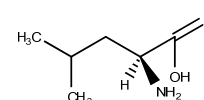


2. L-Isoleucine

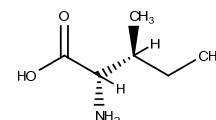

3. L-Valine

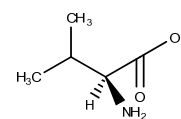

Temperature: 40 °C

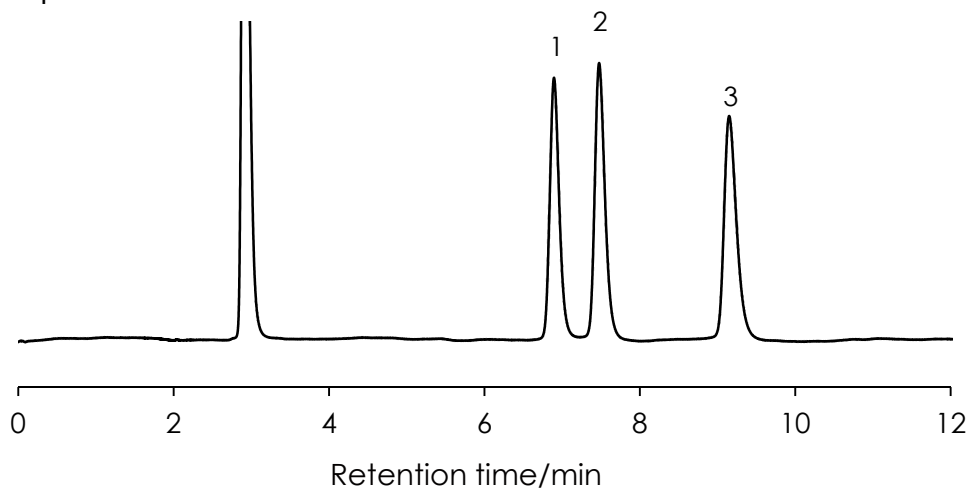

These samples have little temperature dependence.

Column: SunArmor NH2 5 µm, 250 x 4.6 mm


Mobile phase: Acetonitrile:10mM ammonium acetate = 70:30

Flow rate: 1.0 mL/min

Temperature: 25, 40 °C

Detection: RI

Sample: 1 = L-Leucine, 2 = L-Isoleucine, 3 = L-Valine

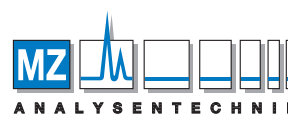

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de