


OPAとFMOC誘導体化アミノ酸の分離(1) Sunniest C18-HT 2 μ m, 100 x 2.1 mm i.d.

Amino Acids derivatized with OPA and FMOC (1)

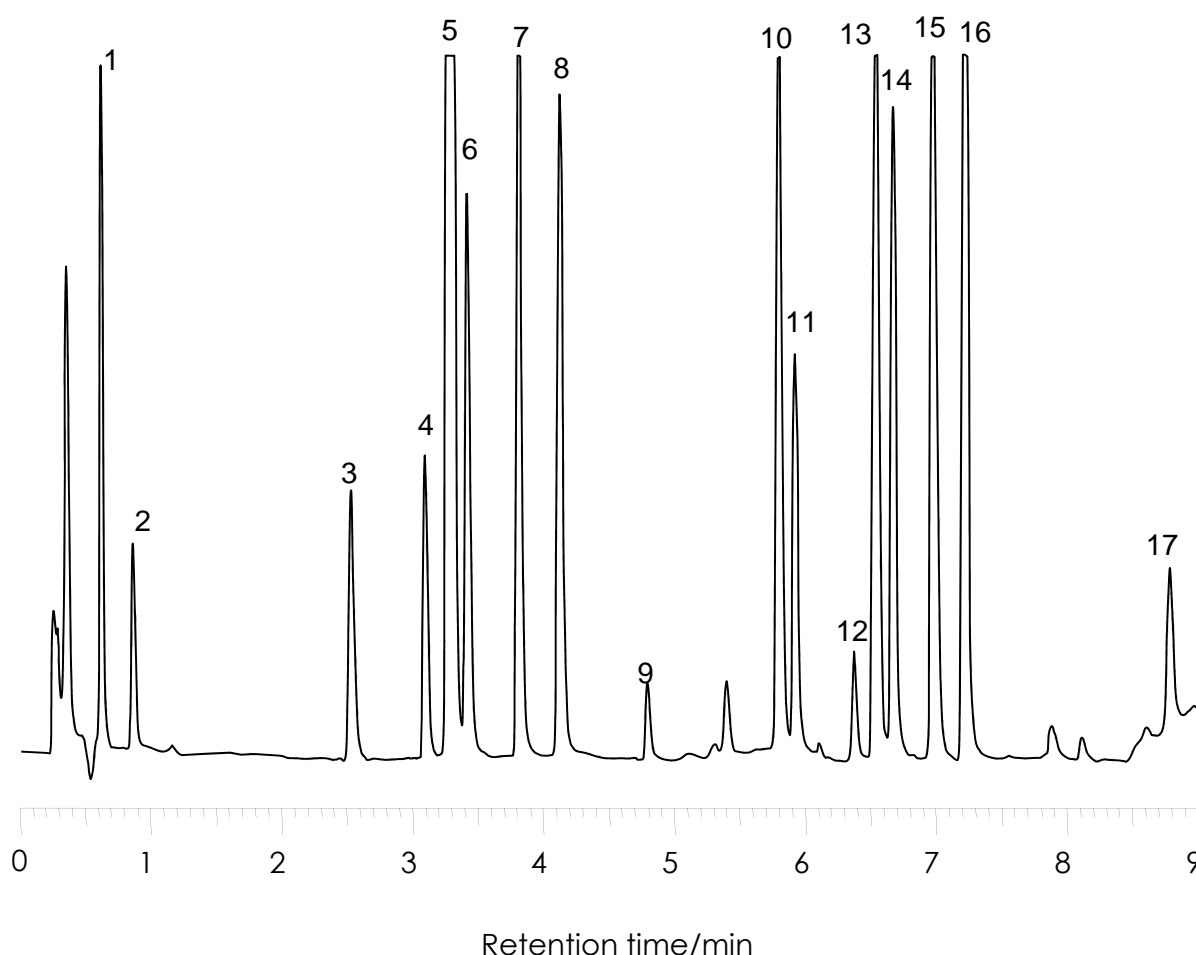


Column: Sunniest C18-HT 2 μ m, 100 x 2.1 mm

Mobile phase: A) 10mM Na₂PO₄ + 10mM Na₂B₄O₇ + 0.5mM NaN₃ (pH7.8)

B) Acetonitrile/Methanol/Water (45/45/10 %V)

Time(min)	0	0.2	7.2	7.8
%B	5	5	50	100

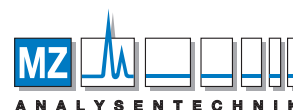

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de

Flow rate: 0.72 mL/min

Temperature: 40 °C

Detection: UV@338 nm

Sample: 1=Aspartic acid, 2=Glutamic acid, 3=Serine, 4=Histidine, 5=Glycine, 6=Threonine, 7=Arginine, 8=Alanine, 9=Tyrosine, 10=Valine, 11=Methionine, 12=Tryptophan, 13=Proline, 14=Isoleucine, 15=Leucine, 16=Lysine, 17=Proline