
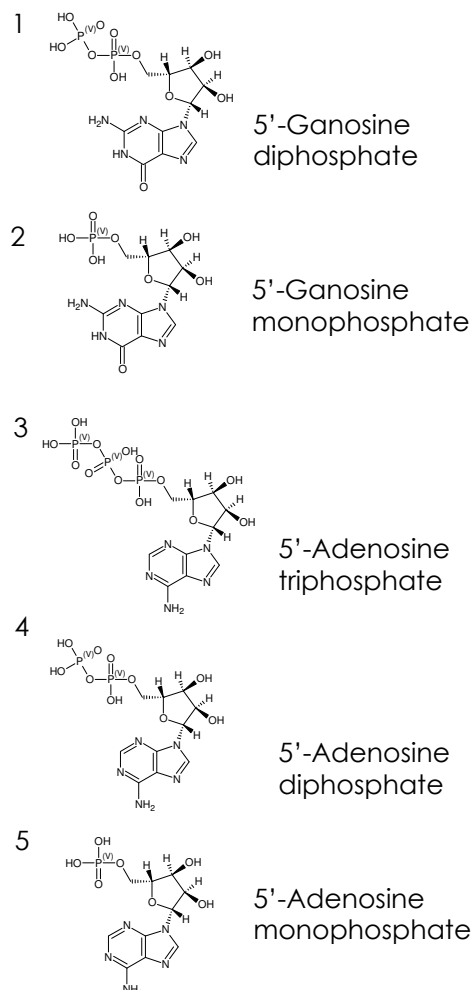


## ヌクレオチドの分離 (1)

### Nucleotides (1)

Sunniest RP-AQUA 3 µm, 150 x 4.6 mm i.d.

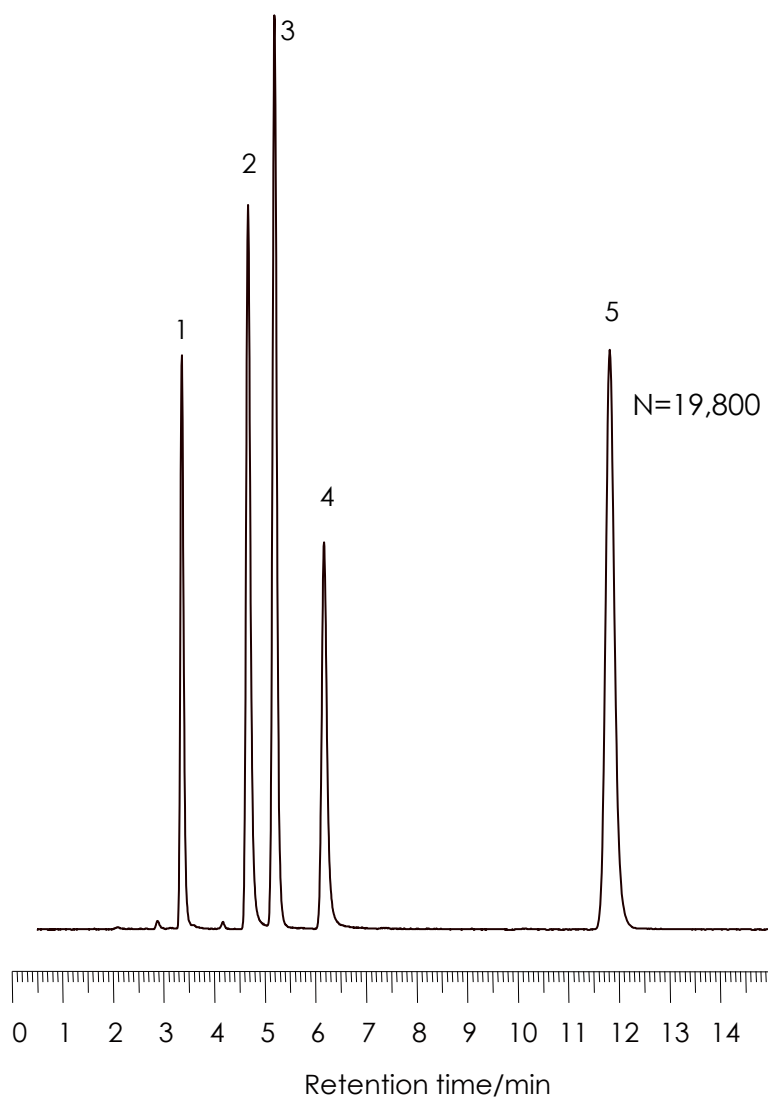

Column: Sunniest RP-AQUA 3 µm, 150 x 4.6 mm


Mobile phase: 20mM Phosphate buffer pH6.0

Flow rate: 1.0 mL/min

Temperature: 40 °C

Detection: UV@250nm

Sample: 1 = 5'-GDP, 2 = 5'-GMP, 3 = 5'-ATP, 4 = 5'-ADP, 5 = 5'-AMP

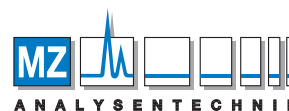

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH  
Barcelona-Allee 17 • D-55129 Mainz  
Tel +49 6131 880 96-0  
Fax +49 6131 880 96-20  
e-mail: info@mz-at.de  
www.mz-at.de