


ヌクレオチドの分離 (2)

Nucleotides (2)

SunShell RP-AQUA 2.6 μ m, 150 x 4.6 mm i.d.

SunShell RP-AQUA

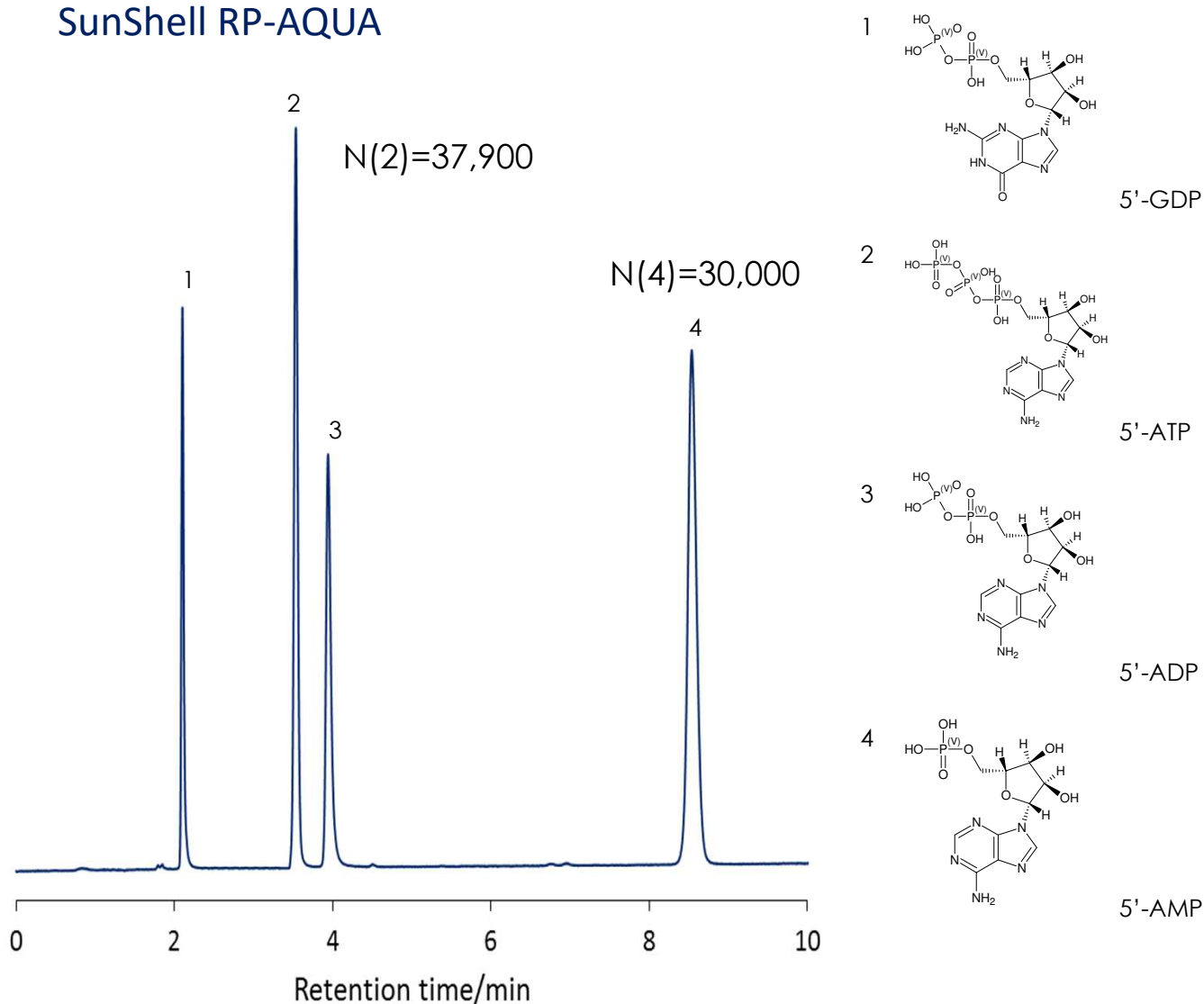

Column: SunShell RP-AQUA 2.6 μ m, 150 x 4.6 mm i.d.

Mobile phase: 20mM Phosphate buffer pH6.0


Flow rate: 1.0 mL/min

Temperature: 25 $^{\circ}$ C

Detection: UV@250nm

Sample: 1 = 5'-GDP, 2 = 5'-ATP, 3 = 5'-ADP, 4 = 5'-AMP

UHPLC: X-LC (Jasco)

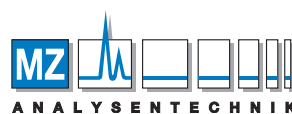

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de