
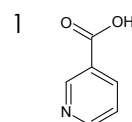


## 水溶性ビタミン類Ⅲ


## Water-soluble vitamins III

SunShell RP-AQUA 2.6  $\mu\text{m}$ , 150 x 4.6 mm i.d.


### SunShell RP-AQUA

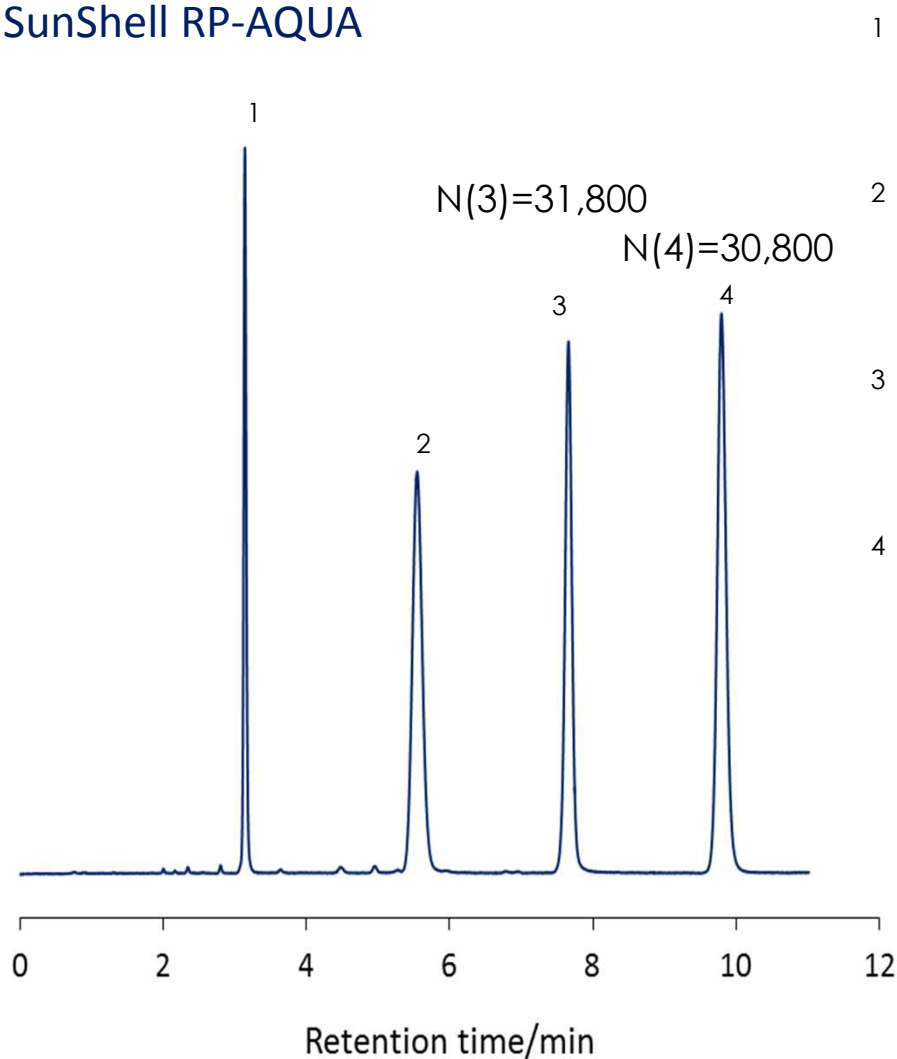


Nicotinic acid

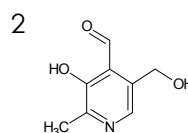

Pyridoxal

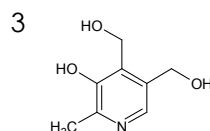

Pyridoxine HCl  
Vitamin B<sub>6</sub>

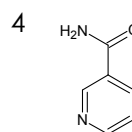

Nicotinamide

Column: SunShell RP-AQUA 2.6  $\mu\text{m}$ , 150 x 4.6 mm i.d.

Mobile phase: 40mM Phosphate buffer pH6.8


Flow rate: 1.0 mL/min

Temperature: 40 °C

Detection: UV@250 nm

Sample: 1 = Nicotinic acid, 2 = Pyridoxal, 3 = Pyridoxine HCl, 4 = Nicotinamide

UHPLC: X-LC (Jasco)

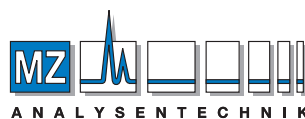

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH  
Barcelona-Allee 17 • D-55129 Mainz  
Tel +49 6131 880 96-0  
Fax +49 6131 880 96-20  
e-mail: info@mz-at.de  
www.mz-at.de