
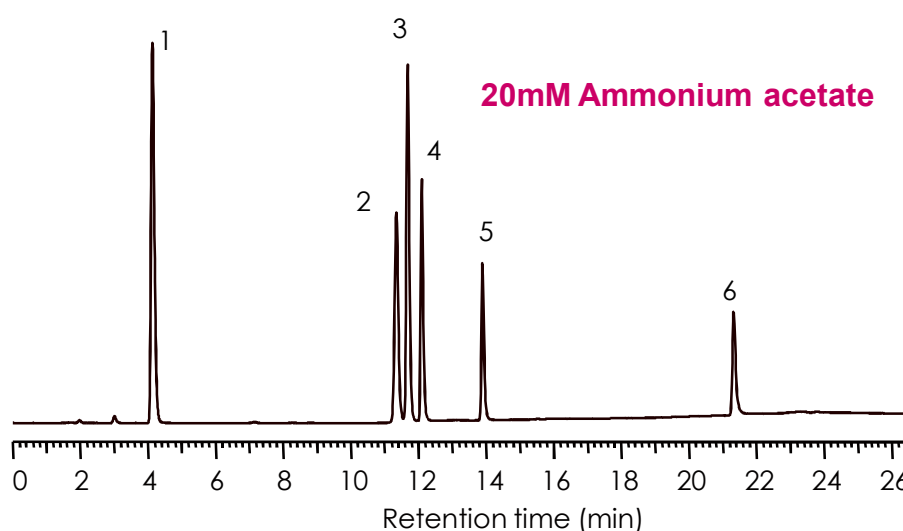
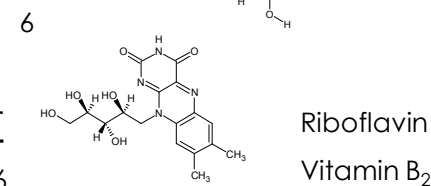
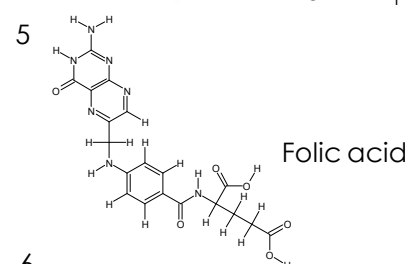
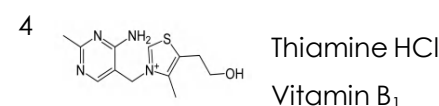
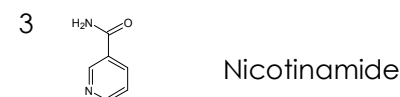
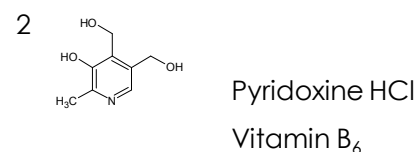


水溶性ビタミン類

Water-soluble vitamins

Sunniest RP-AQUA 5µm, 150 x 4.6 mm i.d.

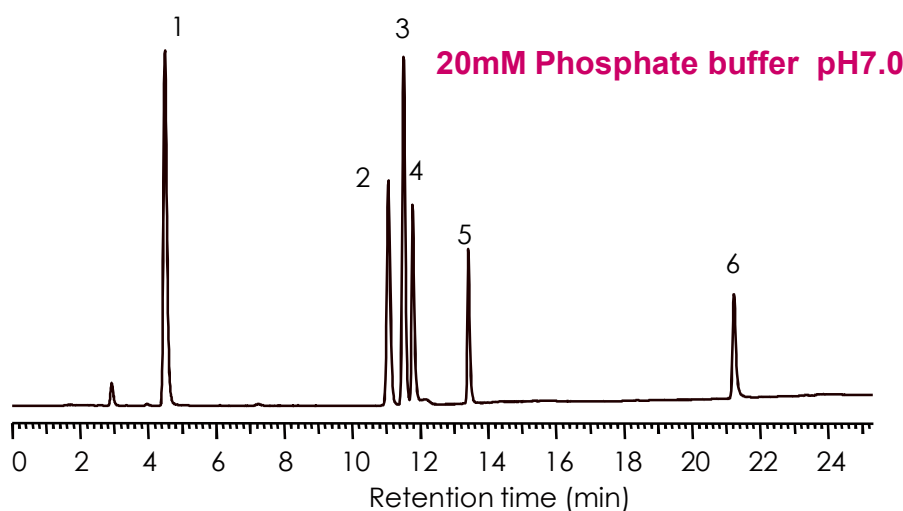

Column: Sunniest RP-AQUA 5 µm, 150 x 4.6 mm i.d.

Mobile phase: A) 20mM Phosphate buffer pH7.0 or 20mM Ammonium acetate

B) Acetonitrile/ A solution (20:80)


| | | | | |
|------------|---|---|-----|-----|
| Time (min) | 0 | 5 | 20 | 25 |
| %B | 0 | 0 | 100 | 100 |

Flow rate: 1.0 mL/min

Temperature: 40 °C

Detection: UV@250 nm

Injection volume: 10 µL

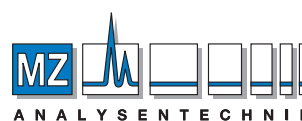

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de