

PART OF THE DURAN® BOTTLE SYSTEM

# SWIRLED

DURAN® BAFFLED BOTTLE



**DURAN  
WHEATON  
KIMBLE**

Excellence in your hands

DURAN® GLS 80® BAFFLED BOTTLE WIDE MOUTH




with GLS 80® thread



During mixing in standard DURAN® GLS 80® laboratory bottles, the liquid moves in a swirling motion that approximates a solid-body rotation. This is a very inefficient flow pattern, and very little mixing actually occurs. By adding three vertical baffles into the inner surface of the DURAN® GLS 80® bottles, the swirling motion is disrupted and an improvement of the top-to-bottom circulation occurs which produces a greater radial and more effective mixing.

Cat. No.	Capacity (mL)	Thread	d (OD) (mm)	h (mm)	Pack Unit
<b>with screw cap and pouring ring from PP (blue)</b>					
21 286 36 58	250	80	95	110	1
21 286 44 57	500	80	101	153	1
21 286 54 53	1 000	80	101	223	1

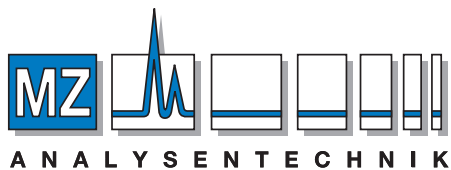

**AUTHORIZED DISTRIBUTOR**

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz  
 Tel +49 6131 880 96-0, Fax +49 6131 880 96-20  
 e-mail: info@mz-at.de, www.mz-at.de



> Find your nearest distributor on our global network:  
[www.DWK-LifeSciences.com/DURAN/distributors](http://www.DWK-LifeSciences.com/DURAN/distributors)

