



CLICK  
AND  
CARRY

THE NEW  
**DURAN**<sup>®</sup>  
GL 45 BOTTLE  
CARRYING SYSTEM

# DURAN®

## GL 45 BOTTLE CARRYING SYSTEM

The new modular Bottle Carrying System fits standard laboratory bottles that have a GL 45 neck finish and is compatible with four bottle sizes 2, 5, 10 and 20 litres. It is ergonomically designed, and offers health & safety benefits to the user and allows easy and safe transport of bottles around the laboratory or production facility. It makes pouring liquids from the larger bottles more comfortable and safer. All materials are fully autoclavable.

### Carrying Handle

- Carrying handle is suitable for bare or gloved hands
- Handle strap is flexible enough to allow pouring

### Tightening Strap

- With indicator arrow
- End modified to prevent unthreading

### Pouring Handle

- Ergonomically shaped
- Anti-slip grip

### GL 45 Bottle Collar


- The collar is securely attached by pushing down over the GL 45 thread with 4 teeth that lock under transfer ring.
- The collar is only removable by cutting – so attachment is 100% secure in use. A cut line indicator for removal of collar when required.

### Horizontal Straps

- Durable stitched construction
- Made from "anti-slip" material



**Lot Specific Retrace Code**  
for traceability and certification

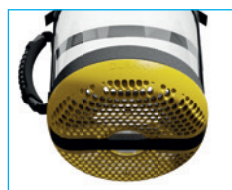


### Silicone Support Base

- Structural support for base of larger 5, 10 and 20 litre bottles
- Protection of glass bottle base from impact during lifting or handling
- Spacer to protect glass from contact with Metal Dolly
- Moulded-in holes minimize liquid traps and allow maximum steam penetration during autoclaving



### Optional Extra: Metal Dolly

- Stainless steel dolly for 10 or 20 litre bottles (2 sizes only)
- Fully compatible with carrying system
- Easy to clean, crevice-free construction
- Optional handle is available for easier handling

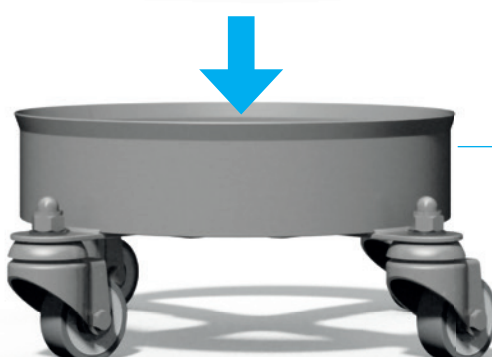

# FOR EASY AND SAFE CARRYING & POURING



## Easy and safe carrying using handle

- Handle attaches at the neck and fully supports the base
- For 2, 5, 10 and 20 litre bottles



## Safe and controlled pouring using handle

- Safer pouring of liquids using a second lower "Pouring Handle" located on the side strap
- Available on the 5 litre carriers only

## Ordering information

Description	Pack Unit	Cat. No.
DURAN® GL 45 Bottle Carrying System for the 2 Litre DURAN® Bottle (Blue)	1	292836307
DURAN® GL 45 Bottle Carrying System for the 2 Litre DURAN® Bottle (Purple)	1	292826306
DURAN® GL 45 Bottle Carrying System for the 2 Litre DURAN® Bottle (Red)	1	292816305
DURAN® GL 45 Bottle Carrying System for the 2 Litre DURAN® Bottle (Yellow)	1	292806304
DURAN® GL 45 Bottle Carrying System for the 5 Litre DURAN® Bottle (Yellow)	1	292807309
DURAN® GL 45 Bottle Carrying System for the 10 Litre DURAN® Bottle (Yellow)	1	292808605
DURAN® GL 45 Bottle Carrying System for the 20 Litre DURAN® Bottle (Yellow)	1	292809104
DURAN® Replacement Handle and Collar for the GL 45 Bottle Carrying System for 5 Litre Bottles	1	292847304
DURAN® Replacement Handle and Collar for the GL 45 Bottle Carrying System for 10 and 20 Litre Bottles	1	292848609
DURAN® Metal Dolly for 10 L DURAN® GL 45 Bottle (Stainless steel)	1	299018607
DURAN® Handle for the 10 L Metal Dolly (Stainless steel)	1	299038609
DURAN® Silicone Bottle Base Protector for the 10 L Metal Dolly (Yellow)	1	292808654
DURAN® Metal Dolly for 20 L DURAN® GL 45 Bottle (Stainless steel)	1	299019106
DURAN® Handle for the 20 L Metal Dolly (Stainless steel)	1	299039108
DURAN® Silicone Bottle Base Protector for the 20 L Metal Dolly (Yellow)	1	292809153

Bottles supplied separately

## Fast and secure tightening

of the polypropylene straps with two side strap adjusters for the larger 5, 10 and 20 litre bottles

4 sizes available:  
2 L / 5 L / 10 L / 20 L



## 2 L bottle carrier in 4 colors

for color coding






**DURAN  
WHEATON  
KIMBLE**


Excellence in your hands

MAX. 1 Litre  
140 °C  
DWK Life Science GmbH  
Otto-Schott-Str. 21  
97877 Wertheim/Main  
Germany  
www.DWK-lifesciences.com

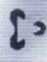RETRACE CODE  
 DURAN  
WHEATON  
KIMBLE

**DWK Life Sciences GmbH**  
Hattenbergstraße 10  
55122 Mainz, Germany

Phone: +49 6131 - 1445 4131  
Fax: +49 6131 - 1445 4016  
sales@DWK.com  
www.DWK.com

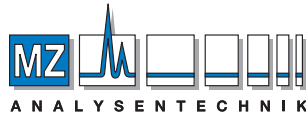

**AUTHORIZED DISTRIBUTOR**

MZ-Analysentechnik GmbH  
Barcelona-Allee 17 • D-55129 Mainz  
Tel +49 6131 880 96-0  
Fax +49 6131 880 96-20  
e-mail: info@mz-at.de  
www.mz-at.de