


InertSearch™ for LC

Inertsil® Applications

Free Amino acids and Organic acids

Data No. LA010-0000

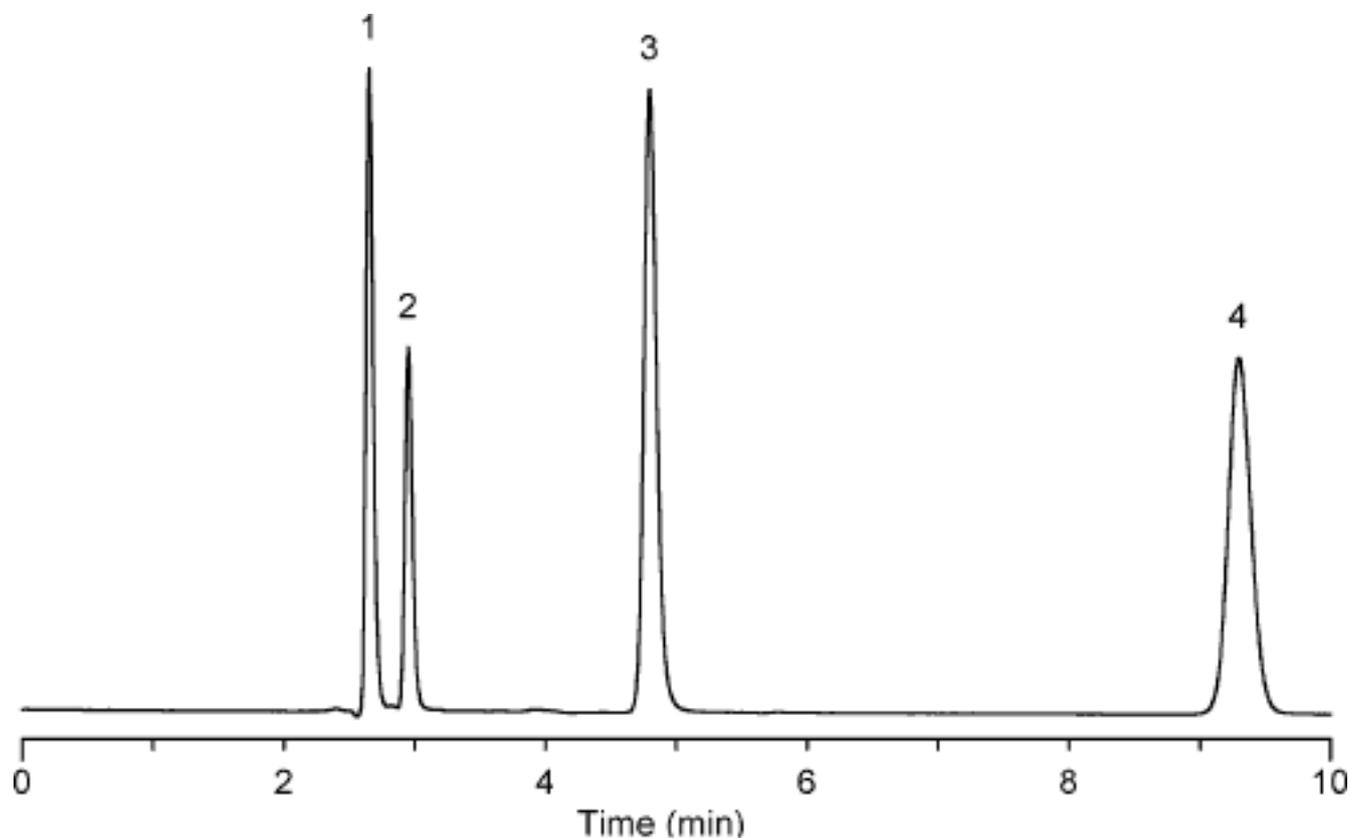

Column : Inertsil CN-3
5µm
250 × 4.6 mm I.D.

Column Cat. No. :5020-01941

Eluent : 20 mM KH₂PO₄ (pH 4.0; H₃PO₄)

Flow rate : 1.0 mL/min.


Col. Temp. : 40°C

Detection : UV 210 nm

Injection Vol. : 5 µL

Sample : Free amino acids

Solutes : 1. Asparagine · H₂O (0.75 mg/mL)
2. Aspartic acid (0.75 mg/mL)
3. Fumaric acid (0.01 mg/mL)
4. Maleic acid (0.01 mg/mL)

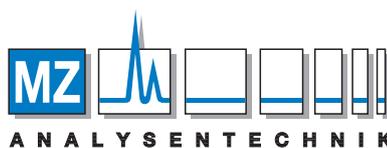

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de