



PTH-Amino acids

Data No. LA002-0000

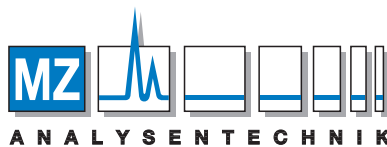

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de

AUTHORIZED DISTRIBUTOR

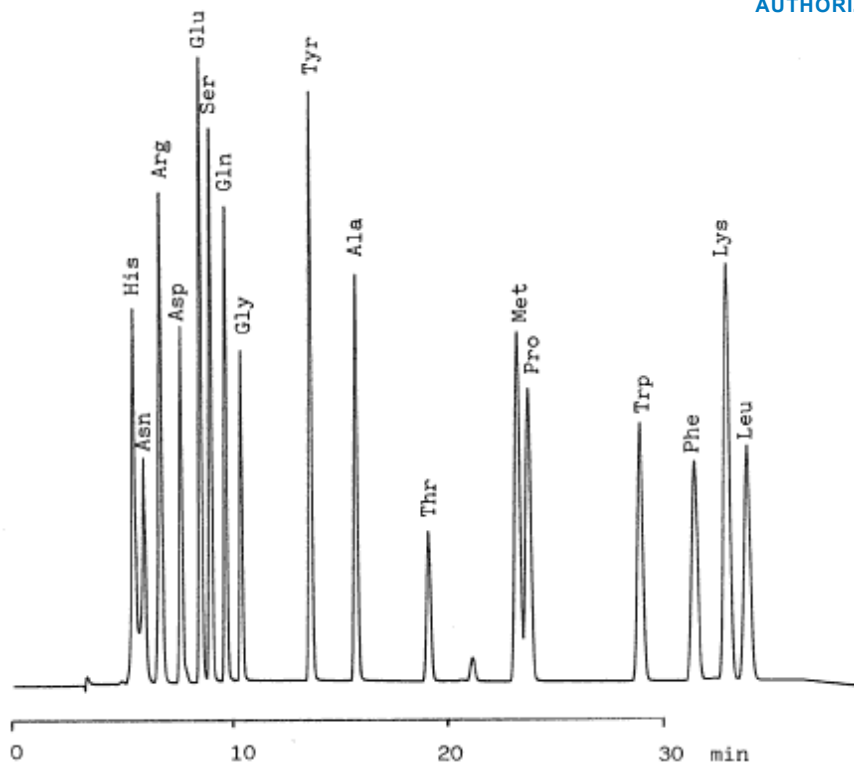

Column : Inertsil ODS-2
5µm
250 × 4.6mm I.D.

Column Cat. No. : 5020-01102

Eluent : A: CH₃CN
B: H₂O

A/B = 20/80 - 5 min - 30/70 - 25 min - 40/60

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254nm

Sample : PTH-amino acids

Solutes

- : 1. Histidine (His)
2. Asparagine (Asn)
3. Arginine (Arg)
4. Aspartic acid (Asp)
5. Glutamic acid (Glu)
6. Serine (Ser)
7. Glutamine (Gln)
8. Glycine (Gly)
9. Tyrosine (Tyr)
10. Alanine (Ala)
11. Threonine (Thr)
12. Methionine (Met)
13. Proline (Pro)
14. Tryptophan (Trp)
15. Phenylalanine (Phe)
16. Lysine (Lys)
17. Leucine (Leu)