


InertSearch™ for LC

Inertsil® Applications

Analysis of 5 Nucleotides (Inertsil AX)

Data No. LA735-0489

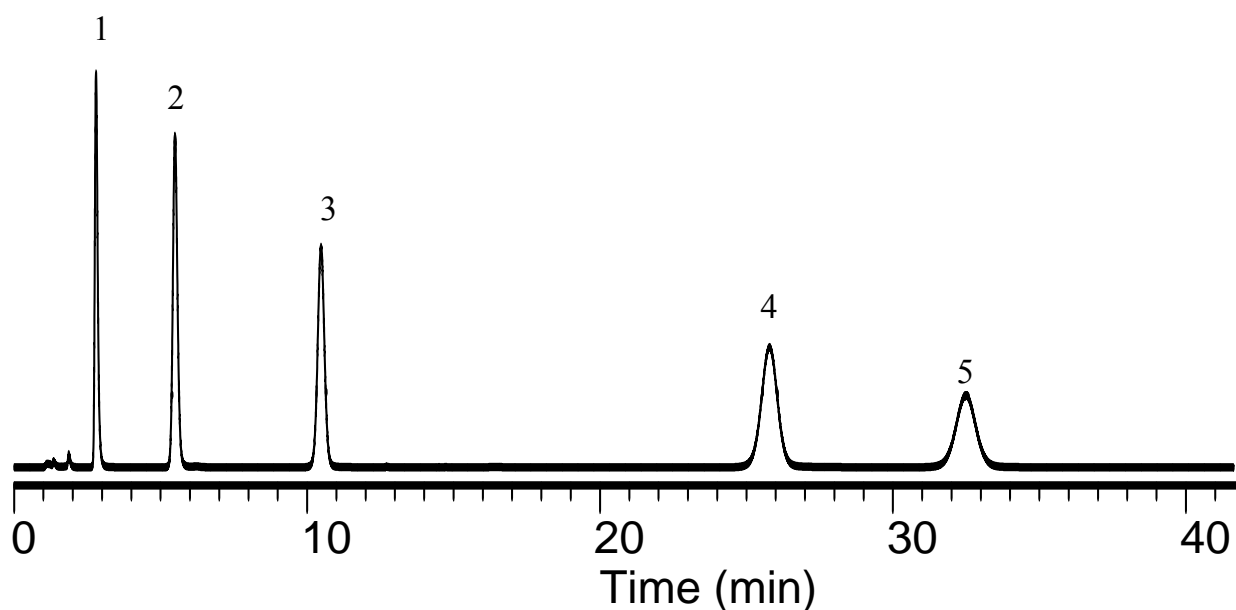


Conditions

Column : Inertsil AX (5 µm, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-07245
Eluent : 0.03 M NH₄H₂PO₄ (pH 3.0, H₃PO₄)
Flow Rate : 1.0 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm
Injection Vol. : 10 µL

Analyte

1. 5'-CMP (0.2 mg/mL)
2. 5'-AMP (0.1 mg/mL)
3. 5'-UMP (0.2 mg/mL)
4. 5'-GMP (0.2 mg/mL)
5. 5'-XMP (0.2 mg/mL)

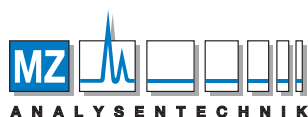

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de