


InertSearch™ for LC

Inertsil® Applications

Analysis of Nucleotides

Data No. LA690-0000

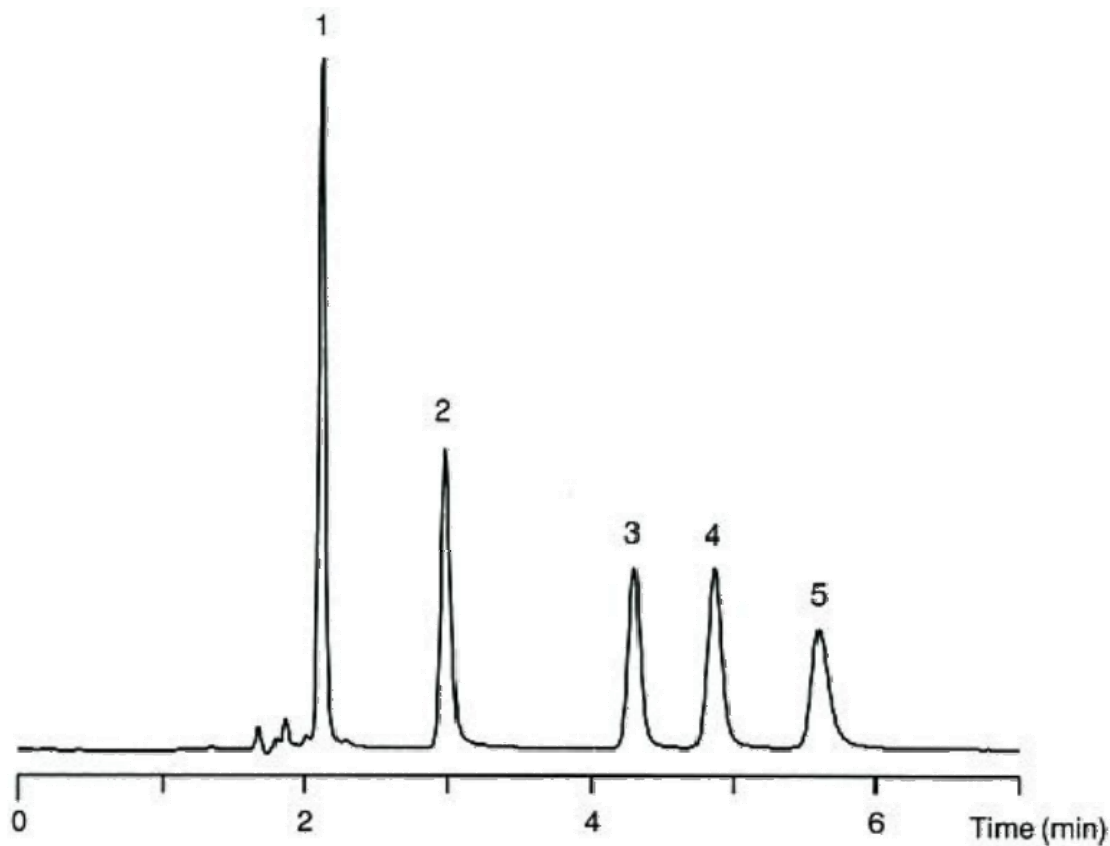


Conditions

Column : Inertsil NH2 (5 µm, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-05545
Eluent : 25 mM Phosphate bufer (pH 2)
Flow Rate : 1.0 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm

Analyte

1. 5'-AMP
2. 5'-GMP
3. 5'-UMP
4. 5'-IMP
5. 5'-TMP

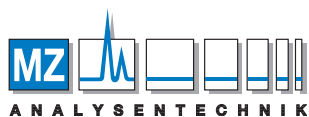

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0, Fax +49 6131 880 96-20
e-mail: info@mz-at.de, www.mz-at.de