


InertSearch™ for LC


Inertsil® Applications

Nucleotides

Data No. LA226-0000

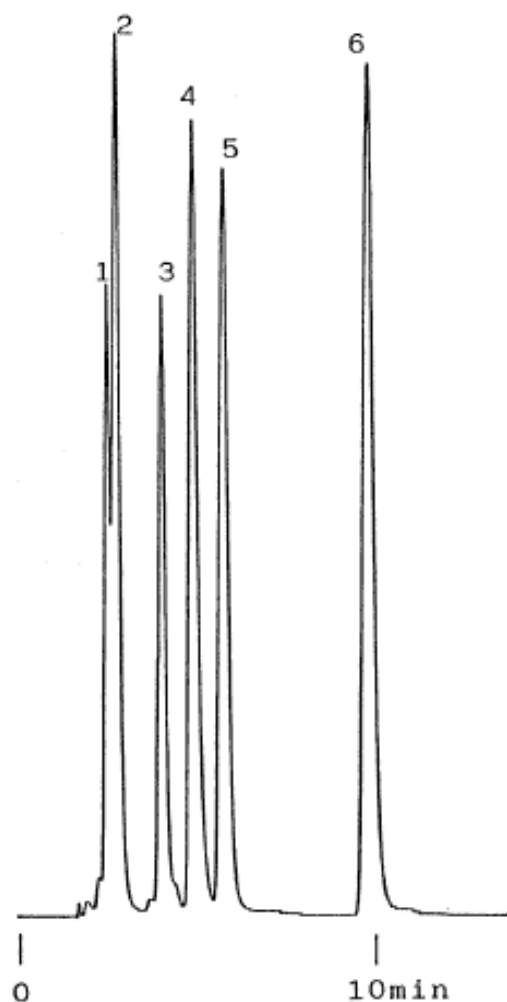

Column	:Inertsil ODS-2 5µm 150 × 0.7 mmI.D.	Solutes	: 1. GTP 2. GDP 3. GMP 4. ATP 5. ADP 6. AMP
Column Cat. No.	:5020-13073		
Eluent	:0.1M Na ₂ HPO ₄ (pH6 by H ₃ PO ₄)		
Flow rate	:30µL/min.		
Col. Temp.	:		
Detection	:UV 260 nm		
Injection Vol.	:		
Sample	:		

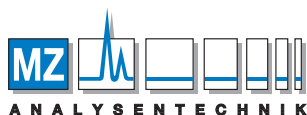

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de

GL Sciences