


InertSearch™ for LC


Inertsil® Applications

Analysis of Neonicotinoids


Data No. LB506-0848

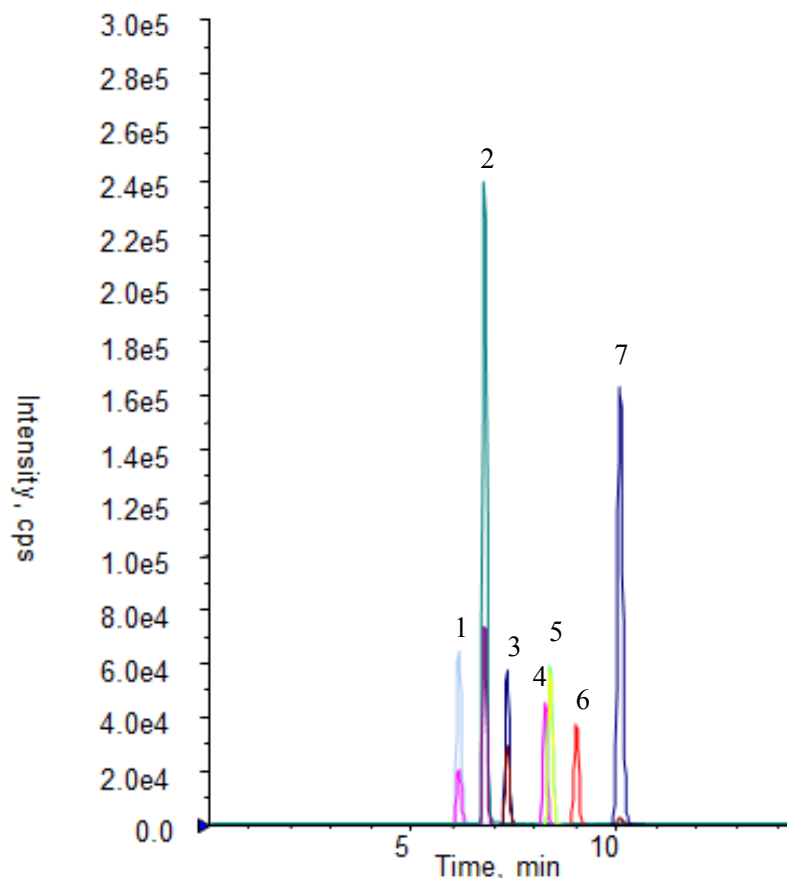


1. Dinotefuran 2. Nitenpyram 3. Thiamethoxam

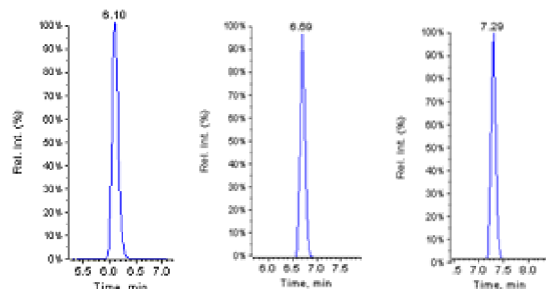

4. Imidacloprid 5. Clothianidin 6. Acetamiprid

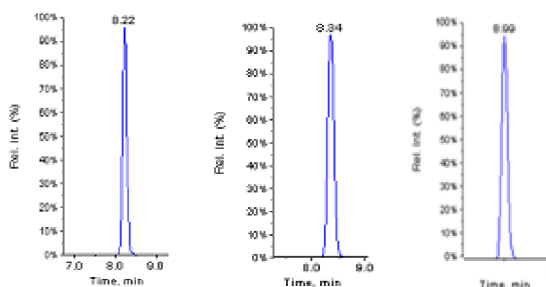

7. Thiacloprid

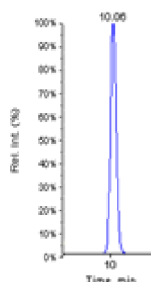

Conditions

Column : InertSustain C18 (3 μm, 75 × 2.1 mm I.D.)
Column Cat. No. : 5020-07413
Eluent : A) 5 mM CH₃COONH₄ in CH₃OH
 B) 5 mM CH₃COONH₄ in H₂O
 A/B = 5/95 – 4 min – 40/60 – 10 min – 50/50

Flow Rate : 0.15 mL/min.


Col. Temp. : 40 °C

Detection : LC/MS/MS (4000QTRAP®: ESI Positive, MRM)

Injection Vol. : 100 μL

No	Compound	Mw	m/z		
			Q1	Q3	
1	Dinotefuran	202.2	203	129	87
2	Nitenpyram	270.7	271	126	237
3	Thiamethoxam	291.7	292	211	181
4	Imidacloprid	255.7	256	209	175
5	Clothianidin	249.7	250	169	132
6	Acetamiprid	222.7	223	126	90
7	Thiacloprid	252.7	253	126	90

Conc. 1.0 μg/L

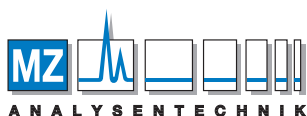

AUTHORIZED DISTRIBUTOR
 MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz
 Tel +49 6131 880 96-0, Fax +49 6131 880 96-20
 e-mail: info@mz-at.de, www.mz-at.de