


# InertSearch™ for LC


Inertsil® Applications

## Organic acids

Data No. LA126-0000

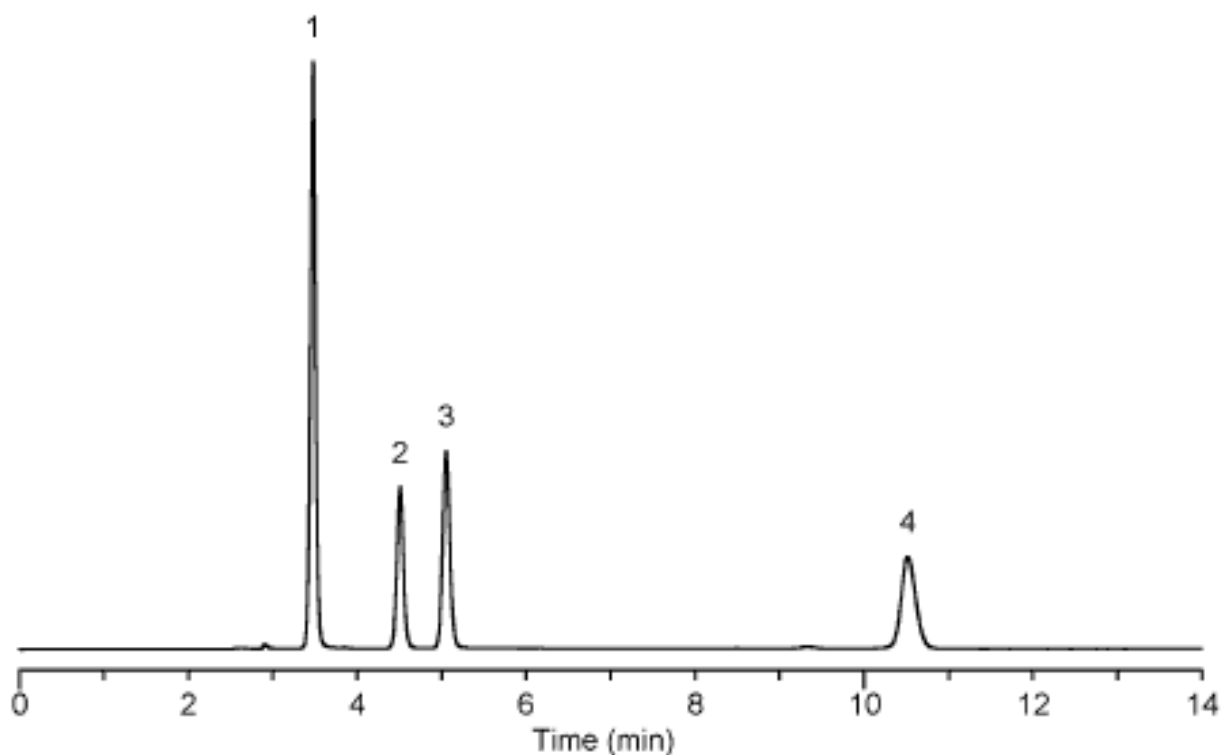

Column	: Inertsil C8-3 5µm 250×4.6mmI.D.	Solutes	: 1. Formic acid (0.1 %) 2. Lactic acid (0.1 %) 3. Acetic acid (0.1 %) 4. Propionic acid (0.1 %)
Column Cat. No.	: 5020-01901		
Eluent	: 20mM (Potassium) Phosphate buffer (pH3.0)		
Flow rate	: 1.0mL/min		
Col. Temp.	: 40°C		
Detection	: UV 210 nm		
Injection Vol.	: 20µL		
Sample	: Organic acids		

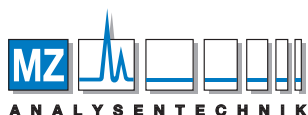

**AUTHORIZED DISTRIBUTOR**

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz  
Tel +49 6131 880 96-0, Fax +49 6131 880 96-20  
e-mail: info@mz-at.de, www.mz-at.de

**GL Sciences**