

DACAPO DX-C18

150 x 3 mm

Application



LC-MS: Nucleosides and bases (mobile phase pH)

核酸塩基およびヌクレオシドのLC-MS分析 (pHの影響)

10mM NH₄OH / MeOH = 100 / 0.5

10mM HCOOH / MeOH = 100 / 0.5

AUTHORIZED DISTRIBUTOR
 MZ-Analysetechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz
 Tel. +49 6131 880 96-0, Fax +49 6131 880 96-20
 e-mail: info@mz-at.de, www.mz-at.de
 MZ ANALYSETECHNIK

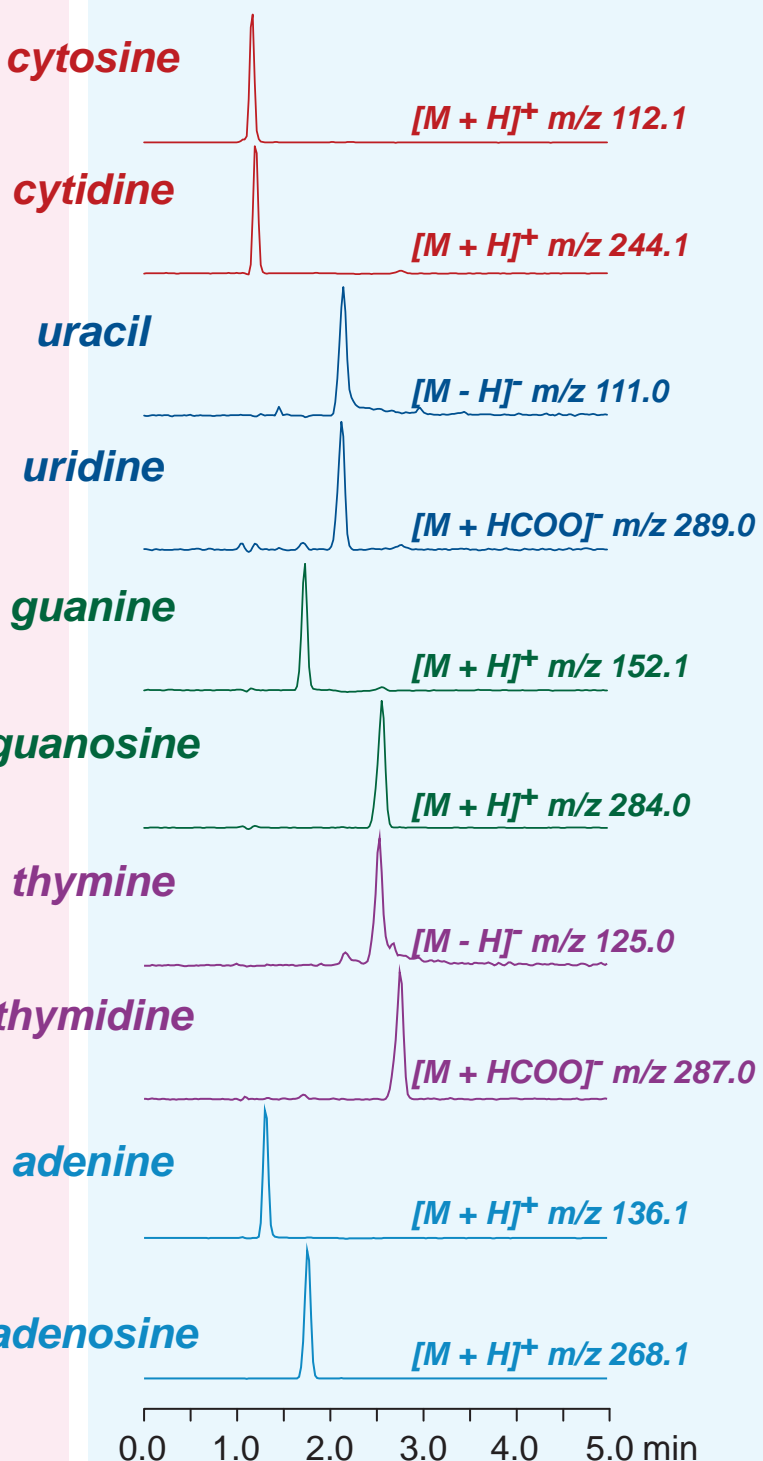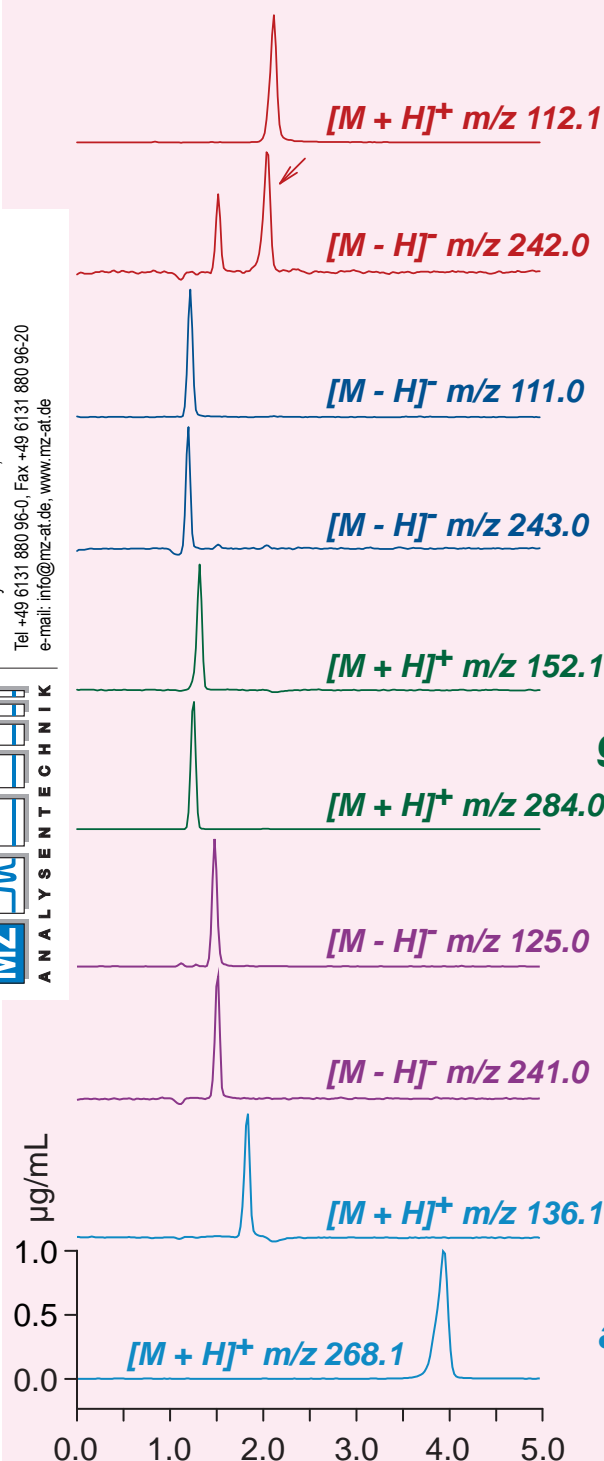

DACAPO DX-C18, 150 x 3 mm
 0.4 mL/min (13 MPa), 37 °C, 10 μL (1 ug/mL, eluent), Single Quad. MS (ESI)