

Intrada Amino Acid

50 x 3 mm

Application

55 Amino acids (standard)

55種アミノ酸(標準試料)

Intrada Amino Acid, 50 x 3 mm

**A: acetonitrile /tetrahydrofuran /25mM ammonium formate /formic acid
= 9 / 75 / 16 / 0.3**

B: acetonitrile / 100mM ammonium formate = 20 / 80



0 %B (0-2.5min)

0-17 %B (2.5-6.5 min)

100 %B (6.5-10 min)

0.6 mL/min (6MPa), 35 deg.C, 5 uL (1000nmol/mL)

Shimadzu LCMS-2020, ESI (SIM, positive)

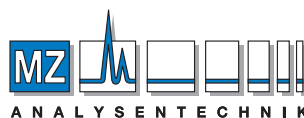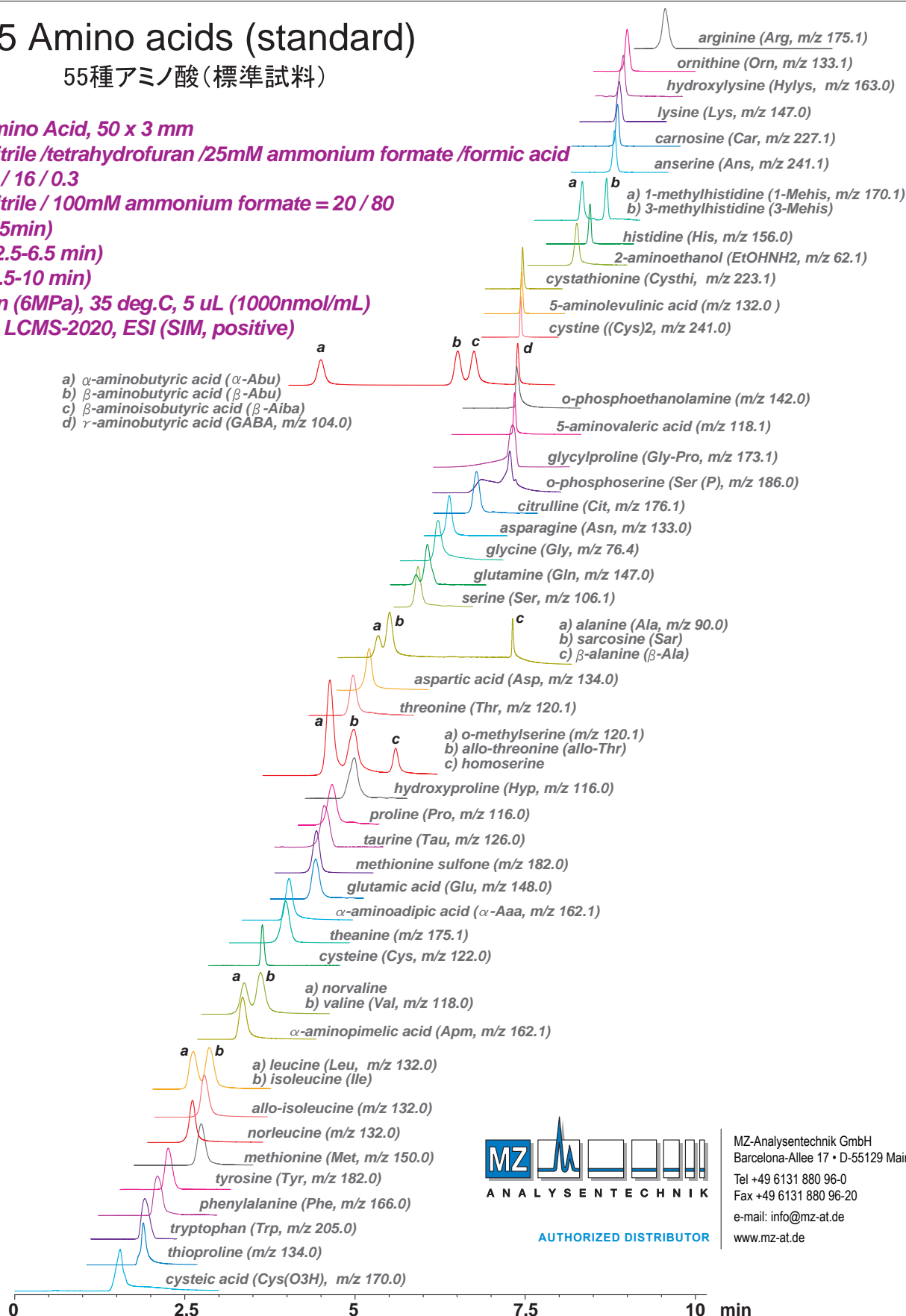

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de