

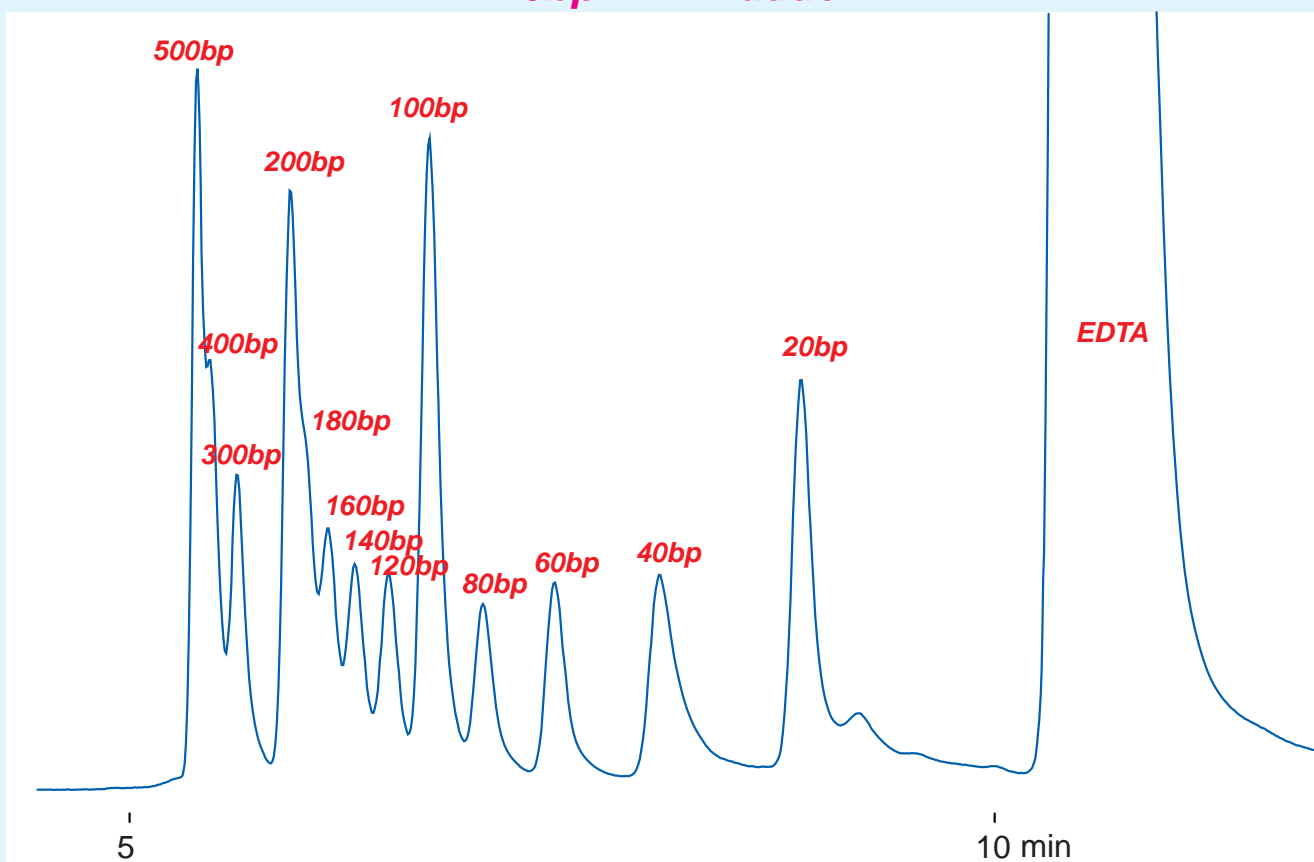
Intrada SEC

250 x 4.6 mm

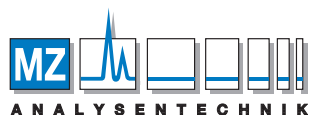
Application

SEC: 20bp DNA Ladder

20bp DNA ラダーのSEC分離

**SEC: DNA Fragments****Intrada SEC****20bp DNA Ladder**

*Intrada SEC, 250 x 4.6 mm
100mM ammonium formate
0.3 mL/min, 37°C, 260nm
20bp DNA Ladder (200µg/mL), 5 µL*

**AUTHORIZED DISTRIBUTOR**

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0, Fax +49 6131 880 96-20
e-mail: info@mz-at.de, www.mz-at.de