

Presto FF-C18


150 x 4.6 mm

Application

LC-MS: Deoxyoligonucleotides

デオキシオリゴヌクレオチドのLC-MS分析

5'-GTTTTCCCAGTCACGACGTT-3'

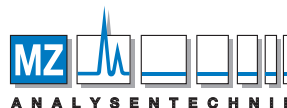

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
 Barcelona-Allee 17 • D-55129 Mainz
 Tel +49 6131 880 96-0
 Fax +49 6131 880 96-20
 e-mail: info@mz-at.de
 www.mz-at.de

$[M-7H]^{7-}$ m/z 864.4

5'-GGAAACAGCTATGACCATGA-3'

$[M-7H]^{8-}$ m/z 768.5

0 1 2 3 4 5 6 7 8 min

Presto FF-C18, 150 x 4.6 mm
 A: 10mM ammonium formate
 B: acetonitrile / 20mM ammonium formate = 50 / 50
 0-100%B (0-8min)
 0.4mL/min (13 MPa), 60°C, MS (ESI, SIM)
 100µL (0.8nmol/mL, 50mM HCOONH₄)