


Presto FF-C18

150 x 4.6 mm

Application

**DNA from salmon sperm**

DNA (サケ精巢由来), MW: 20,000,000 ~ 30,000,000

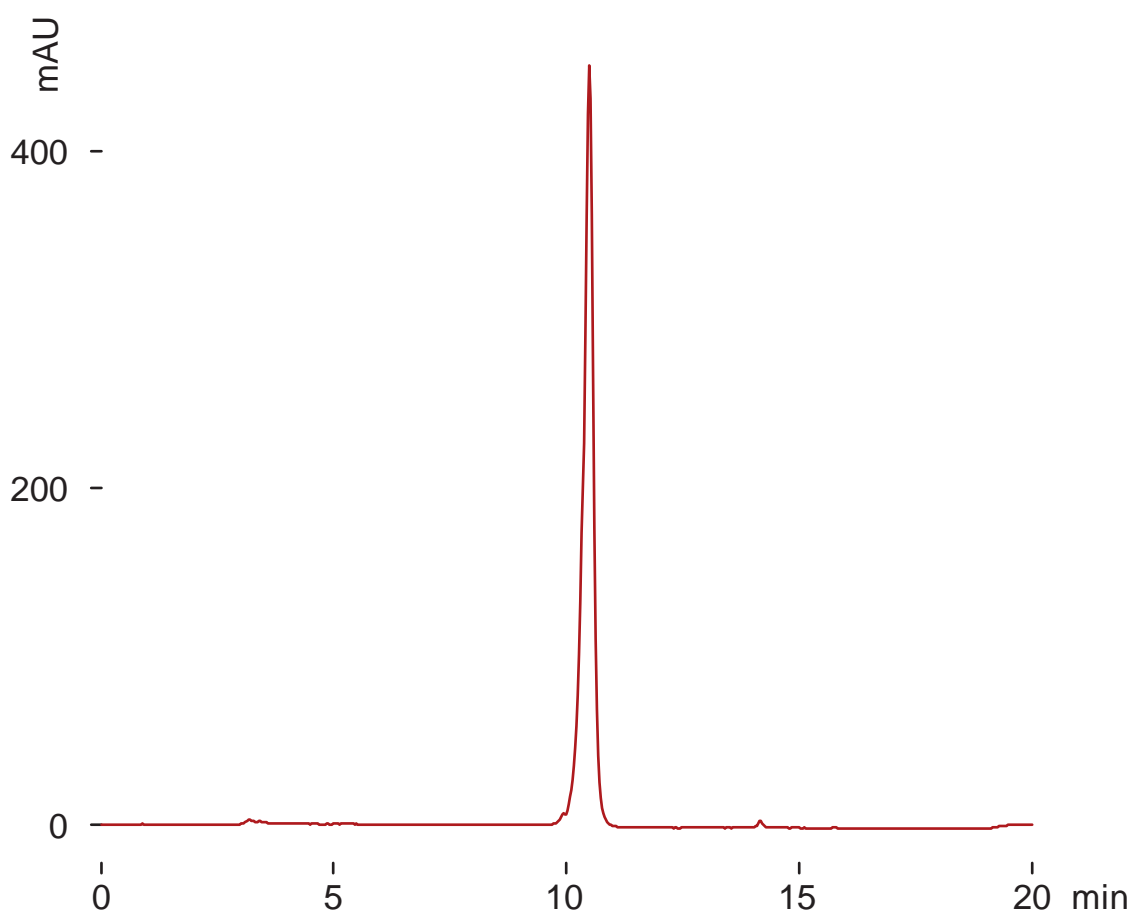

Presto FF-C18, 150 x 4.6 mm

A: 50mM ammonium acetate, B: acetonitrile, 0-50 %B (0-15min)

0.3 mL/min (12MPa), 37 deg.C, 260 nm

1 uL (4ug, 50mM ammonium acetate)

**AUTHORIZED DISTRIBUTOR**  
MZ-Analyse-technik GmbH, Barcelona-Allee 17 • D-55129 Mainz  
Tel +49 6131 880 96-0, Fax +49 6131 880 96-20  
e-mail: info@mz-at.de, www.mz-at.de

