


Presto FF-C18

250 x 4.6 mm

Application

Double-stranded DNA fragments

二重鎖DNAフラグメント

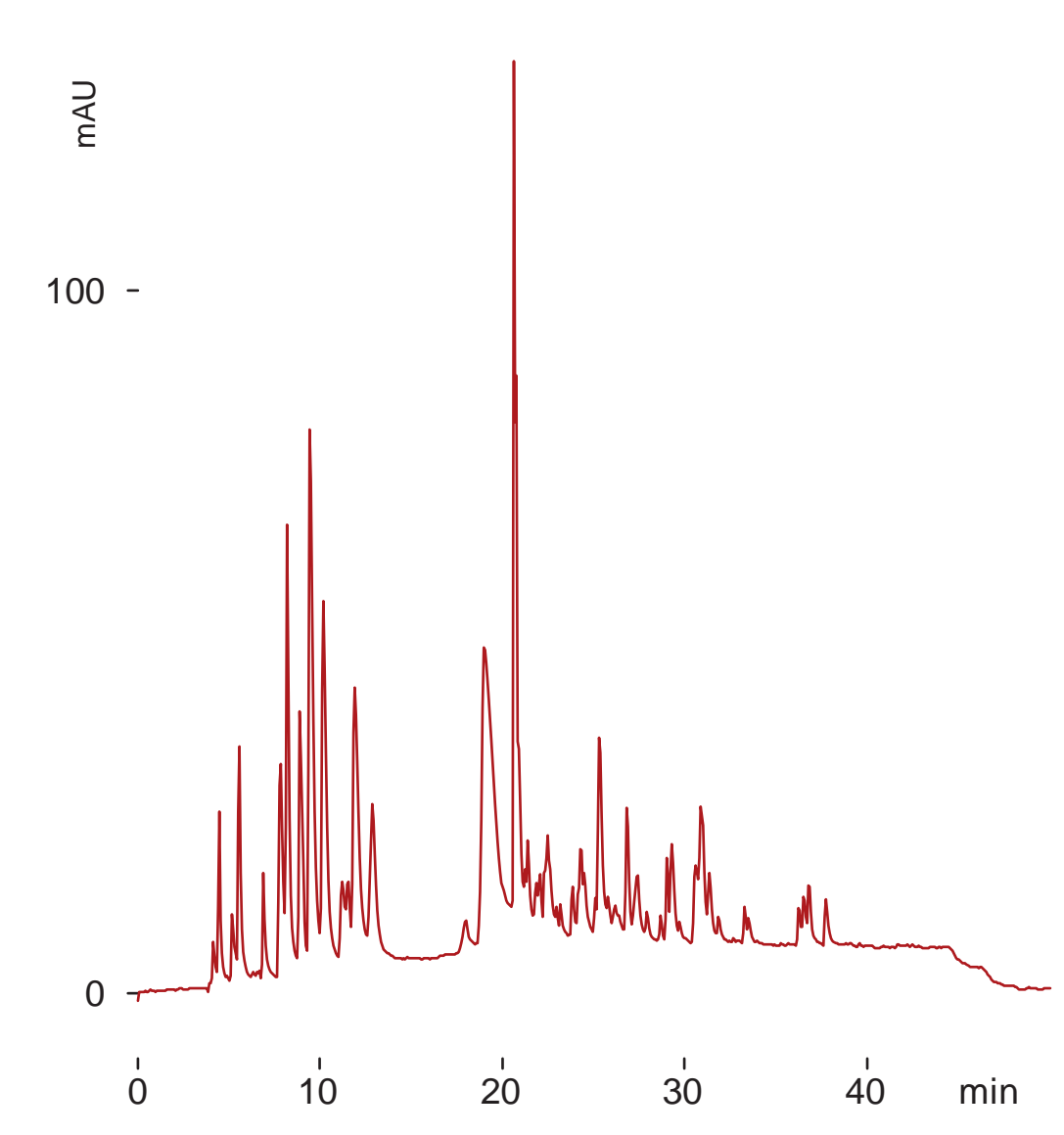

Presto FF-C18, 250 x 4.6 mm

A: water / triethylamine / acetic acid = 100 / 1.1 / 0.5

B: eluent A / acetonitrile = 90 / 10

0-40 %B (0-40min), 0.8 mL/min (24MPa), 50 deg.C, 260 nm

2 uL (salmon DNA / DNase)

AUTHORIZED DISTRIBUTOR
 MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz
 Tel +49 6131 880 96-0, Fax +49 6131 880 96-20
 e-mail: info@mz-at.de, www.mz-at.de

