


Presto FF-C18

150 x 4.6 mm


**HPLC COLUMN APPLICATION**

## Oligonucleotides ( Oligothymidylic acid )

オリゴヌクレオチド ( オリゴチミジル酸 )

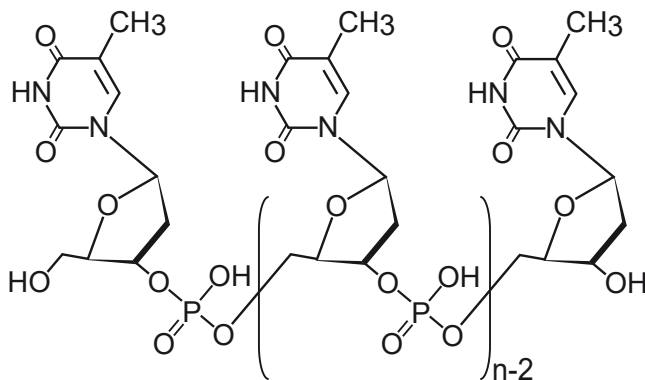


*oligothymidylic acid*  
*d(pT)<sub>12-18</sub>*

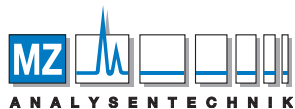

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH  
Barcelona-Allee 17 • D-55129 Mainz  
Tel +49 6131 880 96-0  
Fax +49 6131 880 96-20  
e-mail: info@mz-at.de  
www.mz-at.de

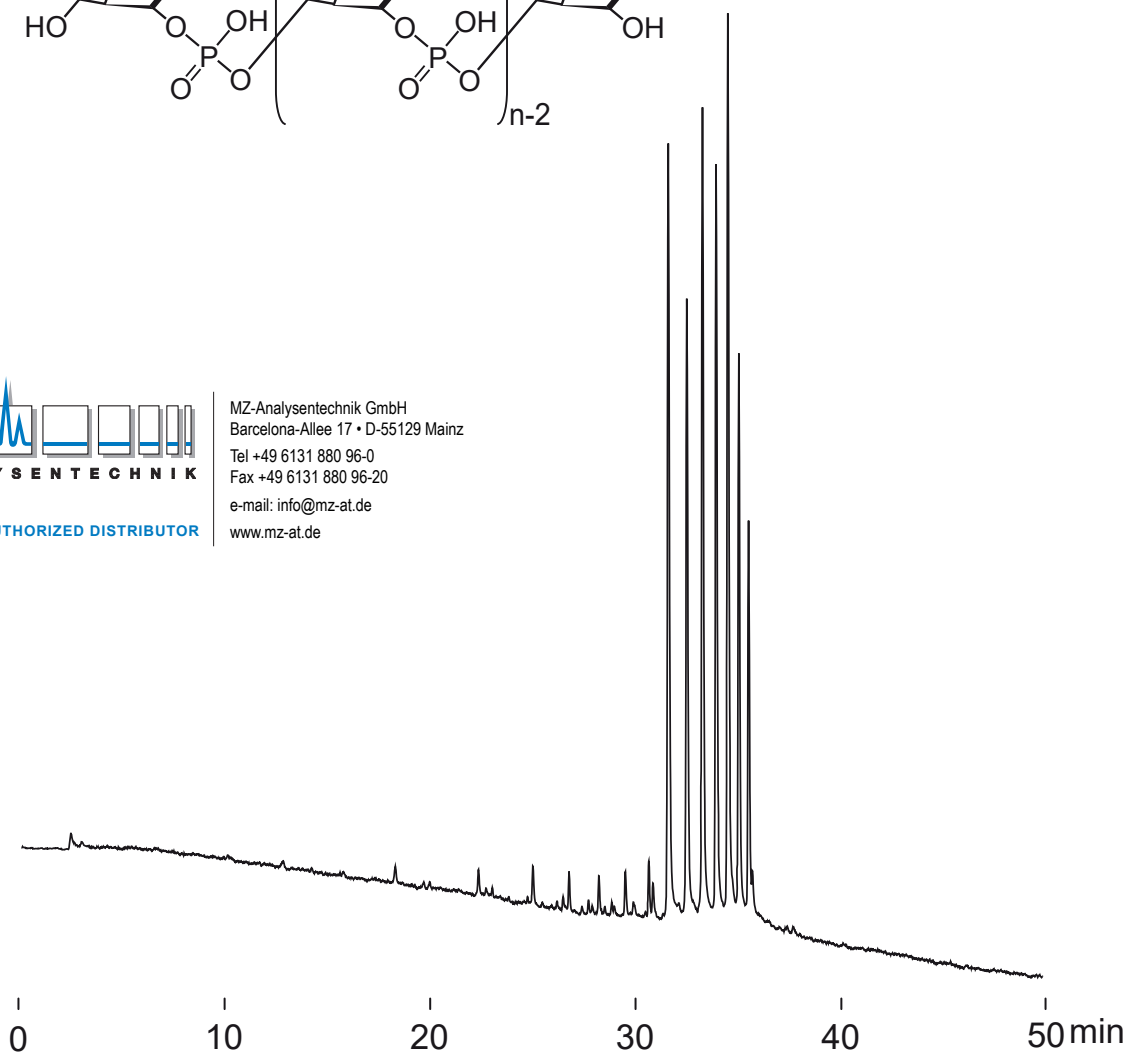

150 x 4.6 mm

A: 100mM ammonium formate

B: methanol / 100mM ammonium formate = 90 / 10

5-35%B (0-50min)

0.4mL/min (15MPa), 60°C, 260nm

4µL (20µg/mL, 100mM ammonium formate)