


Presto FF-C18

250 x 4.6 mm


Application

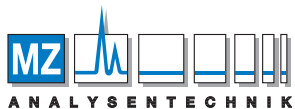

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de

Plasmid DNA プラスミド DNA

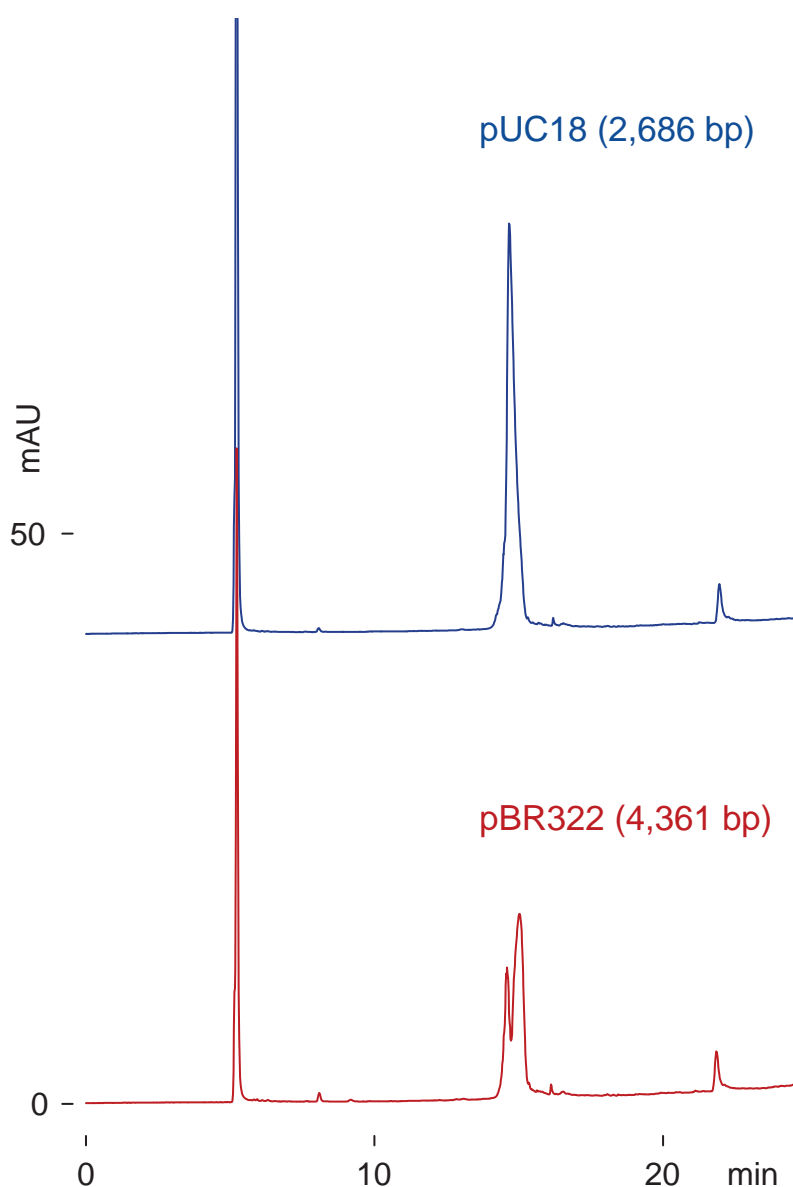

Presto FF-C18, 250 x 4.6 mm
A: 200mM ammonium acetate
B: acetonitrile
0-30 %B (0-20min)
0.3 mL/min (17MPa), 50 deg.C, 260 nm, 1 uL (4ug)