







Scherzo SM-C18

250 x 3 mm

Application

Organophosphorous pesticides,
quaternary ammonium herbicides

有機リン系農薬, 四級アンモニウム系農薬

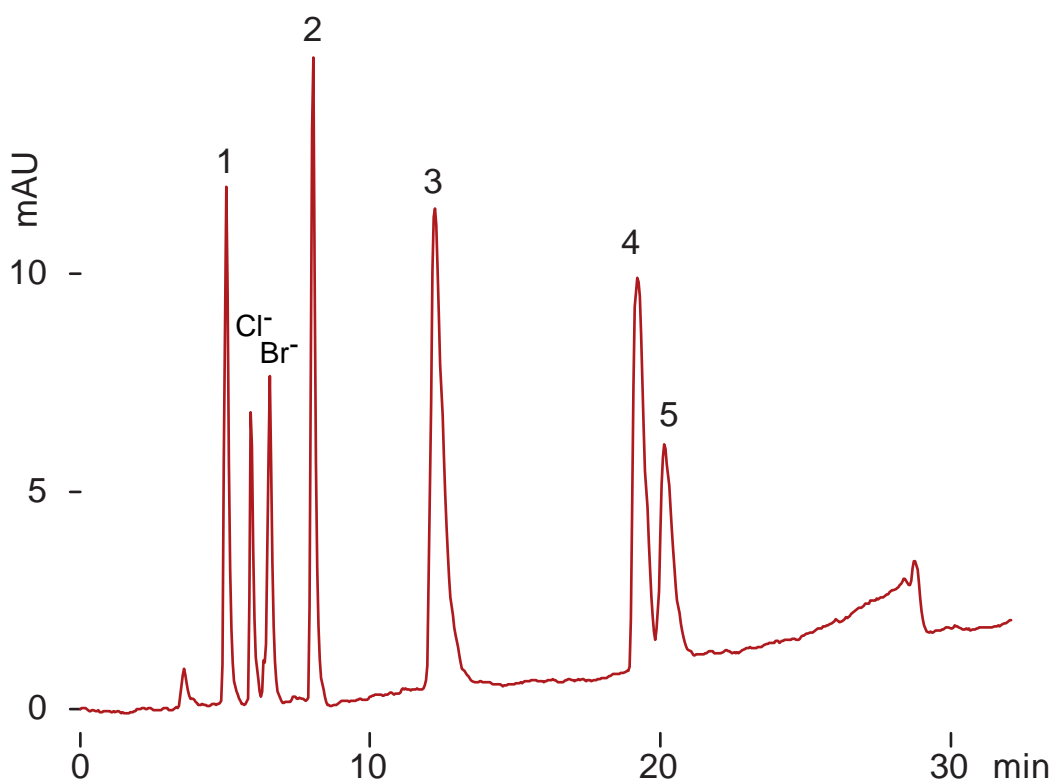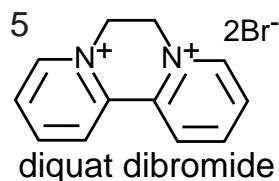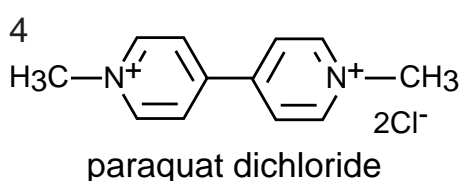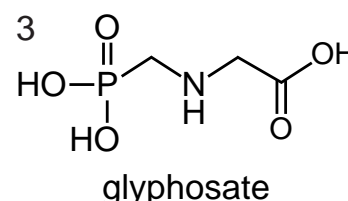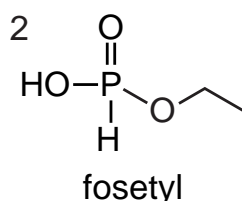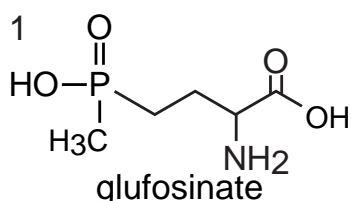

Scherzo SM-C18, 250 x 3 mm

A: 5 mM ammonium formate

B: 100 mM ammonium formate / acetonitrile = 25 / 75

0-100 %B (0-25min), 0.4 mL/min (12MPa)

37 deg.C, ELSD (spray chamber 10 deg.C, drift tube 50 deg.C)

5 uL (0.5-2ug)