


Scherzo SW-C18

150 x 3 mm


Application

siRNA

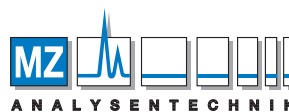

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
 Barcelona-Allee 17 • D-55129 Mainz
 Tel +49 6131 880 96-0
 Fax +49 6131 880 96-20
 e-mail: info@mz-at.de
 www.mz-at.de

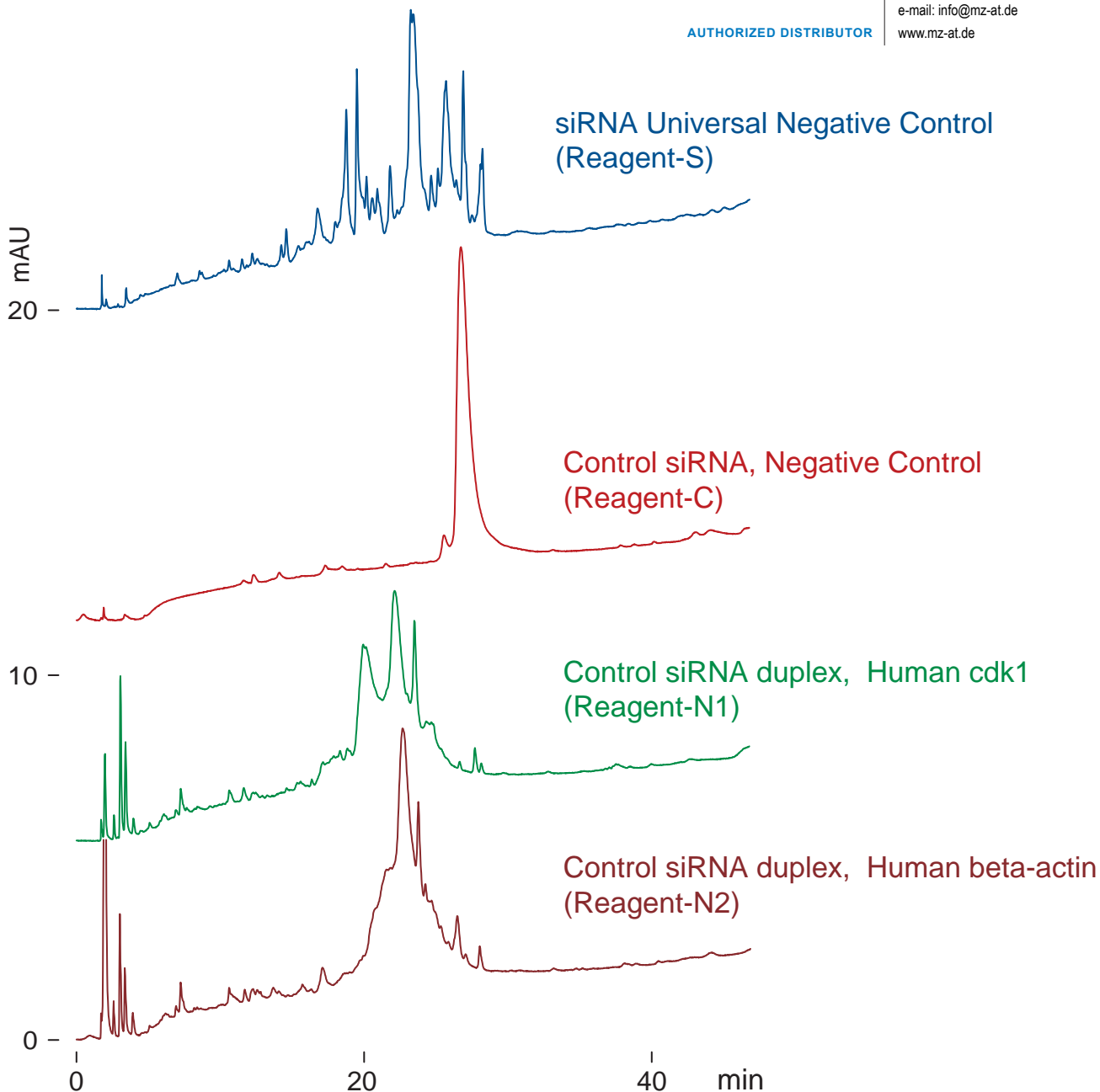

Scherzo SW-C18, 150 x 3 mm
 A: 10 mM ammonium acetate
 B: 200 mM ammonium acetate / acetonitrile = 85 / 15
 0-100 %B (0-45 min)
 0.4 mL/min (10 MPa), 37 deg.C, 260 nm, 0.5-1 uL (35 - 70 ng)