


Unison UK-Amino

150 x 2 mm


Application

Dehydroascorbic Acid and Related Compounds


デヒドロアスコルビン酸と関連化合物

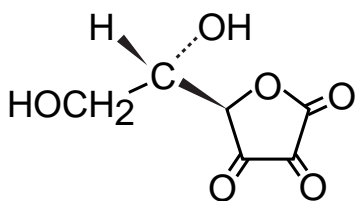


1 dehydroascorbic acid

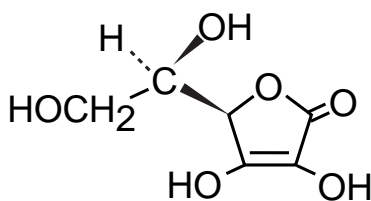

2 isoascorbic acid (erythorbic acid)

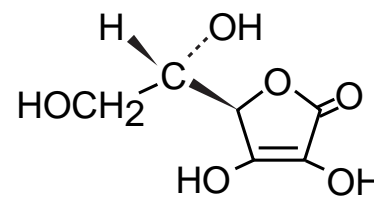


3 L-ascorbic acid (Vitamin C)

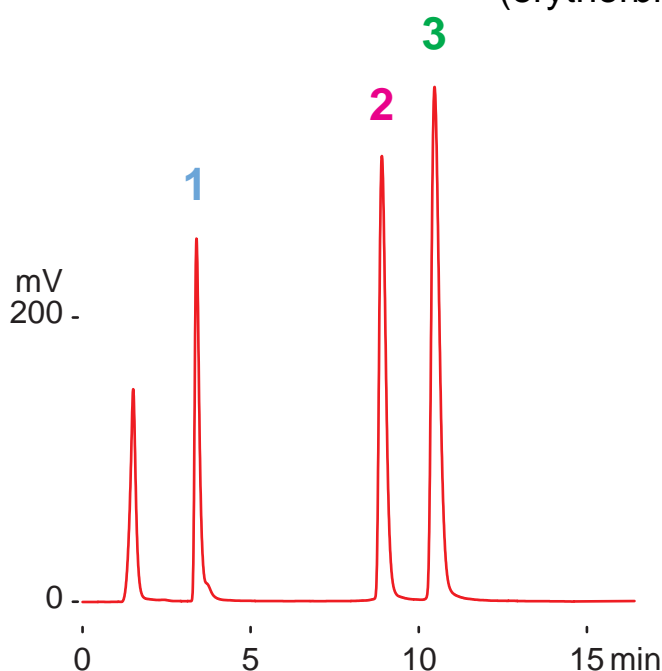


normal phase

Unison UK-Amino
150 x 2 mm
acetonitrile /water /acetic acid =
80 / 20 / 5
0.2mL/min (5MPa),
37deg.C, ELSD,
1uL(2-4ug)

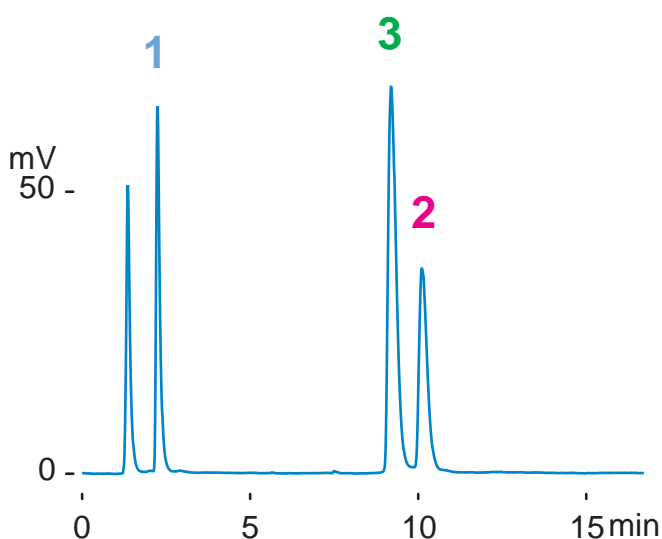

anion exchange

Unison UK-Amino
150 x 2 mm
water /acetic acid = 100 /0.5
0.2mL/min (9MPa),
37deg.C, ELSD,
1uL(2-4ug)

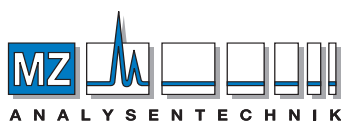

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17• D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de