

Unison UK-C18


50 x 3 mm

Application

LC-MS: Deoxyoligonucleotides

デオキシオリゴヌクレオチドのLC-MS分析


5'-GTTTTCCAGTCACGACGTT-3'

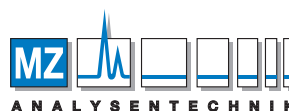

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de

5'-GGAAACAGCTATGACCATGA-3'

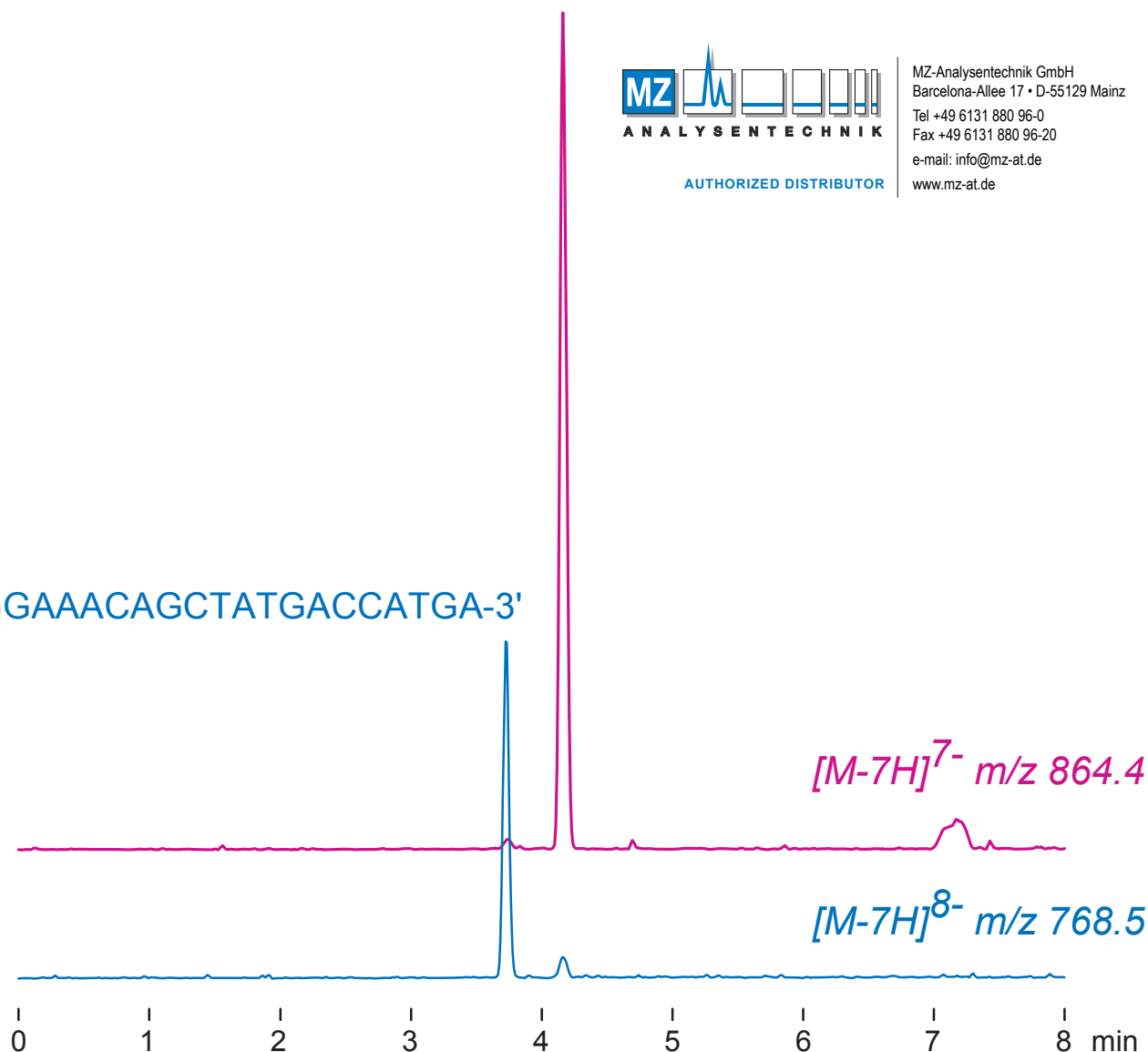

Unison UK-C18, 50 x 3 mm

A: 10mM ammonium formate

B: acetonitrile / 20mM ammonium formate = 50 / 50

0-40%B (0-8min)

0.4mL/min (3 MPa), 60°C, MS (ESI, SIM)

100µL (0.8nmol/mL, 50mM HCOONH₄)