
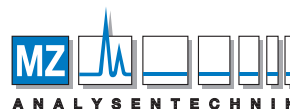


Unison US-C18 150 x 4.6 mm 150 x 3 mm 150 x 2 mm Application

Oligonucleotides


オリゴヌクレオチド

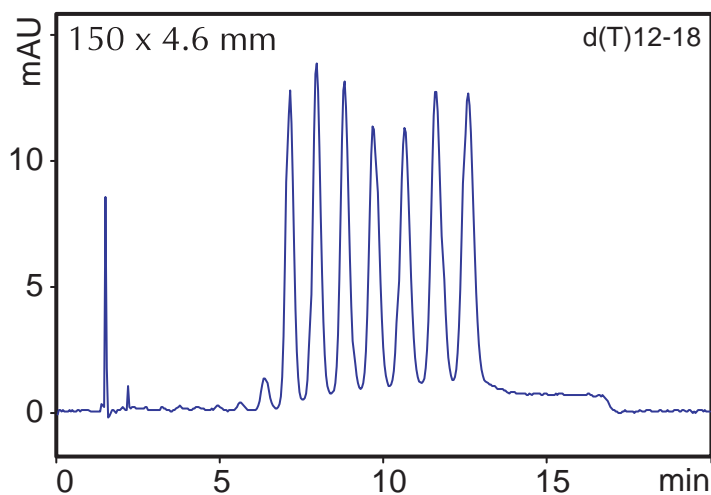


AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
 Barcelona-Allee 17 • D-55129 Mainz
 Tel +49 6131 880 96-0
 Fax +49 6131 880 96-20
 e-mail: info@mz-at.de
 www.mz-at.de

Unison US-C18, 150 x 4.6 mm
 A: 10mM dibutylammonium acetate
 B: acetonitrile
 28-32%B (0-10min)
 1.2 mL/min, 37 °C, 260 nm, 7.8 MPa

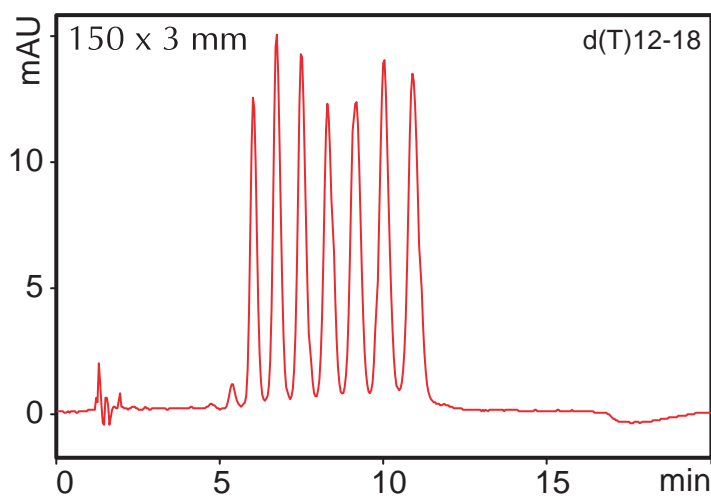

Unison US-C18, 150 x 3 mm
 A: 10mM dibutylammonium acetate
 B: acetonitrile
 28-32%B (0-10min)
 0.6 mL/min, 37 °C, 260 nm, 7.2 MPa

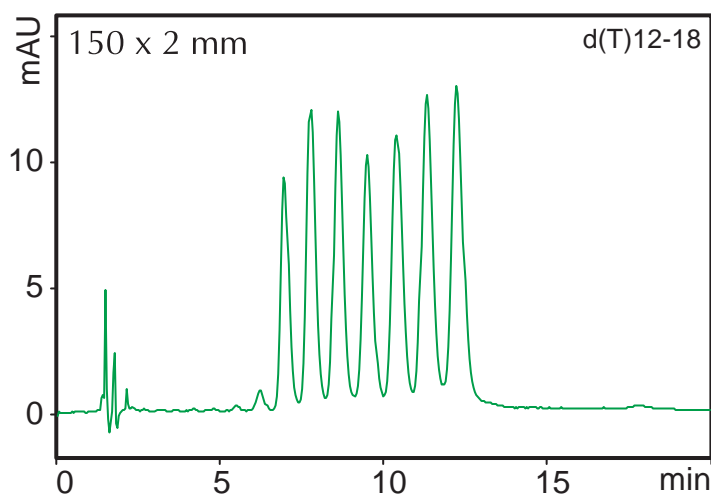

Unison US-C18, 150 x 2 mm
 A: 10mM dibutylammonium acetate
 B: acetonitrile
 28-32%B (0-10min)
 0.24 mL/min, 37 °C, 260 nm, 7.3 MPa