


Separation of Acrylamide, Methacrylamide and Methacrylic Acid

Application Number 123010

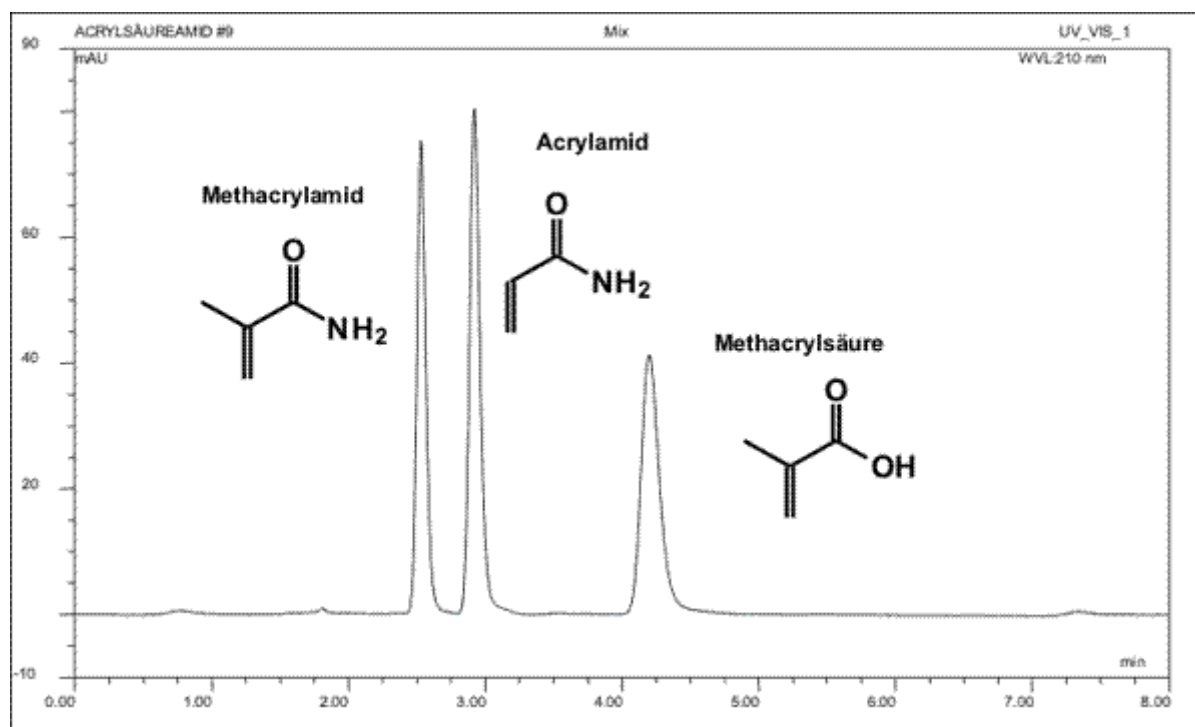

Legend	1 - methacrylamide
	2 - acrylamide
	3 - methacrylic acid
Column	EC 125/4 NUCLEODUR HILIC, 5µm
REF	760551.40
Conditions	Eluent: acetonitril / 0.1 % formic acid (98:2, v/v)
Flow rate	0.6 mL/min
Temperature	22 °C
Injection volume	0.5 µL
Detection	UV, 210 nm
Author	HPLC department
Source	MACHEREY-NAGEL, Germany 2009