
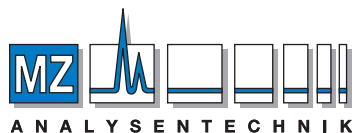


HPLC Analysis of Fipronil Enantiomers on Astec® (R,R) P-CAP™

column Astec (R,R) P-CAP, 25 cm x 4.6 mm I.D., 5 µm particles (31024AST)
 mobile phase (A) heptane; (B) ethanol; (80:20, A:B)
 flow rate 1 mL/min
 column temp. 25 °C
 detector UV, 235 nm
 injection 1 µL
 sample 1 mg/mL in methanol
 Application No. **G004720**

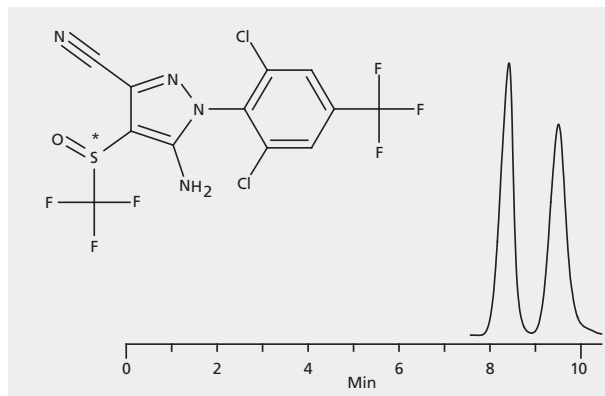

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17• D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de