

# <621> CHROMATOGRAPHY

## INTRODUCTION

Chromatographic separation techniques are multistage separation methods in which the components of a sample are distributed between two phases, of which one is stationary and the other is mobile. The stationary phase may be a solid or a liquid supported on a solid or a gel. The stationary phase may be packed in a column, spread as a layer, distributed as a film, or applied by other techniques. The mobile phase may be in a gaseous or liquid form, or a supercritical fluid. The separation may be based on adsorption, mass distribution (partition), or ion exchange; or it may be based on differences among the physicochemical properties of the molecules, such as size, mass, and volume. This chapter contains general procedures, definitions, and calculations of common parameters and describes general requirements for system suitability. The types of chromatography useful in qualitative and quantitative analyses employed in *USP* procedures are column, gas (GC), paper, thin-layer (TLC) [including high-performance thin-layer chromatography (HPTLC)], and pressurized liquid chromatography [commonly called high-pressure or high-performance liquid chromatography (HPLC)].

## GENERAL PROCEDURES

This section describes the basic procedures used when a chromatographic method is described in a monograph. These procedures are followed unless otherwise indicated in the individual monograph.

### Paper Chromatography

#### STATIONARY PHASE

The stationary phase is a sheet of paper of suitable texture and thickness. Development may be ascending, in which the solvent is carried up the paper by capillary forces, or descending, in which the solvent flow is also assisted by gravitational force. The orientation of paper grain, with respect to solvent flow, is to be kept constant in a series of chromatograms. The machine direction is usually designated by the manufacturer.

#### APPARATUS

The essential equipment for paper chromatography consists of a vapor-tight chamber with inlets for the addition of solvent and a rack of corrosion-resistant material about 5 cm shorter than the inside height of the chamber. The rack serves as a support for solvent troughs and antisiphon rods that, in turn, hold up the chromatographic sheets. The bottom of the chamber is covered with the prescribed solvent system or mobile phase. Saturation of the chamber with solvent vapor is facilitated by lining the inside walls with paper wetted with the prescribed solvent system.

#### SPOTTING

The substance or substances analyzed are dissolved in a suitable solvent. Convenient volumes delivered from suitable micropipettes of the resulting solution, normally containing 1–20 µg of the compound, are placed in 6- to 10-mm spots NLT 3 cm apart.

#### DESCENDING PAPER CHROMATOGRAPHY PROCEDURE

1. A spotted chromatographic sheet is suspended in the apparatus, using the antisiphon rod to hold the upper end of the sheet in the solvent trough. [NOTE—Ensure that the portion of the sheet hanging below the rods is freely suspended in the chamber without touching the rack, the chamber walls, or the fluid in the chamber.]
2. The chamber is sealed to allow equilibration (saturation) of the chamber and the paper with the solvent vapor. Any excess pressure is released as necessary.
3. After equilibration of the chamber, the prepared mobile phase is introduced into the trough through the inlet.
4. The inlet is closed, and the mobile solvent phase is allowed to travel the desired distance down the paper.
5. The sheet is removed from the chamber.
6. The location of the solvent front is quickly marked, and the sheet is dried.
7. The chromatogram is observed and measured directly or after suitable development to reveal the location of the spots of the isolated drug or drugs.

## ASCENDING PAPER CHROMATOGRAPHY PROCEDURE

1. The mobile phase is added to the bottom of the chamber.
2. The chamber is sealed to allow equilibration (saturation) of the chamber and the paper with the solvent vapor. Any excess pressure is released as necessary.
3. The lower edge of the stationary phase is dipped into the mobile phase to permit the mobile phase to rise on the chromatographic sheet by capillary action.
4. When the solvent front has reached the desired height, the chamber is opened, the sheet is removed, the location of the solvent front is quickly marked, and the sheet is dried.
5. The chromatogram is observed and measured directly or after suitable development to reveal the location of the spots of the isolated drug or drugs.

**Thin-Layer Chromatography**

## STATIONARY PHASE

The stationary phase is a relatively thin, uniform layer of dry, finely powdered material applied to a glass, plastic, or metal sheet or plate (typically called the plate). The stationary phase of TLC plates has an average particle size of 10–15  $\mu\text{m}$ , and that of HPTLC plates has an average particle size of 5  $\mu\text{m}$ . Commercial plates with a preadsorbent zone can be used if they are specified in a monograph. The sample applied to the preadsorbent region develops into sharp, narrow bands at the preadsorbent–sorbent interface. The separations achieved may be based on adsorption, partition, or a combination of both effects, depending on the particular type of stationary phase.

## APPARATUS

A chromatographic chamber made of inert, transparent material and having the following specifications is used: a flat-bottom or twin trough, a tightly fitted lid, and a size suitable for the plates. The chamber is lined on at least one wall with filter paper. Sufficient mobile phase or developing solvent is added to the chamber so that, after impregnation of the filter paper, a depth appropriate to the dimensions of the plate used is available. The chromatographic chamber is closed and allowed to equilibrate. [NOTE—Unless otherwise indicated, the chromatographic separations are performed in a saturated chamber.]

## DETECTION/VISUALIZATION

An ultraviolet (UV) light source suitable for observations under short- (254 nm) and long- (365 nm) wavelength UV light and a variety of other spray reagents, used to make spots visible, is often used.

## SPOTTING

Solutions are spotted on the surface of the stationary phase (plate) at the prescribed volume in sufficiently small portions to obtain circular spots of 2–5 mm in diameter (1–2 mm on HPTLC plates) or bands of 10–20 mm  $\times$  1–2 mm (5–10 mm  $\times$  0.5–1 mm on HPTLC plates) at an appropriate distance from the lower edge and sides of the plate. [NOTE—During development, the application position must be at least 5 mm (TLC) or 3 mm (HPTLC) above the level of the mobile phase.] The solutions are applied on a line parallel to the lower edge of the plate with an interval of at least 10 mm (5 mm on HPTLC plates) between the centers of spots or 4 mm (2 mm on HPTLC plates) between the edges of bands, then allowed to dry.

## PROCEDURE

1. Place the plate in the chamber, ensuring that the spots or bands are above the surface of the mobile phase.
2. Close the chamber.
3. Allow the mobile phase to ascend the plate until the solvent front has traveled three-quarters of the length of the plate, or the distance prescribed in the monograph.
4. Remove the plate, mark the solvent front with a pencil, and allow to dry.
5. Visualize the chromatograms as prescribed.
6. Determine the chromatographic *Retardation factor* ( $R_f$ ) values for the principal spots or zones.
7. Presumptive identification can be made by observation of spots or zones of identical  $R_f$  value and about equal magnitude obtained, respectively, with an unknown and a standard chromatographed on the same plate. A visual comparison of the size or intensity of the spots or zones may serve for semiquantitative estimation. Quantitative measurements are possible by means of densitometry (absorbance or fluorescence measurements).

## Column Chromatography

### SOLID SUPPORT

Purified siliceous earth is used for normal-phase separation. Silanized chromatographic siliceous earth is used for reverse-phase partition chromatography.

### STATIONARY PHASE

The solid support is modified by the addition of a stationary phase specified in the individual monograph. If a mixture of liquids is used as the stationary phase, mix the liquids before the introduction of the solid support.

### MOBILE PHASE

The mobile phase is specified in the individual monograph. If the stationary phase is an aqueous solution, equilibrate with water. If the stationary phase is a polar organic fluid, equilibrate with that fluid.

### APPARATUS

Unless otherwise specified in the individual monograph, the chromatographic tube is about 22 mm in its inside diameter and 200–300 mm long. Attached to it is a delivery tube, without stopcock, about 4 mm in its inside diameter and about 50 mm long.

**Apparatus preparation:** Pack a pledget of fine glass wool in the base of the tube. Combine the specified volume of stationary phase and the specified amount of solid support to produce a homogeneous, fluffy mixture. Transfer this mixture to the chromatographic tube, and tamp using gentle pressure to obtain a uniform mass. If the specified amount of solid support is >3 g, transfer the mixture to the column in portions of approximately 2 g, and tamp each portion. If the assay or test requires a multisegment column with a different stationary phase specified for each segment, tamp after the addition of each segment, and add each succeeding segment directly to the previous one. Pack a pledget of fine glass wool above the completed column packing. [NOTE—The mobile phase should flow through a properly packed column as a moderate stream or, if reverse-phase chromatography is applied, as a slow trickle.]

If a solution of the analyte is incorporated into the stationary phase, complete the quantitative transfer to the chromatographic tube by scrubbing the beaker used for the preparation of the test mixture with a mixture of about 1 g of solid support and several drops of the solvent used to prepare the sample solution before adding the final portion of glass wool.

### PROCEDURE

1. Transfer the mobile phase to the column space above the column packing, and allow it to flow through the column under the influence of gravity.
2. Rinse the tip of the chromatographic column with about 1 mL of mobile phase before each change in composition of mobile phase and after completion of the elution.
3. If the analyte is introduced into the column as a solution in the mobile phase, allow it to pass completely into the column packing, then add the mobile phase in several small portions, allowing each to drain completely, before adding the bulk of the mobile phase.
4. Where the procedure indicates the use of multiple chromatographic columns mounted in series and the addition of mobile phase in divided portions is specified, allow each portion to drain completely through each column, and rinse the tip of each with mobile phase before the addition of each succeeding portion.

## Gas Chromatography

### LIQUID STATIONARY PHASE

This type of phase is available in packed or capillary columns.

### PACKED COLUMN GAS CHROMATOGRAPHY

The liquid stationary phase is deposited on a finely divided, inert solid support, such as diatomaceous earth, porous polymer, or graphitized carbon, which is packed into a column that is typically 2–4 mm in internal diameter and 1–3 m in length.

## CAPILLARY COLUMN GAS CHROMATOGRAPHY

In capillary columns, which contain no packed solid support, the liquid stationary phase is deposited on the inner surface of the column and may be chemically bonded to it.

## SOLID STATIONARY PHASE

This type of phase is available only in packed columns. In these columns the solid phase is an active adsorbent, such as alumina, silica, or carbon, packed into a column. Polyaromatic porous resins, which are sometimes used in packed columns, are not coated with a liquid phase. [NOTE—Packed and capillary columns must be conditioned before use until the baseline and other characteristics are stable. The column or packing material supplier provides instructions for the recommended conditioning procedure.]

## APPARATUS

A gas chromatograph consists of a carrier gas source, injection port, column, detector, and recording device. The injection port, column, and detector are temperature controlled and may be varied as part of the analysis. The typical carrier gas is helium, nitrogen, or hydrogen, depending on the column and detector in use. The type of detector used depends on the nature of the compounds analyzed and is specified in the individual monograph. Detector output is recorded as a function of time, and the instrument response, measured as peak area or peak height, is a function of the amount present.

## TEMPERATURE PROGRAM

The length and quality of a GC separation can be controlled by altering the temperature of the chromatographic column. When a temperature program is necessary, the individual monograph indicates the conditions in table format. The table indicates the initial temperature, rate of temperature change (ramp), final temperature, and hold time at the final temperature.

## PROCEDURE

1. Equilibrate the column, injector, and detector with flowing carrier gas until a constant signal is received.
2. Inject a sample through the injector septum, or use an autosampler.
3. Begin the temperature program.
4. Record the chromatogram.
5. Analyze as indicated in the monograph.

**Liquid Chromatography**

LC, as used in the compendia, is synonymous with HPLC (both high-pressure and high-performance). LC is a separation technique based on a solid stationary phase and a liquid mobile phase.

## STATIONARY PHASE

Separations are achieved by partition, adsorption, or ion-exchange processes, depending on the type of stationary phase used. The most commonly used stationary phases are modified silica or polymeric beads. The beads are modified by the addition of long-chain hydrocarbons. The specific type of packing needed to complete an analysis is indicated by the “L” designation in the individual monograph (see also *Chromatographic Columns*). The size of the beads is often described in the monograph as well. Changes in the packing type and size are covered in *System Suitability* in this chapter.

## CHROMATOGRAPHIC COLUMN

The term “column” includes stainless steel, lined stainless steel, and polymeric columns, packed with a stationary phase. The length and inner diameter of the column affects the separation, and therefore typical column dimensions are included in the individual monograph. Changes to column dimensions are discussed in *System Suitability*. Compendial monographs do not include the name of appropriate columns; this omission avoids the appearance of endorsement of a vendor’s product and natural changes in the marketplace. See *Chromatographic Columns* for more information.

In LC procedures, a guard column may be used with the following requirements, unless otherwise indicated in the individual monograph: (a) the length of the guard column must be NMT 15% of the length of the analytical column, (b) the inner diameter must be the same or smaller than that of the analytical column, and (c) the packing material should be the same as the analytical column (e.g., silica) and contain the same bonded phase (e.g., C18). In any case, all system suitability requirements specified in the official procedure must be met with the guard column installed.

## MOBILE PHASE

The mobile phase is a solvent or a mixture of solvents, as defined in the individual monograph.

## APPARATUS

A liquid chromatograph consists of a reservoir containing the mobile phase, a pump to force the mobile phase through the system at high pressure, an injector to introduce the sample into the mobile phase, a chromatographic column, a detector, and a data collection device.

## GRADIENT ELUTION

The technique of continuously changing the solvent composition during the chromatographic run is called gradient elution or solvent programming. The gradient elution profile is presented in the individual monograph as a gradient table, which lists the time and proportional composition of the mobile phase at the stated time.

## PROCEDURE

1. Equilibrate the column and detector with mobile phase at the specified flow rate until a constant signal is received.
2. Inject a sample through the injector, or use an autosampler.
3. Begin the gradient program.
4. Record the chromatogram.
5. Analyze as directed in the monograph.

**Change to read:****CHROMATOGRAPHIC COLUMNS**

A complete list of packings (L), phases (G), and supports (S) used in USP–NF tests and assays is located in USP–NF, 15 (USP40) *Reagents, Indicators, and Solutions—Chromatographic Columns*. This list is intended to be a convenient reference for the chromatographer in identifying the pertinent chromatographic column specified in the individual monograph.

**DEFINITIONS AND INTERPRETATION OF CHROMATOGRAMS**

**Chromatogram:** A graphical representation of the detector response, concentration of analyte in the effluent, or other quantity used as a measure of effluent concentration versus effluent volume or time. In planar chromatography, chromatogram may refer to the paper or layer with the separated zones.

Figure 1 represents a typical chromatographic separation of two substances, 1 and 2.  $t_{R1}$  and  $t_{R2}$  are the respective retention times;  $h$  is the height,  $h/2$  is the half-height, and  $W_{h/2}$  is the width at half-height, for peak 1.  $W_1$  and  $W_2$  are the respective widths of peaks 1 and 2 at the baseline. Air peaks are a feature of gas chromatograms and correspond to the solvent front in LC. The retention time of these air peaks, or unretained components, is designated as  $t_M$ .

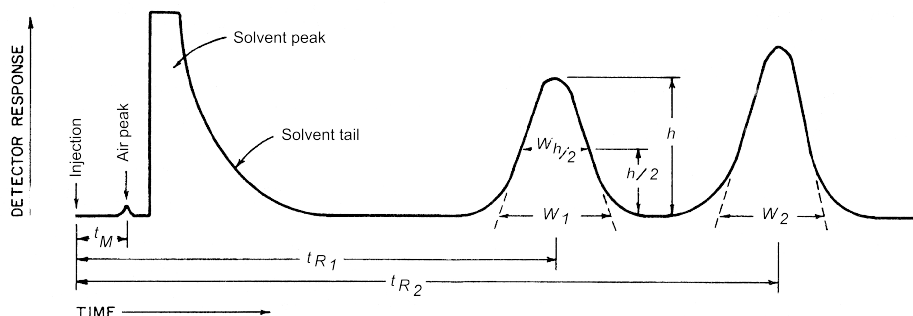


Figure 1. Chromatographic separation of two substances.

**Dwell volume (D):** Also known as “gradient delay volume”, is the volume between the point at which the eluents meet and the top of the column.

**Hold-up time ( $t_M$ ):** The time required for elution of an unretained component (see Figure 1, shown as an air or unretained solvent peak, with the baseline scale in minutes).

**Hold-up volume ( $V_M$ ):** The volume of mobile phase required for elution of an unretained component. It may be calculated from the hold-up time and the flow rate,  $F$ , in mL/min:

$$V_M = t_M \times F$$

In size-exclusion chromatography, the symbol  $V_0$  is used.

**Number of theoretical plates ( $N$ ):**<sup>1</sup> A measure of column efficiency. For Gaussian peaks, it is calculated by:

$$N = 16(t_R/W)^2$$

where  $t_R$  is the retention time of the substance, and  $W$  is the peak width at its base, obtained by extrapolating the relatively straight sides of the peak to the baseline. The value of  $N$  depends upon the substance being chromatographed as well as the operating conditions, such as the flow rate and temperature of the mobile phase or carrier gas, the quality of the packing, the uniformity of the packing within the column, and, for capillary columns, the thickness of the stationary phase film and the internal diameter and length of the column.

Where electronic integrators are used, it may be convenient to determine the number of theoretical plates, by the equation:

$$N = 5.54 \left( \frac{t_R}{W_{h/2}} \right)^2$$

where  $W_{h/2}$  is the peak width at half-height. However, in the event of dispute, only equations based on peak width at baseline are to be used.

**Peak:** The portion of the chromatographic recording of the detector response when a single component is eluted from the column. If separation is incomplete, two or more components may be eluted as one unresolved peak.

**Peak-to-valley ratio ( $p/v$ ):**  $p/v$  may be employed as a system suitability criterion in a test for related substances when baseline separation between two peaks is not achieved. *Figure 2* represents a partial separation of two substances, where  $H_p$  is the height above the extrapolated baseline of the minor peak and  $H_v$  is the height above the extrapolated baseline at the lowest point of the curve separating the minor and major peaks:

$$p/v = H_p/H_v$$

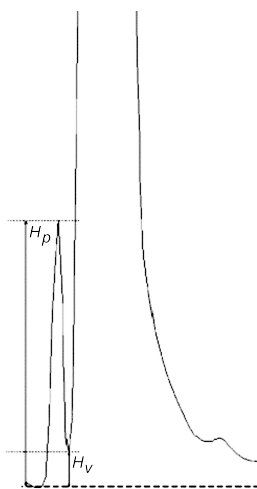


Figure 2. Peak-to-valley ratio determination.

**Relative retardation ( $R_{ret}$ ):** The ratio of the distance traveled by the analyte to the distance simultaneously traveled by a reference compound (see *Figure 3*) and is used in planar chromatography.

$$R_{ret} = b/c$$

<sup>1</sup> The parameters  $k$ ,  $N$ ,  $r$ , and  $r_G$  were developed for isothermal GC separations and isocratic HPLC separations. Because these terms are thermodynamic parameters, they are only valid for separations made at a constant temperature, mobile phase composition, and flow rate. However, for separations made with a temperature program or solvent gradient, these parameters may be used simply as comparative means to ensure that adequate chromatographic conditions exist to perform the methods as intended in the monographs.

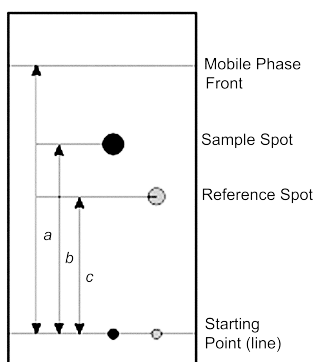


Figure 3. Typical planar chromatography.

**Relative retention ( $r$ ):**<sup>1</sup> The ratio of the adjusted retention time of a component relative to that of another used as a reference, obtained under identical conditions:

$$r = (t_{R2} - t_M)/(t_{R1} - t_M)$$

where  $t_{R2}$  is the retention time measured from the point of injection of the compound of interest;  $t_{R1}$  is the retention time measured from the point of injection of the compound used as reference; and  $t_M$  is the retention time of a nonretained marker defined in the procedure, all determined under identical experimental conditions on the same column.

**Relative retention time (RRT):** Also known as the “unadjusted relative retention”. Comparisons in *USP–NF* are normally made in terms of unadjusted relative retention, unless otherwise indicated.

$$\text{RRT} = t_{R2}/t_{R1}$$

The symbol  $r_c$  is also used to designate unadjusted relative retention values.

**Relative standard deviation in percentage (%RSD):**

$$\%RSD = \frac{100}{\bar{x}} \left( \frac{\sum_{i=1}^N (x_i - \bar{x})^2}{N - 1} \right)^{1/2}$$

**Resolution ( $R_S$ ):** The resolution is the separation of two components in a mixture, calculated by:

$$R_S = 2 \times (t_{R2} - t_{R1})/(W_1 + W_2)$$

where  $t_{R2}$  and  $t_{R1}$  are the retention times of the two components; and  $W_2$  and  $W_1$  are the corresponding widths at the bases of the peaks obtained by extrapolating the relatively straight sides of the peaks to the baseline.

Where electronic integrators are used, it may be convenient to determine the resolution, by the equation:

$$R_S = 1.18 \times (t_{R2} - t_{R1})/(W_{1,h/2} + W_{2,h/2})$$

**Retardation factor ( $R_f$ ):** The ratio of the distance traveled by the center of the spot to the distance simultaneously traveled by the mobile phase and is used in planar chromatography. Using the symbols in *Figure 3*:

$$R_f = b/a$$

**Retention factor ( $k$ ):**<sup>1</sup> Also known as the “capacity factor ( $k'$ )”. Defined as:

$$k = \frac{\text{amount of substance in stationary phase}}{\text{amount of substance in mobile phase}}$$

or

$$k = \frac{\text{time spent by substance in stationary phase}}{\text{time spent by substance in mobile phase}}$$

The  $k$  of a component may be determined from the chromatogram:

$$k = (t_R - t_M)/t_M$$

**Retention time ( $t_R$ ):** In LC and GC, the retention time,  $t_R$ , is defined as the time elapsed between the injection of the sample and the appearance of the maximum peak response of the eluted sample zone.  $t_R$  may be used as a parameter for identifi-

cation. Chromatographic retention times are characteristic of the compounds they represent but are not unique. The coincidence of retention times of a sample and a reference substance can be used as a partial criterion in construction of an identity profile, but may not be sufficient on its own to establish identity. Absolute retention times of a given compound may vary from one chromatogram to the next.

**Retention volume ( $V_R$ ):** The volume of mobile phase required for elution of a component. It may be calculated from the retention time and the flow rate in mL/min:

$$V_R = t_R \times F$$

**Separation factor ( $\alpha$ ):** The relative retention calculated for two adjacent peaks (by convention, the value of the separation factor is always >1):

$$\alpha = k_2/k_1$$

**Symmetry factor ( $A_s$ ):**<sup>2</sup> Also known as the “tailing factor”, of a peak (see Figure 4) is calculated by:

$$A_s = W_{0.05}/2f$$

where  $W_{0.05}$  is the width of the peak at 5% height and  $f$  is the distance from the peak maximum to the leading edge of the peak, the distance being measured at a point 5% of the peak height from the baseline.

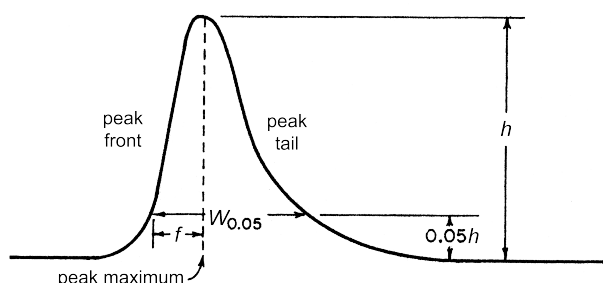


Figure 4. Asymmetrical chromatographic peak.

**Tailing factor ( $T$ ):** See *Symmetry factor*.

### Change to read:

## SYSTEM SUITABILITY

System suitability tests are an integral part of GC and LC methods. These tests are used to verify that the chromatographic system is adequate for the intended analysis.

The tests are based on the concept that the equipment, electronics, analytical operations, and samples analyzed constitute an integral system that can be evaluated as such.

Factors that may affect chromatographic behavior include the following:

- Composition, ionic strength, temperature, and apparent pH of the mobile phase
- Flow rate, column dimensions, column temperature, and pressure
- Stationary phase characteristics, including type of chromatographic support (particle-based or monolithic), particle or macropore size, porosity, and specific surface area
- Reverse-phase and other surface modification of the stationary phases, the extent of chemical modification (as expressed by end-capping, carbon loading, and others)

$R_s$  is a function of the number of theoretical plates,  $N$  (also referred to as efficiency),  $\alpha$ , and  $k$ . [NOTE—All terms and symbols are defined in *Definitions and Interpretation of Chromatograms*.] For a given stationary phase and mobile phase,  $N$  may be specified to ensure that closely eluting compounds are resolved from each other, to establish the general resolving power of the system, and to ensure that internal standards are resolved from the drug. This is a less reliable means to ensure resolution, as opposed to direct measurement. Column efficiency is, in part, a reflection of peak sharpness, which is important for the detection of trace components.

Replicate injections of a standard preparation or other standard solutions are compared to ascertain whether requirements for precision are met. Unless otherwise specified in the individual monograph, data from five replicate injections of the analyte are used to calculate the relative standard deviation (RSD), if the requirement is  $\leq 2.0\%$ ; data from six replicate injections are used if the RSD requirement is  $> 2.0\%$ .

<sup>2</sup> It is also a common practice to measure the “Asymmetry factor” as the ratio of the distance between the vertical line connecting the peak apex with the interpolated baseline and the peak front, and the distance between that line and the peak back measured at 10% of the peak height (see Figure 4), which would be  $(W_{0.10} - f_{0.10})/f_{0.10}$ . However, for the purposes of USP, only the formula ( $A_s$ ) as presented here is valid.



For the assay in a drug substance monograph, where the value is 100% for the pure substance, and no maximum RSD is stated, the maximum permitted %RSD is calculated for a series of injections of the reference solution:

$$\%RSD = KB\sqrt{n}/t_{90\%,n-1}$$

where  $K$  is a constant (0.349), obtained from the expression  $K = (0.6/\sqrt{2}) \times (t_{90\%,5}/\sqrt{6})$ , in which  $0.6/\sqrt{2}$  represents the required %RSD after six injections for  $B = 1.0$ ;  $B$  is the upper limit given in the definition of the individual monograph – 100%;  $n$  is the number of replicate injections of the reference solution ( $3 \leq n \leq 6$ ); and  $t_{90\%,n-1}$  is the Student's  $t$  at the 90% probability level (double sided) with  $n-1$  degrees of freedom.

Unless otherwise prescribed, the maximum permitted RSD does not exceed the appropriate value given in *Table 1* of repeatability requirements. This requirement does not apply to tests for related substances.

**Table 1. RSD Requirements**

B (%)	Number of Individual Injections			
	3	4	5	6
	Maximum Permitted RSD			
2.0	0.41	0.59	0.73	0.85
2.5	0.52	0.74	0.92	1.06
3.0	0.62	0.89	1.10	1.27

$A_s$ , a measure of peak symmetry, is unity for perfectly symmetrical peaks; its value increases as tailing becomes more pronounced (see *Figure 4*). In some cases, values less than unity may be observed. As peak symmetry moves away from values of 1, integration, and hence precision, become less reliable.

The signal-to-noise (S/N) ratio is a useful system suitability parameter. The S/N ratio is calculated as follows:

$$S/N \text{ ratio} = 2H/h$$

where  $H$  is the height of the peak measured from the peak apex to a baseline extrapolated over a distance  $\geq 5$  times the peak width at its half-height; and  $h$  is the difference between the largest and smallest noise values observed over a distance  $\geq 5$  times the width at the half-height of the peak and, if possible, situated equally around the peak of interest (see *Figure 5*).

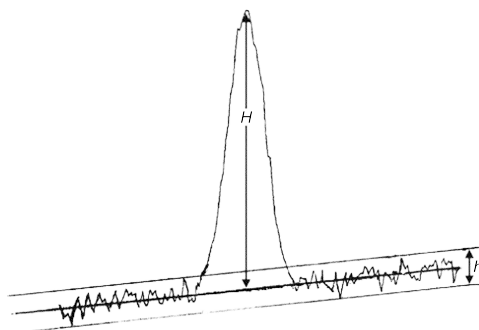


Figure 5. Noise and chromatographic peak, components of the S/N ratio.

These system suitability tests are performed by collecting data from replicate injections of standard or other solutions as specified in the individual monograph.

The specification of definitive parameters in a monograph does not preclude the use of other suitable operating conditions.

Adjustments to the specified chromatographic system may be necessary in order to meet system suitability requirements. Adjustments to chromatographic systems performed in order to comply with system suitability requirements are not to be made in order to compensate for column failure or system malfunction. Adjustments are permitted only when suitable standards (including Reference Standards) are available for all compounds used in the suitability test, and the adjustments or column change yields a chromatogram that meets all the system suitability requirements specified in the official procedure.

If adjustments of operating conditions are necessary in order to meet system suitability requirements, each of the items in the following list is the maximum variation that can be considered, unless otherwise directed in the monograph; these changes may require additional verification data. To verify the suitability of the method under the new conditions, assess the relevant analytical performance characteristics potentially affected by the change. Multiple adjustments can have a cumulative effect on the performance of the system and are to be considered carefully before implementation. In some circumstances, it may be desirable to use an HPLC column with different dimensions to those prescribed in the official procedure (different length, internal diameter, and/or particle size). In either case, changes in the chemical characteristics ("L" designation) of the stationary phase will be considered a modification to the method and will require full validation. Adjustments to the composition of the mobile phase in gradient elution may cause changes in selectivity and are not recommended. If adjustments are necessary, a change in column packing (maintaining the same chemistry), the duration of an initial isocratic hold (when prescribed), and/or the dwell volume are allowed. Additional allowances for gradient adjustment are noted in the text below.

**pH of mobile phase (HPLC):** The pH of the aqueous buffer used in the preparation of the mobile phase can be adjusted to within  $\pm 0.2$  units of the value or range specified. Applies to both gradient and isocratic separations.

**Concentration of salts in buffer (HPLC):** The concentration of the salts used in the preparation of the aqueous buffer employed in the mobile phase can be adjusted to within  $\pm 10\%$  if the permitted pH variation (see above) is met. Applies to both gradient and isocratic separations.

**Ratio of components in mobile phase (HPLC):** The following adjustment limits apply to minor components of the mobile phase (specified as  $\leq 50\%$ ). The amounts of these components can be adjusted by  $\pm 30\%$  relative. However, the change in any component cannot exceed  $\pm 10\%$  absolute (i.e., in relation to the total mobile phase). Adjustment can be made to one minor component in a ternary mixture. Examples of adjustments for binary and ternary mixtures are given below.

#### Binary mixtures

*Specified ratio of 50:50:* 30% of 50 is 15% absolute, but this exceeds the maximum permitted change of  $\pm 10\%$  absolute in either component. Therefore, the mobile phase ratio may be adjusted only within the range of 40:60–60:40.

*Specified ratio of 2:98:* 30% of 2 is 0.6% absolute. Therefore, the maximum allowed adjustment is within the range of 1.4:98.6–2.6:97.4.

#### Ternary mixtures

*Specified ratio of 60:35:5:* For the second component, 30% of 35 is 10.5% absolute, which exceeds the maximum permitted change of  $\pm 10\%$  absolute in any component. Therefore, the second component may be adjusted only within the range of 25%–45% absolute. For the third component, 30% of 5 is 1.5% absolute. In all cases, a sufficient quantity of the first component is used to give a total of 100%. Therefore, mixture ranges of 50:45:5–70:25:5 or 58.5:35:6.5–61.5:35:3.5 would meet the requirement.

**Wavelength of UV-visible detector (HPLC):** Deviations from the wavelengths specified in the procedure are not permitted. The procedure specified by the detector manufacturer, or another validated procedure, is used to verify that error in the detector wavelength is, at most,  $\pm 3$  nm.

#### Stationary phase

*Column length (GC):* Can be adjusted by as much as  $\pm 70\%$ .

*Column length (HPLC):* See *Particle size (HPLC)* below.

*Column inner diameter (HPLC):* Can be adjusted if the linear velocity is kept constant. See *Flow rate (HPLC)*.

*Column inner diameter (GC):* Can be adjusted by as much as  $\pm 50\%$ .

*Film thickness (capillary GC):* Can be adjusted by as much as  $-50\%$  to  $100\%$ .

**Particle size (HPLC):** For isocratic separations, the particle size and/or the length of the column may be modified provided that the ratio of the column length ( $L$ ) to the particle size ( $dp$ ) remains constant or into the range between  $-25\%$  and  $50\%$  of the prescribed  $L/dp$  ratio. Alternatively (as for the application of particle-size adjustment to superficially porous particles), other combinations of  $L$  and  $dp$  can be used provided that the number of theoretical plates ( $N$ ) is within  $-25\%$  to  $50\%$ , relative to the prescribed column. Caution should be used when the adjustment results in a higher number of theoretical plates that generate smaller peak volumes, which may require adjustments to minimize extra-column band broadening by factors such as instrument plumbing, detector cell volume and sampling rate, and injection volume. ■<sub>1S</sub> (USP40) For gradient separations, changes in length, column inner diameter, and particle size are not allowed.

**Particle size (GC):** Changing from a larger to a smaller or from a smaller to a larger particle size GC mesh support is acceptable if the chromatography meets the requirements of system suitability and the same particle size range ratio is maintained. The particle size range ratio is defined as the diameter of the largest particle divided by the diameter of the smallest particle.

**Flow rate (GC):** The flow rate can be adjusted by as much as  $\pm 50\%$ . ■[NOTE—When the monograph specifies a linear velocity parameter, the allowed velocity adjustment is between  $+50\%$  and  $-25\%$ , provided the carrier gas system can be maintained under control at the desired set points.] ■<sub>1S</sub> (USP40)

**Flow rate (HPLC):** When the particle size is changed, the flow rate may require adjustment, because smaller-particle columns will require higher linear velocities for the same performance (as measured by reduced plate height). Flow rate changes for both a change in column diameter and particle size can be made by:

$$F_2 = F_1 \times [(dc_2^2 \times dp_1)/(dc_1^2 \times dp_2)]$$

where  $F_1$  and  $F_2$  are the flow rates for the original and modified conditions,  $dc_1$  and  $dc_2$  are the respective column diameters, and  $dp_1$  and  $dp_2$  are the particle sizes.

When a change is made from  $\geq 3\text{-}\mu\text{m}$  to  $< 3\text{-}\mu\text{m}$  particles in isocratic separations, an additional increase in linear velocity (by adjusting flow rate) may be justified, provided that the column efficiency does not drop by  $> 20\%$ . Similarly, a change from  $< 3\text{-}\mu\text{m}$  to  $\geq 3\text{-}\mu\text{m}$  particles may require additional reduction of linear velocity (flow rate) to avoid reduction in column efficiency by  $> 20\%$ . Changes in  $F$ ,  $dc$ , and  $dp$  are not allowed for gradient separations.

Additionally, the flow rate can be adjusted by  $\pm 50\%$  (isocratic only).

*Examples:* Adjustments in column length, internal diameter, particle size, and flow rate can be used in combination to give equivalent conditions (same  $N$ ), but with differences in pressure and run time. Table 2 lists some of the more popular column configurations to give equivalent efficiency ( $N$ ), by adjusting these variables.

Table 2. Column Configurations

Length (L, mm)	Column Diameter ( $d_c$ , mm)	Particle Size ( $d_p$ , $\mu\text{m}$ )	Relative Values				
			$L/d_p$	$F$	$N$	Pressure	Run Time
250	4.6	10	25,000	0.5	0.8	0.2	3.3
150	4.6	5	30,000	1.0	1.0	1.0	1.0
150	2.1	5	30,000	0.2	1.0	1.0	1.0
100	4.6	3.5	28,600	1.4	1.0	1.9	0.5
100	2.1	3.5	28,600	0.3	1.0	1.9	0.5
75	4.6	2.5	30,000	2.0	1.0	4.0	0.3
75	2.1	2.5	30,000	0.4	1.0	4.0	0.3
50	4.6	1.7	29,400	2.9	1.0	8.5	0.1
50	2.1	1.7	29,400	0.6	1.0	8.5	0.1

For example, if a monograph specifies a 150-mm  $\times$  4.6-mm; 5- $\mu\text{m}$  column operated at 1.5 mL/min, the same separation may be expected with a 75-mm  $\times$  2.1-mm; 2.5- $\mu\text{m}$  column operated at 1.5 mL/min  $\times$  0.4 = 0.6 mL/min, along with a pressure increase of about four times and a reduction in run time to about 30% of the original.

**Injection volume (HPLC):** The injection volume can be adjusted as far as it is consistent with accepted precision, linearity, and detection limits. Note that excessive injection volume can lead to unacceptable band broadening, causing a reduction in  $N$  and resolution, which applies to both gradient and isocratic separations.

**Injection volume and split volume (GC):** The injection volume and split volume may be adjusted if detection and repeatability are satisfactory.

**Column temperature (HPLC):** The column temperature can be adjusted by as much as  $\pm 10^\circ$ . Column thermostating is recommended to improve control and reproducibility of retention time, which applies to both gradient and isocratic separations.

**Oven temperature (GC):** The oven temperature can be adjusted by as much as  $\pm 10\%$ .

**Oven temperature program (GC):** Adjustment of temperatures is permitted as stated above. When the specified temperature must be maintained or when the temperature must be changed from one value to another, an adjustment of up to  $\pm 20\%$  is permitted.

Unless otherwise directed in the monograph, system suitability parameters are determined from the analyte peak.

Measured values of  $R_r$ ,  $R_p$ , or  $t_r$  for the sample substance do not deviate from the values obtained for the reference compound and mixture by more than the statistically determined reliability estimates from replicate assays of the reference compound. RRT may be provided in monographs for informational purposes only to aid in peak identification. There are no acceptance criteria applied to RRT.

Suitability testing is used to ascertain the effectiveness of the final operating system, which should be subjected to this testing. Make injections of the appropriate preparation(s), as required, in order to demonstrate adequate system suitability (as described in the *Chromatographic system* section of the method in a monograph) throughout the run.

The preparation can be a standard preparation or a solution containing a known amount of analyte and any additional materials (e.g., excipients or impurities) useful in controlling the analytical system. Whenever there is a significant change in the chromatographic system (equipment, mobile phase component, or other components) or in a critical reagent, system suitability is to be reestablished. No sample analysis is acceptable unless the suitability of the system has been demonstrated.

## QUANTITATION

During quantitation, disregard peaks caused by solvents and reagents or arising from the mobile phase or the sample matrix.

In the linear range, peak areas and peak heights are usually proportional to the quantity of compound eluting. The peak areas and peak heights are commonly measured by electronic integrators but may be determined by more classical approaches. Peak areas are generally used but may be less accurate if peak interference occurs. The components measured are separated from any interfering components. Peak tailing and fronting is minimized, and the measurement of peaks on tails of other peaks are avoided when possible.

Although comparison of impurity peaks with those in the chromatogram of a standard at a similar concentration is preferred, impurity tests may be based on the measurement of the peak response due to impurities and expressed as a percentage of the area of the drug peak. The standard may be the drug itself at a level corresponding to, for example, 0.5% impurity, assuming similar peak responses. When impurities must be determined with greater certainty, use a standard of the impurity itself or apply a correction factor based on the response of the impurity relative to that of the main component.

### External Standard Method

The concentration of the component(s) quantified is determined by comparing the response(s) obtained with the sample solution to the response(s) obtained with a standard solution.

### Internal Standard Method

Equal amounts of the internal standard are introduced into the sample solution and a standard solution. The internal standard is chosen so that it does not react with the test material and does not contain impurities with the same retention time as that of the analytes, and is stable and resolved from the component(s) quantified (analytes). The concentrations of the analytes are determined by comparing the ratios of their peak areas or peak heights and the internal standard in the sample solution with the ratios of their peak areas or peak heights and the internal standard in the standard solution.

### Normalization Procedure


The percent content of a component of the test material is calculated by determining the area of the corresponding peak as a percentage of the total area of all the peaks, excluding those due to solvents or reagents or arising from the mobile phase or the sample matrix and those at or below the limit at which they can be disregarded.

### Calibration Procedure

The relationship between the measured or evaluated signal  $y$  and the quantity (e.g., concentration or mass) of substance  $x$  is determined, and the calibration function is calculated. The analytical results are calculated from the measured signal or evaluated signal of the analyte and its position on the calibration curve.

In tests for impurities for both the external standard method, when a dilution of the sample solution is used for comparison, and the normalization procedure, any correction factors indicated in the monograph are applied (e.g., when the relative response factor is outside the range of 0.8–1.2).

When the impurity test prescribes the total of impurities or there is a quantitative determination of an impurity, choice of an appropriate threshold setting and appropriate conditions for the integration of the peak areas is important. In such tests the limit at or below which a peak is disregarded is generally 0.05%. Thus, the threshold setting of the data collection system corresponds to at least half of this limit. Integrate the peak area of any impurity that is not completely separated from the principal peak, preferably by valley-to-valley extrapolation (tangential skim).

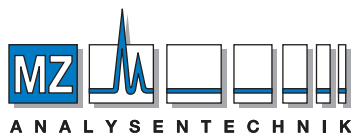

MZ-Analysentechnik GmbH, Barcelona-Allee 17• D-55129 Mainz  
Tel +49 6131 880 96-0, Fax +49 6131 880 96-20  
e-mail: [info@mz-at.de](mailto:info@mz-at.de), [www.mz-at.de](http://www.mz-at.de)