


## **COSMOSIL Application Data**

Column: 5Diol-\*\*\*-II  
Column size: 7.5mmI.D.-300mm  
Mobile phase: 20mmol/l Phosphate buffer(pH7.0),  
100mmol/l Na<sub>2</sub>SO<sub>4</sub>  
Flow rate: 1.0 ml/min  
Temperature: 30°C  
Detection: UV260nm

Sample: 1; dsDNA, 1000bp  
2; dsDNA, 100bp  
3; ssDNA, 20mer  
4; Uracil

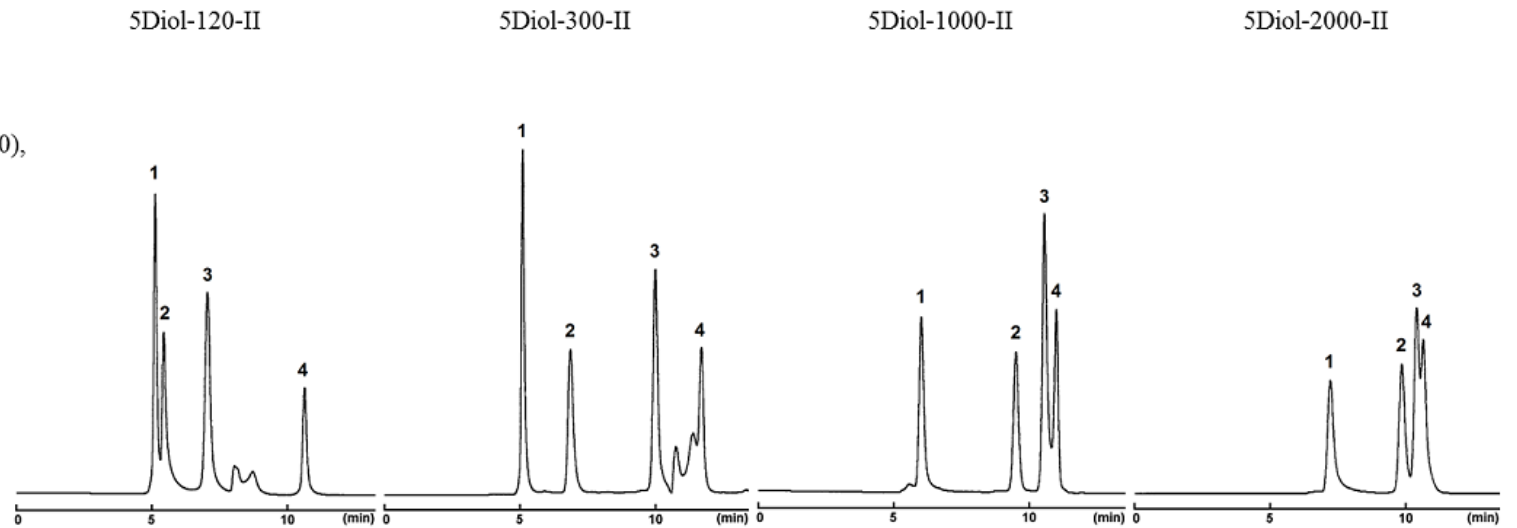

NACALAI TESQUE, INC

AP-1481

**AUTHORIZED DISTRIBUTOR**

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: [info@mz-at.de](mailto:info@mz-at.de), [www.mz-at.de](http://www.mz-at.de)