



COSMOSIL Application Data

Column: COSMOSIL RNA-RP1
Column size: 2.0 mm I.D.-100 mm
Mobile phase: A; 100 mmol/L Triethylamine, 100 mmol/L Acetic Acid (pH7)
B; Acetonitrile
B conc. 9→14 % (0→20 min), 14 % (20-24 min)
Flow rate: 0.2 ml/min
Temperature: 65 °C
Detection: UV260 nm

Sample: mRNA (0.8 knt)
(413 ng/μl)
Inj. Vol: 0.25 μl

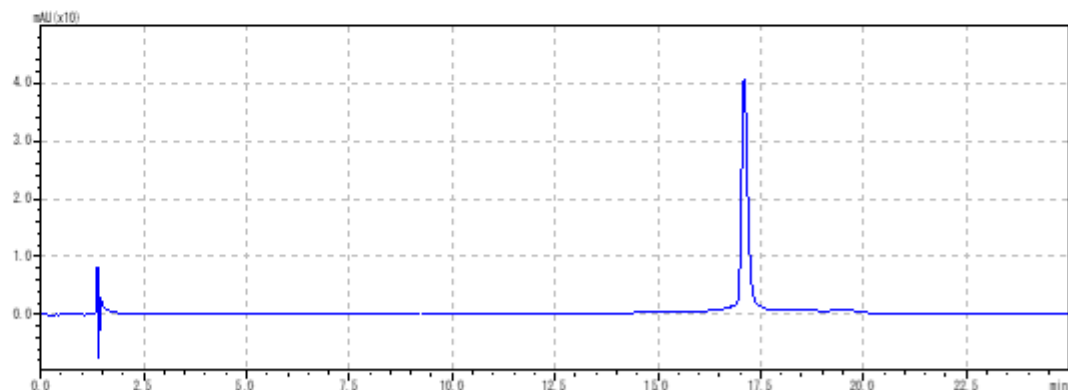

NACALAI TESQUE, INC

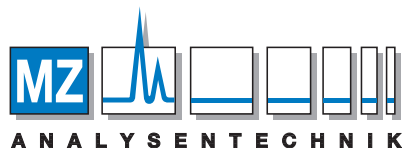

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de

AP-1873