


COSMOSIL Application Data

Column: COSMOSIL RNA-RP1
Column size: 4.6 mm I.D.-100 mm
Mobile phase: A; 100 mmol/L Triethylamine, 100 mmol/L Acetic Acid (pH7)
B; Acetonitrile
B conc. 5→13 % (0→20 min), 13 % (20-24 min)
Flow rate: 1.0 ml/min
Temperature: 65 °C
Detection: UV260 nm

Sample: ssRNA Ladder,
mixture
Inj. Vol: 4.0 µl

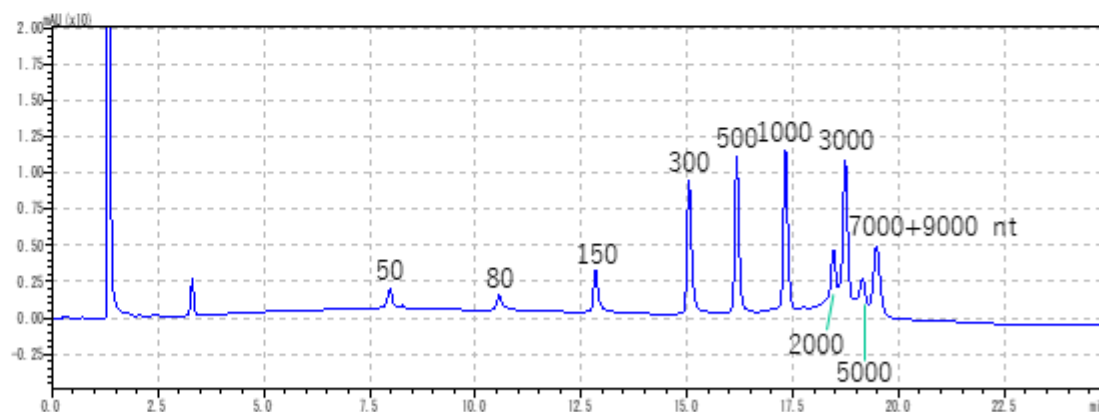

NACALAI TESQUE, INC