


COSMOSIL Application Data


Column: COSMOSIL RNA-SEC-****
Column size: 4.6 mmI.D.-250 mm
Mobile phase: 20 mmol/L Phosphate buffer (pH7.0)
Flow rate: 0.3 ml/min
Temperature: 40 °C
Detection: UV260 nm

Sample: 100 bp DNA Ladder (130 ng/μl)
Inj. Vol: 4.0 μl

COSMOSIL
RNA-SEC-1000

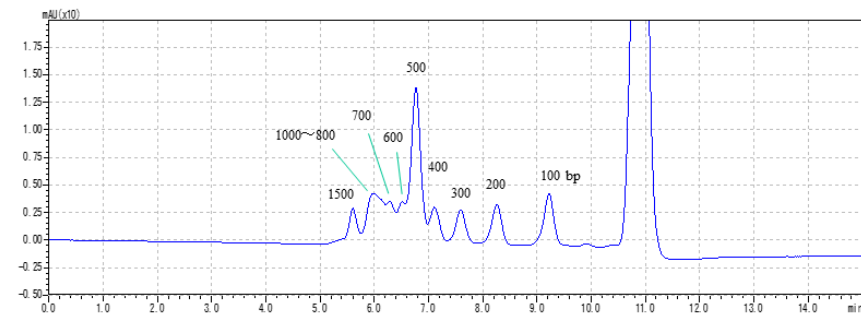

COSMOSIL
RNA-SEC-2000

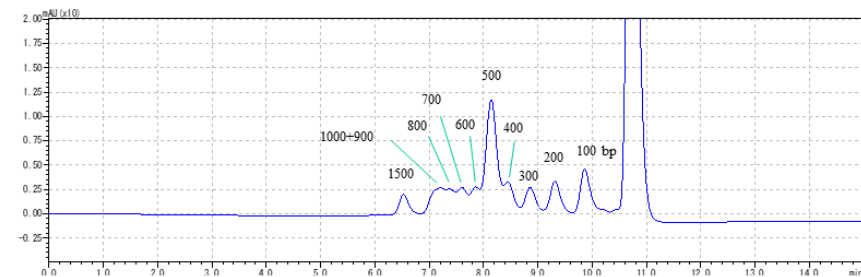

NACALAI TESQUE, INC

AP-1853