


## COSMOSIL Application Data


Column: COSMOSIL RNA-SEC-\*\*\*\*  
Column size: 4.6 mmI.D.-250 mm  
Mobile phase: 20 mmol/L Phosphate buffer (pH7.0)  
Flow rate: 0.3 ml/min  
Temperature: 40 °C  
Detection: UV260 nm

Sample: 20 bp DNA Ladder (200 ng/μl)  
Inj. Vol: 2.5 μl

COSMOSIL  
RNA-SEC-1000

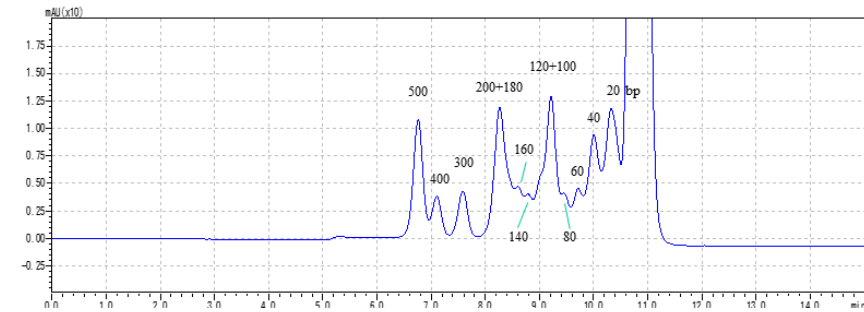


COSMOSIL  
RNA-SEC-2000

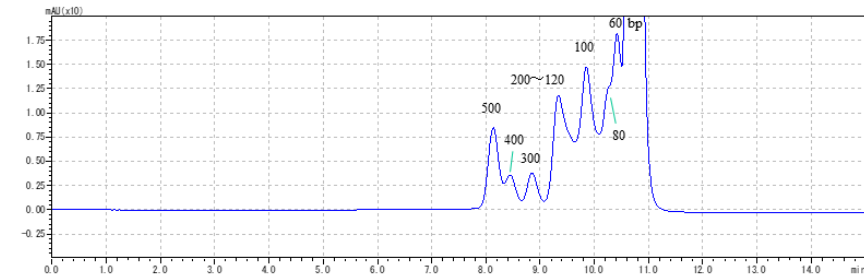

NACALAI TESQUE, INC

AP-1854

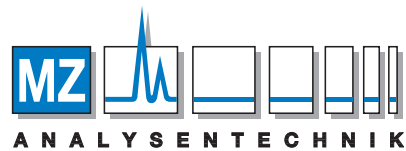

### AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: [info@mz-at.de](mailto:info@mz-at.de), [www.mz-at.de](http://www.mz-at.de)