


COSMOSIL Application Data


Column: COSMOSIL RNA-SEC-****
Column size: 4.6 mm I.D. x 250 mm
Mobile phase: 20 mmol/L Phosphate buffer (pH 7.0)
Flow rate: 0.3 ml/min
Temperature: 65 °C
Detection: UV260 nm

Sample: ssRNA Ladder (200 ng/μl)
Inj. Vol: 2.0 μl

COSMOSIL
RNA-SEC-1000

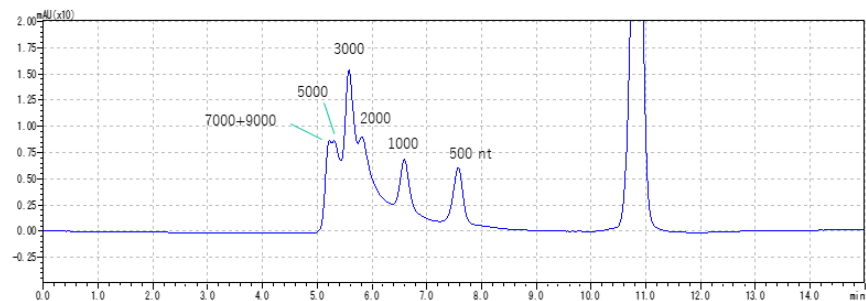


COSMOSIL
RNA-SEC-2000

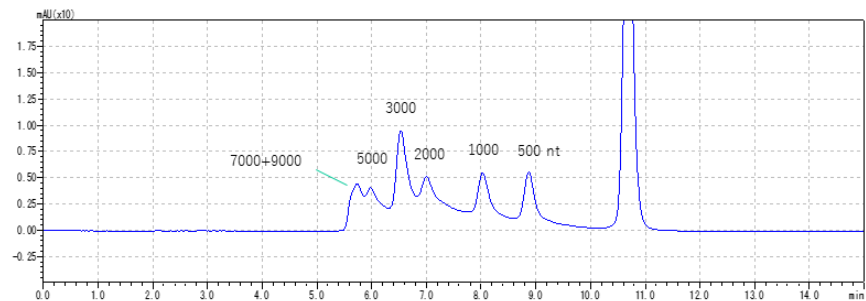

NACALAI TESQUE, INC

AP-1861

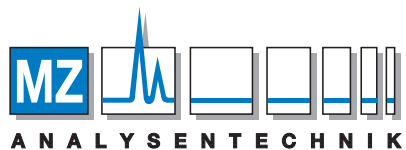

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de