



Fat-soluble vitamins

Separation of vitamins A, D and E on Kromasil EternityXT C18


1 = retinol

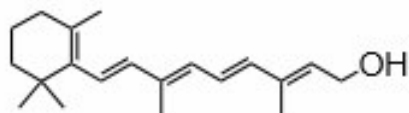



2 = ergocalciferol (vitamin D2)

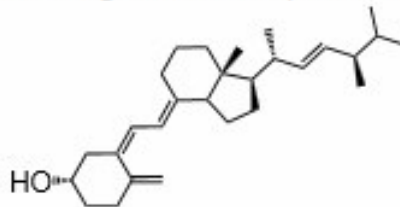


3 = cholecalciferol (vitamin D3)

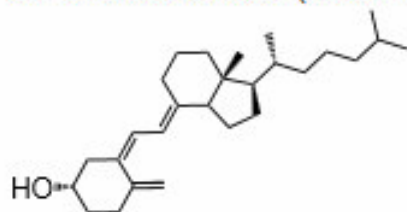

4 = δ -tocopherol (vitamin E)

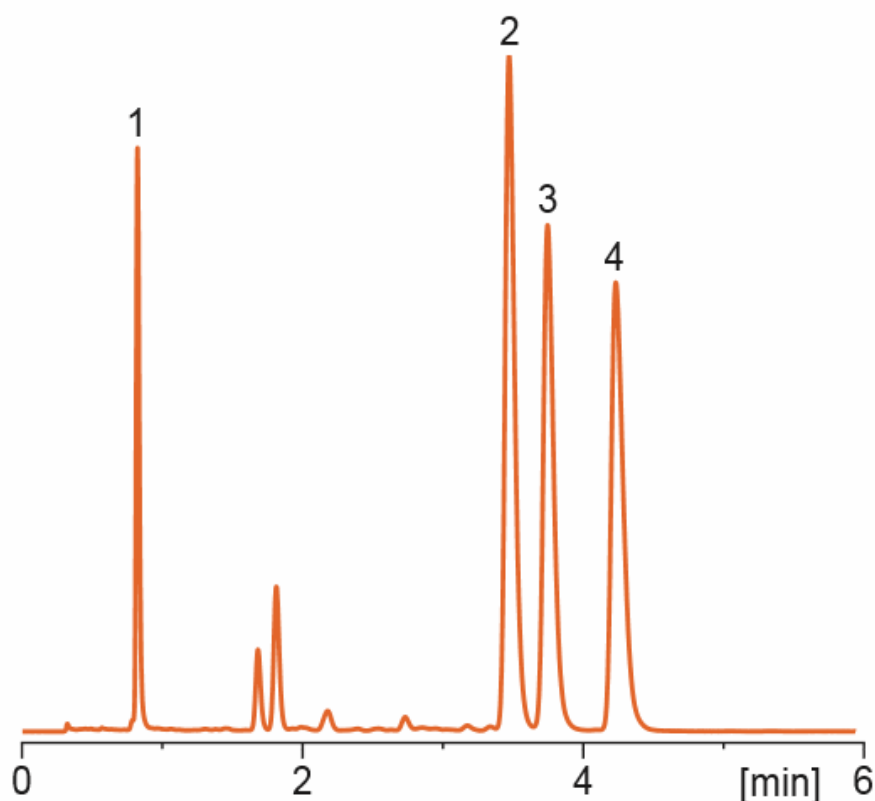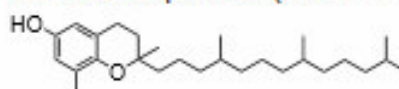

Column: Kromasil EternityXT-1.8-C18 2.1 × 100 mm 
Part number: XF1CLD10
Eluent: acetonitrile
Flow rate: 0.7 ml/min
Temperature: 35 °C
Detection: UV 265 nm

AUTHORIZED DISTRIBUTOR

MZ-Analyse-Technik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de

