

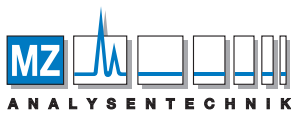


Data provided by : Dr.Oka, Aichi Prefectural Institute of Public Health

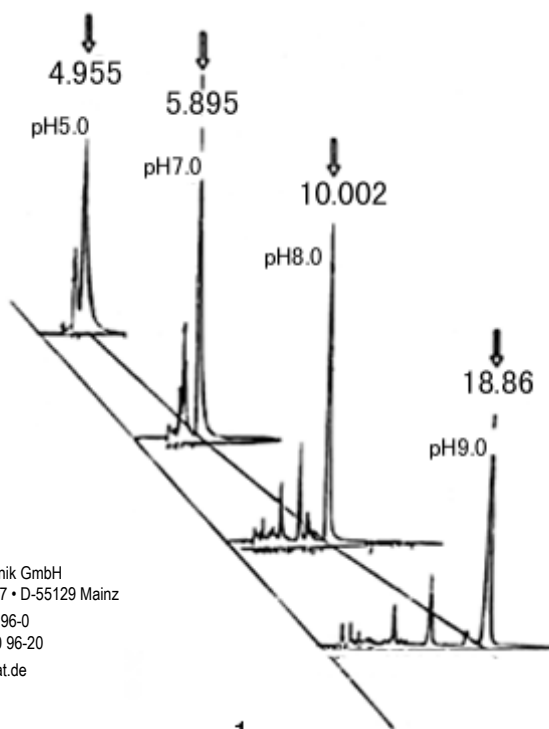

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de

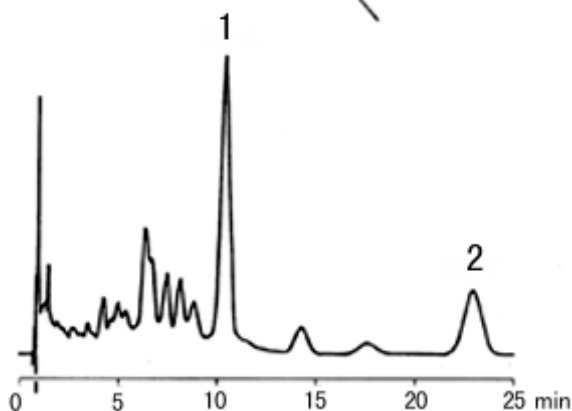

1 Bacitracin-A

2 Bacitracin-F

HPLC Conditions 1

Column CAPCELL PAK C₁₈ AG120 S5 : 4.6 mm i.d.×250 mm

Mobile phase 50mmol/L NaH₂PO₄ / CH₃CN

Detect UV 210 nm

HPLC Conditions 2

Column CAPCELL PAK C₁₈ AG120 S5 : 4.6 mm i.d.×250 mm

Mobile phase CH₃OH / 40mmol/L Na₂HPO₄ = 60 / 40, pH9-10

Detect

UV 234 nm