


- | | | | | | |
|---|----------------------------|---|--------------------|---|-----------------|
| 1 | Ribonuclease-A (14K) | 2 | Cytochrome-C (12K) | 3 | Lysozyme (14K) |
| 4 | Bovine serum albumin (67K) | 5 | Myoglobin (17K) | 6 | Ovalbumin (43K) |

HPLC Conditions

Column	CAPCELL PAK C ₁₈ SG300 S5 : 4.6 mm i.d.×150 mm
Mobile phase	(A) 0.1vol% TFA, H ₂ O (B) 0.1vol% TFA, CH ₃ CN B 15% - 60% (15min) Gradient
Flow rate	1.5 mL/min
Temperature	40 °C
Detect	UV 214 nm

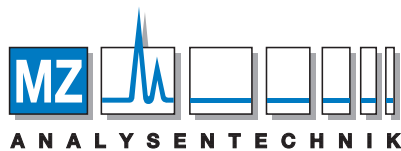

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de