


1	Ribonuclease-A	2	Cytochrome-C	3	Lysozyme
	2 µg		2 µg		2 µg
4	Bovine serum albumin	5	Myoglobin	6	Ovalbumin
	2 µg		2 µg		2 µg

HPLC Conditions

Column	CAPCELL PAK C18 UG300 S5 : 1.5 mm i.d.×150 mm
Mobile phase	(A) 0.1vol% TFA, H <sub>2</sub> O (B) 0.1vol% TFA, CH <sub>3</sub> CN B 15% - 60% (15min) - 60% (25min) Gradient
Flow rate	280 µL/min
Temperature	40 °C
Detect	UV 214 nm
Inj. vol.	1.0 µL

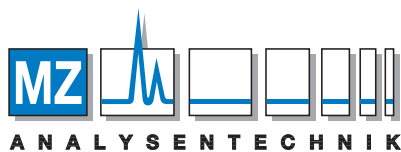

**AUTHORIZED DISTRIBUTOR**

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz  
 Tel +49 6131 880 96-0, Fax +49 6131 880 96-20  
 e-mail: info@mz-at.de, www.mz-at.de