


## Nitrosamines Stabilwax®-DB

AUTHORIZED DISTRIBUTOR

1. solvent
2. N-nitrosodimethylamine
3. N-nitrosodiethylamine
4. N-nitrosodiisopropylamine
5. N-nitrosodipropylamine
6. N-nitrosobutylamine
7. N-nitrosopiperidine
8. N-nitrosopyrrolidine
9. N-nitrosomorpholine

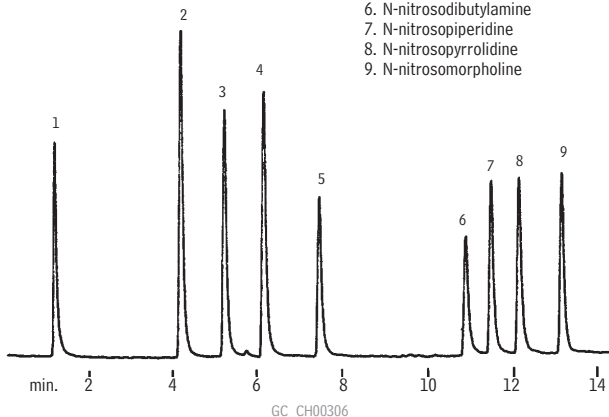

60m, 0.53mm ID, 1.0µm Stabilwax®-DB (cat.# 10858)

Direct injection of nitrosamines

Conc.: 1.0µg/mL

Oven temp.: 100°C (hold 1 min.) to 170°C @ 5°C/min.

Inj. & det. temp.: 200°C

Carrier gas: helium

Linear velocity: 100cm/sec. (flow rate: 15cc/min.)

Detector: TSD