

Perfluorinated Compounds (PFCs) on Raptor C18 (1.8 µm)


Column Raptor C18 (cat.# 9304252)
Dimensions: 50 mm x 2.1 mm ID
Particle Size: 1.8 µm
Pore Size: 90 Å
Temp.: 40 °C
Sample Diluent: 80:20 Water:methanol
Conc.: 5-250 ng/mL
Inj. Vol.: 5 µL
Mobile Phase
A: 5 mM Ammonium acetate in water
B: Methanol

Time (min)	Flow (mL/min)	%A	%B
0.00	0.4	80	20
6.00	0.4	5	95
6.50	0.4	5	95
6.51	0.4	80	20
8.00	0.4	80	20

Detector MS/MS
Ion Mode: ESI-
Mode: MRM
Instrument UHPL


Peaks	tr (min)	Conc. (ng/mL)	Precursor Ion	Product Ion
1. Perfluorobutanoic acid (PFBA)	0.90	5	213.07	169.07
2. Perfluoropentanoic acid (PFPeA)	1.83	5	263.10	219.11
3. Perfluorobutane sulfonic acid (PFBS)	2.09	5	299.10	79.99
4. 1H, 1H, 2H, 2H-perfluorohexane sulfonic acid (4:2 FTS)	2.63	5	327.10	307.08
5. Perfluorohexanoic acid (PFHxA)	2.70	5	313.10	269.12
6. Perfluoropentane sulfonic acid (PFPeS)	2.84	5	349.10	79.98
7. Perfluoroheptanoic acid (PFHpA)	3.33	5	363.16	319.09
8. Perfluorohexane sulfonic acid (PFHxS)	3.39	5	399.13	79.98
9. 2-Perfluorohexyl ethanoic acid (FHEA)	3.46	250	377.13	293.07
10. 1H, 1H, 2H, 2H-perfluorooctane sulfonic acid (6:2 FTS)	3.78	5	427.17	407.18
11. Perfluorooctanoic acid (PFOA)	3.81	5	413.16	369.10
12. Perfluoroheptane sulfonic acid (PFHpS)	3.84	5	449.17	79.98
13. Perfluorononanoic acid (PFNA)	4.19	5	463.16	419.19
14. Perfluorooctane sulfonic acid (PFOS)	4.21	5	499.17	79.98
15. 2-Perfluorooctyl ethanoic acid (FOEA)	4.31	250	477.17	393.15
16. 1H, 1H, 2H, 2H-perfluorodecane sulfonic acid (8:2 FTS)	4.51	5	527.17	507.16
17. Perfluorodecanoic acid (PFDA)	4.52	5	513.17	469.16
18. Perfluorononane sulfonic acid (PFNS)	4.53	5	549.17	79.98
19. N-methyl perfluorooctane sulfonamidoacetic acid (N-MeFOSAA)	4.66	5	570.20	419.17
20. Perfluorodecane sulfonic acid (PFDS)	4.79	5	599.17	79.98
21. N-ethyl perfluorooctane sulfonamidoacetic acid (N-EtFOSAA)	4.79	5	584.20	419.18
22. Perfluorooctane sulfonamide (FOSA)	4.79	5	498.17	77.97
23. Perfluoroundecanoic acid (PFUnA)	4.80	5	563.23	519.24
24. 2-Perfluorodecyl ethanoic acid (FDEA)	4.91	250	577.23	493.18
25. Perfluorododecanoic acid (PFDoA)	5.04	5	613.23	569.19
26. 1H, 1H, 2H, 2H-perfluorododecane sulfonic acid (10:2 FTS)	5.05	5	627.23	607.20
27. Perfluorotridecanoic acid (PFTrDA)	5.25	5	663.23	619.21
28. N-methyl perfluorooctane sulfonamide (N-MeFOSA)	5.31	5	512.17	169.07
29. N-methyl perfluorooctane sulfonamidoethanol (N-MeFOSE)	5.31	50	616.27	59.11
30. Perfluorotetradecanoic acid (PFTeDA)	5.44	5	713.23	669.23
31. N-ethyl perfluorooctane sulfonamide (N-EtFOSA)	5.47	5	526.23	169.06
32. N-ethyl perfluorooctane sulfonamidoethanol (N-EtFOSE)	5.47	50	630.33	59.11
33. Perfluorohexadecanoic acid (PFHxDA)	5.74	5	813.30	769.28
34. Perfluorooctadecanoic acid (PFODa)	5.98	5	913.30	869.32

- High chromatographic efficiency
- Fast 8-minute analysis time

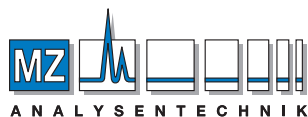

MZ-Analysentechnik GmbH
 Barcelona-Allee 17 • D-55129 Mainz
 Tel +49 6131 880 96-0
 Fax +49 6131 880 96-20
 e-mail: info@mz-at.de
 www.mz-at.de


AUTHORIZED DISTRIBUTOR

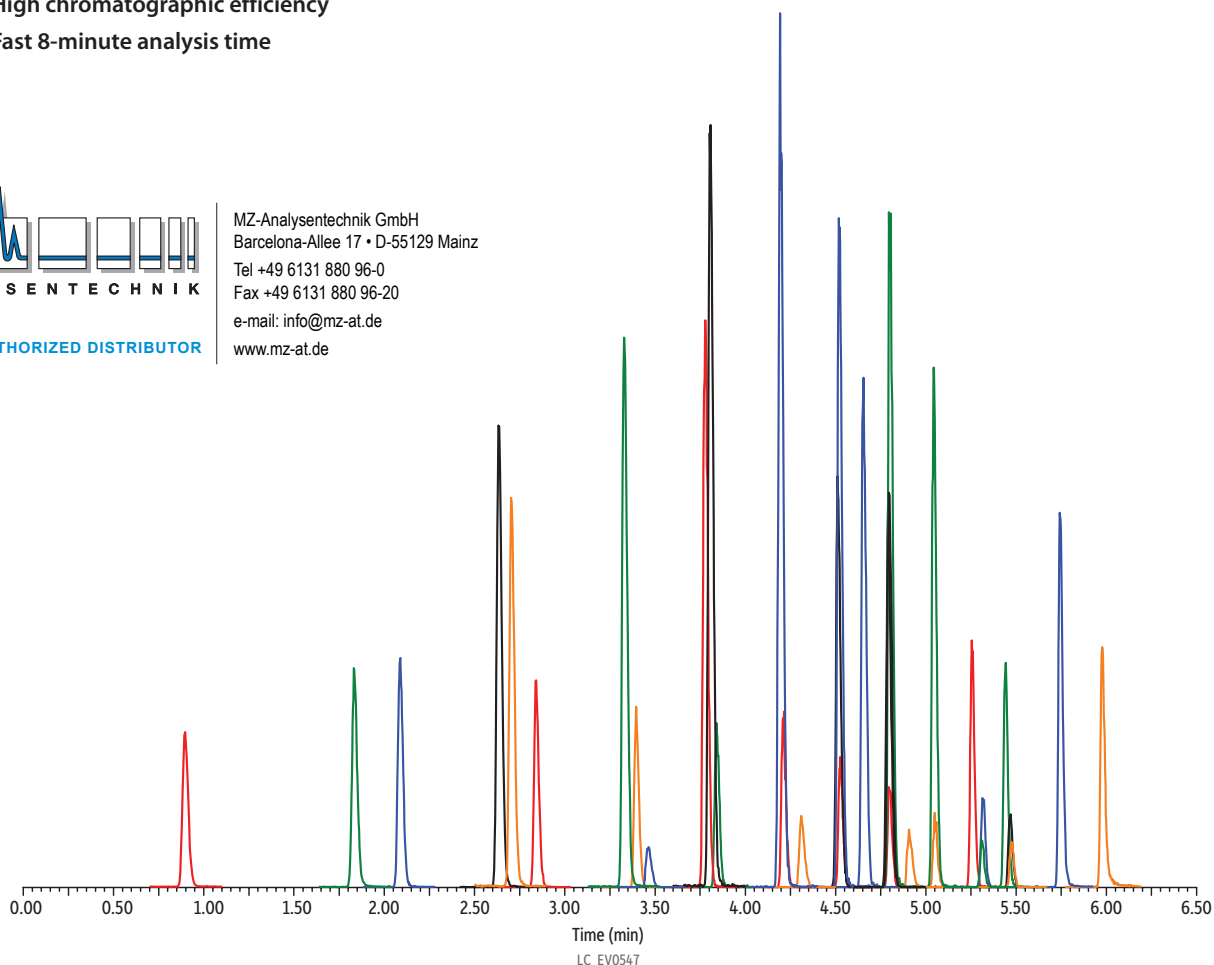


© 2017 Restek Corporation. All rights reserved.
 www.restek.com

Restek patents and trademarks are the property of Restek Corporation. (See www.restek.com/Patents-Trademarks for full list.) Other trademarks in Restek literature or on its website are the property of their respective owners. Restek registered trademarks are registered in the U.S. and may also be registered in other countries.


Perfluoroalkylcarboxylic Acid

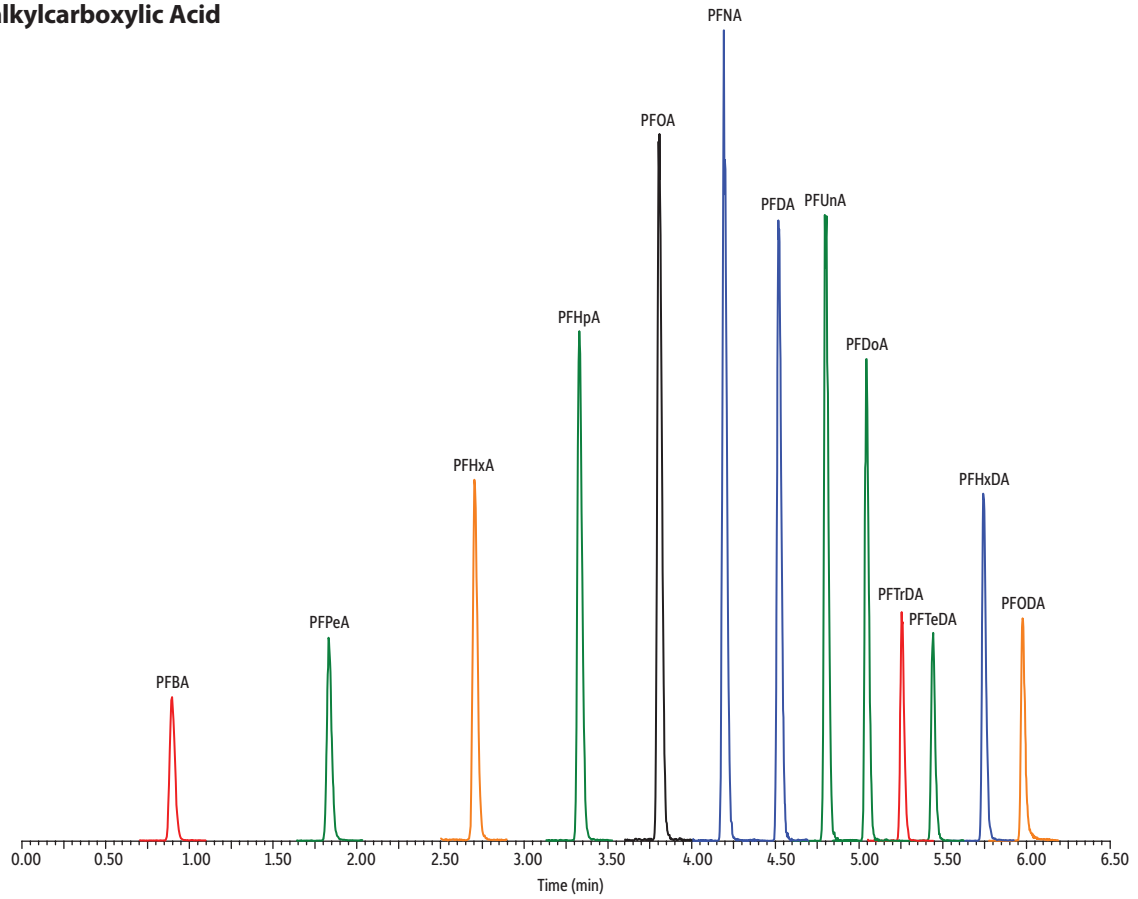

Perfluoroalkyl Sulfonate

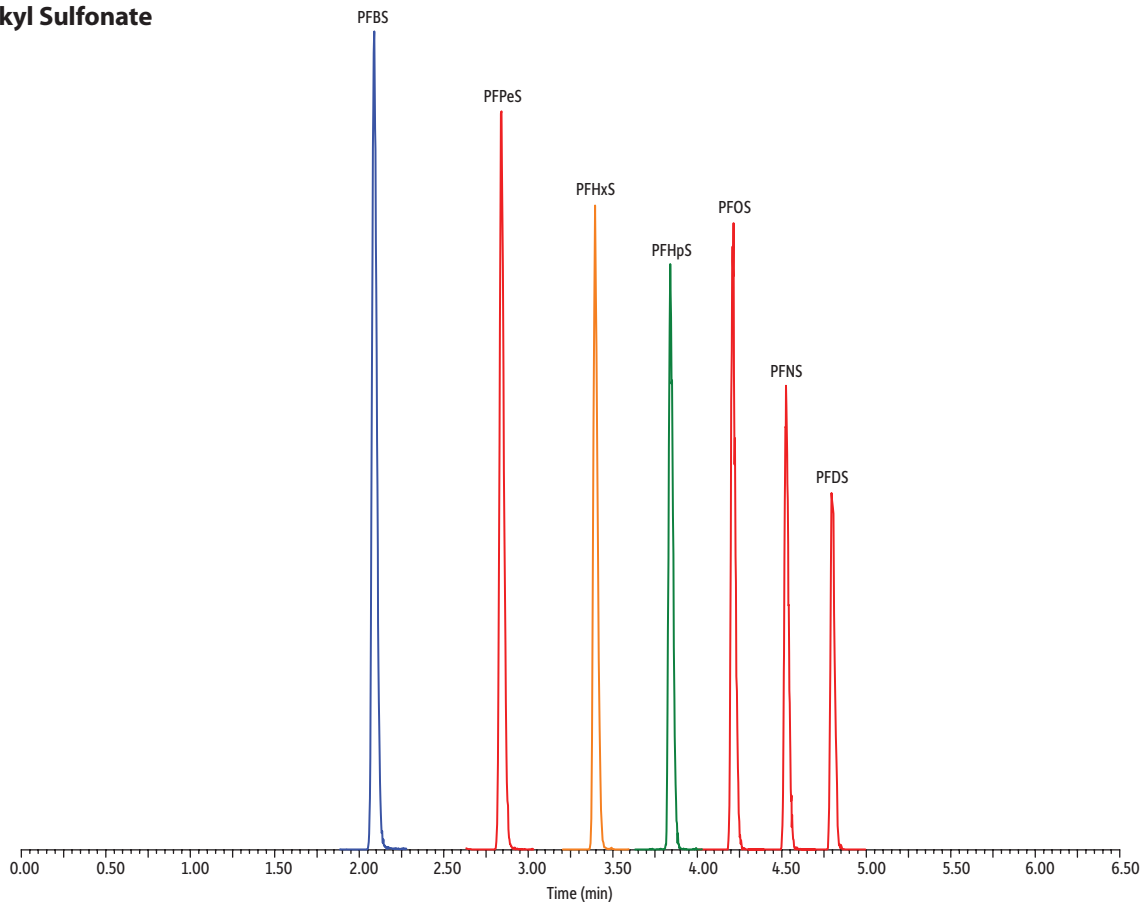

Fluorotelomer Sulfonate
Perfluoroalkyl Sulfonamidoethanol
Perfluoroalkyl Sulfonamidoacetic Acid

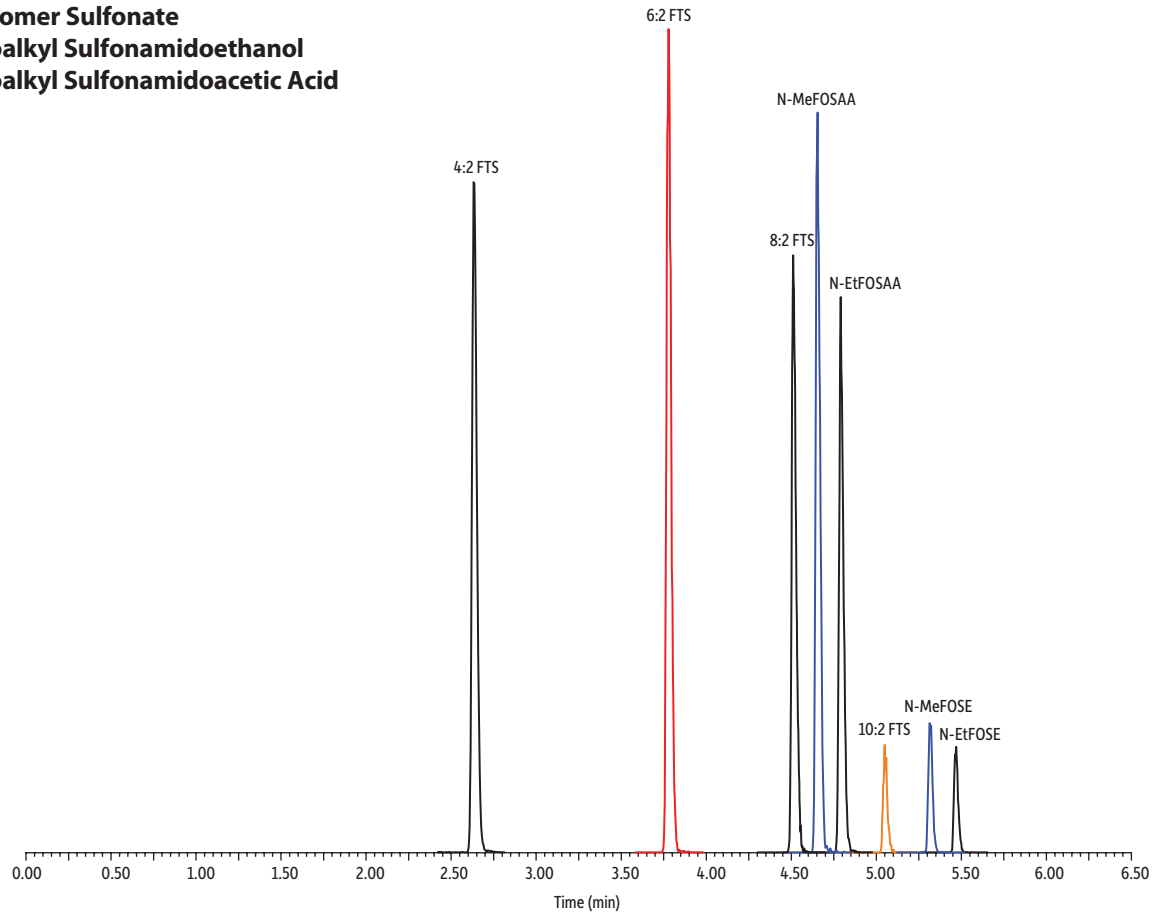

Perfluoroalkyl Sulfonamide
Fluorotelomer Carboxylic Acid

