



## Aromatic Amino Acids on Ultra Aqueous C18 (Fast Analysis)

MZ-Analysetechnik GmbH  
Barcelona-Allee 17 • D-55129 Mainz  
Tel. +49 6131 880 96-0  
Fax +49 6131 880 96-20  
e-mail: info@mz-at.de  
www.mz-at.de

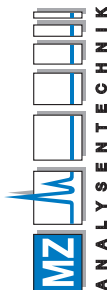

AUTHORIZED DISTRIBUTOR

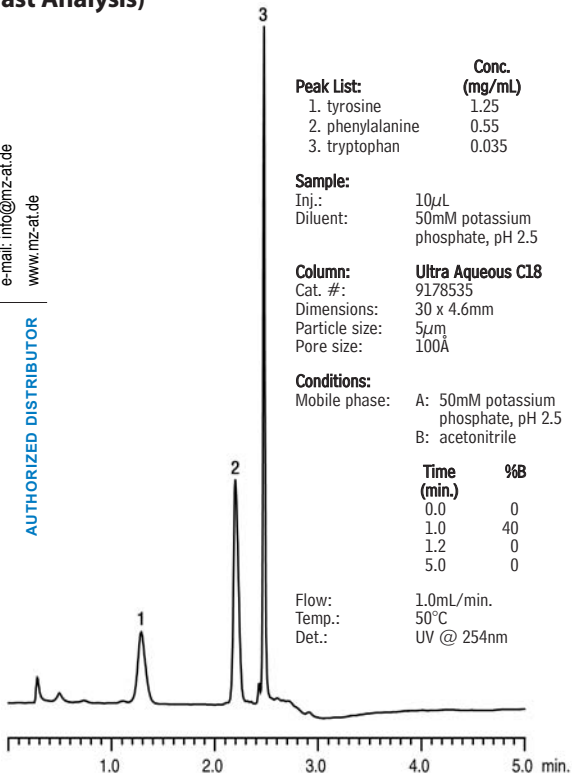

Peak List:	Conc. (mg/mL)
1. tyrosine	1.25
2. phenylalanine	0.55
3. tryptophan	0.035

**Sample:**  
Inj.: 10 $\mu$ L  
Diluent: 50mM potassium phosphate, pH 2.5

**Column:** Ultra Aqueous C18  
Cat. #: 9178535  
Dimensions: 30 x 4.6mm  
Particle size: 5 $\mu$ m  
Pore size: 100Å

**Conditions:**  
Mobile phase: A: 50mM potassium phosphate, pH 2.5  
B: acetonitrile

Time (min.)	%B
0.0	0
1.0	40
1.2	0
5.0	0

Flow: 1.0mL/min.  
Temp.: 50°C  
Det.: UV @ 254nm

LC\_0307