




Proteomix-RP Polymer Based Reverse Phase

Proteomix RP-300, 10 μ m, 300 A 4.6 x 150 mm

Part Number: 469300-4615

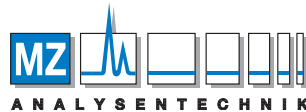

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0, Fax +49 6131 880 96-20
e-mail: info@mz-at.de, www.mz-at.de

Contact Us:

Website: Sepax-Tech.com

Phone: 1-877-SEPAX-US

Email: sales@sepax-tech.com

LinkedIn: [Sepax-Technologies-Inc](https://www.linkedin.com/company/Sepax-Technologies-Inc)



Separation of DNA by Proteomix RP Column

Column: Proteomix RP-300, 10 μm , 300 A 4.6 x 150 mm; **Part Number:** 469300-4615

Mobile Phase: A: 0.1% TEAA, pH=7.0; B: ACN; **Flow Rate:** 1.0 mL/min; **Gradient:** 0% - 12% - 30%B (0 - 30 - 50 min); **Wavelength:** 260 nm; **Column Temperature:** RT; **Sample:** DNA (21 bp, MW ~ 6,000) (1.0 mg/mL)

